

CADD™

**CADD®-Solis VIP
Ambulatory Infusion Pump**

Operator's Manual
Model 2120

smiths medical

The CADD®-Solis VIP (variable infusion profile) Ambulatory Infusion Pump is designed to promote patient care and safety for a variety of adult and pediatric patients, and clinical care areas.

This manual describes only the CADD®-Solis VIP Ambulatory Infusion Pump. Refer to the CADD®-Solis Medication Safety Software installation guide and online help for information specific to that program.

This pump can be programmed with a pump protocol configuration consisting of a therapy, qualifier, and drug. The pump can deliver medication via patient controlled analgesia (PCA), continuous, intermittent, variable stepped rate, and tapered infusions.

This manual is intended for use by clinicians and system administrators only. Do not permit patients to have access to it. The pump has three security levels designed to limit overall access to certain pump features. Only disclose the pump's security codes to those who are authorized. Access to the pump key should also be restricted.

The issue date of this Operator's Manual is included on the back cover. In the event one year has elapsed between the issue date and product use, contact Smiths Medical to see if an updated version of this manual is available.

Technical Assistance

If you have comments or questions concerning the operation of the CADD®-Solis VIP Ambulatory Infusion Pump, call the number given below. When calling, specify the pump software version number. This information is located in the device information report (See *Device Information* on page 98 for more information).

Smiths Medical ASD is available to help with the programming and operation of the CADD®-Solis VIP Ambulatory Infusion Pump.

Smiths Medical ASD, Inc.

St. Paul, MN 55112 USA
1 800.258.5361 (USA)
1 214.618.0218
www.smiths-medical.com

European Distribution:

Smiths Medical International Ltd.

WD24 4LG, UK
+44 (0) 1923 246434

Australian Distribution:

Smiths Medical Australasia Pty. Ltd.

61 Brandl Street
Eight Mile Plains, QLD 4113, Australia
+61 (0) 7 3340 1300

Read this entire operator's manual before operating the CADD®-Solis VIP ambulatory infusion pump. Failure to properly follow warnings, cautions, and instructions could result in death or serious injury to the patient.

Contraindications

- The pump is not to be used for delivery of blood or cellular blood products, as blood and cellular blood products will be damaged by the pumping mechanism.
- This pump is not to be used in any intra-articular space infusion.

Warnings

- This operator's manual should be used by authorized clinicians and system administrators only. Do not permit patients to have access to this manual, as the information contained allows complete access to all programming and operating functions.
- To avoid explosion hazard, do not use the pump in the presence of flammable anesthetics or explosive gases.
- For patients likely to be adversely affected by unintended operations and failures, including interrupted medication or fluid delivery from the device, close supervision and provision for immediate corrective action should be provided in order to assure minimum medication delivery interruption. Pump failure suspends medication delivery, and unintended pump operations could lead to a variety of consequences for the patient.
- If the pump is used to deliver life-sustaining medication, an additional pump must be available, and close supervision and provision for immediate corrective action should be provided to assure minimum medication delivery interruption in the event of a pump failure. Pump failure suspends medication delivery.
- The user should ensure that the performance offered by the pump is fit for the intended use and that the pump is not used in any way or for any purpose other than its intended use.
- If the pump is dropped or hit, inspect it for damage. Do not use a pump that is damaged or not functioning properly. Contact Smiths Medical Customer Service to return a pump for service.
- Use of a syringe with the CADD® administration set may result in UNDERDELIVERY of medication. Syringe function can be adversely affected by variations in plunger dimension and lubricity, which can result in greater force required to move the syringe plunger. A syringe plunger loses lubrication as it ages and as a result, the amount of underdelivery increases, which could on occasion, be significant. Therefore, the type of medication and delivery accuracy required must be considered when using a syringe with the CADD®-Solis VIP ambulatory infusion pump.

The volume remaining in the syringe must be regularly compared to the displayed values on the pump (for example, Reservoir Vol. and Total Given) in order to determine whether underdelivery of medication is occurring. If necessary, appropriate action must be taken.

- Do not administer drugs to the epidural space or subarachnoid space unless the drug is indicated for administration to those spaces.
- To prevent the infusion of drugs that are not indicated for epidural space or subarachnoid space infusion, DO NOT use administration sets that incorporate injection sites.
- If a CADD™ medication cassette reservoir, CADD® extension set, or CADD® administration set is used for epidural space or subarachnoid space drug delivery, it is strongly recommended that it be clearly differentiated from those used for other routes of infusion (for example, by color coding or other means of identification).
- If the air detector is turned off, the pump does not detect air in the fluid path. It is recommended that you periodically inspect the fluid path and remove any air to prevent air embolism.
- When the upstream occlusion sensor is turned off, the pump does not detect occlusions upstream between the pump and the reservoir. Periodic inspection of the fluid path for kinks, a closed clamp,

or other upstream obstructions is recommended. Upstream occlusions may result in underdelivery of medication.

- Do not disclose pump security codes or any other information that would allow the patient or unauthorized clinician complete access to programming and operating functions.
- Do not leave the pump unattended while unlocked. All programming functions are accessible while the pump is unlocked.
- The manual program mode does not contain programming limits. Carefully review each parameter to ensure it accurately matches the prescription.
- Always carefully review the program on the pump after it has been programmed to verify that the pump is programmed correctly.
- The remote dose cord is for *patient use only*. Operation by anyone other than the patient may cause overdelivery of medication.
- Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected.
- Exercise care when using the clinician bolus function. Because there are no limits to the frequency of delivering a bolus, and because the bolus amount may be set as high as 50 mL (or the mg or mcg equivalent), you should not permit the patient or unauthorized clinician to become familiar with the procedure for giving a clinician bolus.
- To prevent the patient from accessing the clinician bolus function, do not let the patient know the clinician or administrator security codes.
- Never leave the pump unattended while on the clinician bolus edit screen. You must deliver the programmed value, or cancel and leave the screen.
- Do not use rechargeable NiCd or nickel metal hydride (NiMH) batteries. Do not use carbon zinc (“heavy duty”) batteries. They do not provide sufficient power for the pump to operate properly.
- Always have new batteries available for replacement. If power is lost, nondelivery of drug occurs.
- A rechargeable battery pack that has reached the end of its useful life must be replaced with either another CADD®-Solis rechargeable battery pack or with 4 AA batteries.
- There is no pump alarm to alert users that a battery has not been properly installed. An improperly installed battery could result in loss of power and nondelivery of drug.
- Always check the battery compartment for fluid or debris before inserting the batteries, and do not allow any fluid or debris to fall into the battery compartment. Fluid or debris in the battery compartment may damage the battery contacts and could result in loss of power and nondelivery of drug.
- If the pump is dropped or hit, the battery door may become broken or damaged. Do not use the pump if the battery door is damaged because the batteries will not be properly secured. This may result in loss of power and non-delivery of drug.
- Follow the instructions for use provided with the CADD™ medication cassette reservoir, CADD® extension set, or CADD® administration set, paying particular attention to all warnings and cautions associated with their use.
- Per general rules of safe practice, always clamp the tubing before removing the cassette from the pump. Removing the cassette without closing the clamp could potentially cause unregulated gravity infusion.
- Attach the cassette properly. A detached or improperly attached cassette could result in unregulated gravity infusion of medication from the reservoir or a reflux of blood.

If you are using a CADD® administration set or CADD™ medication cassette reservoir that does not have the Flow Stop feature, you must use a CADD® extension set with anti-siphon valve, or a CADD®

administration set with either an integral or add-on anti-siphon valve to protect against unregulated gravity infusion that can result from an improperly attached cassette.

- Do not prime the fluid path with the tubing connected to a patient as this could result in overdelivery of medication or air embolism.
- To prevent air embolism, ensure that the entire fluid path is free of all air bubbles before connecting to the patient.
- System delivery inaccuracies may occur as a result of back pressure or fluid resistance, which depends upon drug viscosity, catheter size, and extension set tubing (for example, microbore tubing), and placing the infusion reservoir and/or pump above or below the level of the patient. System delivery inaccuracy may result in under or overdelivery of medication.
- The use of power supplies and a remote dose cord other than those listed in the electromagnetic emissions declaration may result in increased emissions or decreased immunity of the pump.
- The pump should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, verify normal operation of the pump in the configuration in which it is to be used.
- There are potential health hazards associated with improper disposal of batteries, electronics, and contaminated (used) reservoirs and extension sets. Dispose of used batteries, reservoirs, extension sets, and other used accessories, or a pump that has reached the end of its useful life, in an environmentally safe manner, and according to any regulations that may apply.

Cautions

- Do not operate the pump at temperatures below 2°C (36°F) or above 40°C (104°F) to avoid damaging the electronic circuitry.
- Do not store the pump at temperatures below -20°C (-4°F) or above 60°C (140°F) to avoid damaging the electronic circuitry. Do not store the pump with a CADD™ medication cassette reservoir or CADD® administration set attached.
- Do not expose the pump to humidity levels below 20% or above 90% relative humidity to avoid damaging the electronic circuitry.
- CADD® pumps are sealed units. A broken or damaged seal will therefore be considered conclusive evidence that the pump has been misused and/or altered, which voids any and all warranties. All service and repair of CADD® pumps must be performed by Smiths Medical or its authorized agents.
- Use only Smiths Medical accessories that are specified for use with the CADD®-Solis ambulatory infusion pump, as other brands may adversely affect pump operation.
- Inspect the AA batteries for damage or wear to the metal or plastic insulation prior to use, or after the pump has been dropped or hit. Replace the batteries if any damage is noted.
- Do not store the pump for prolonged periods with the batteries installed. Battery leakage could damage the pump.
- If the power up results in an error message indicating that the protocol library was lost, do not proceed with using the pump. Follow your facility's procedures for downloading protocol libraries.
- If you are using a CADD™ medication cassette reservoir in which the medication is frozen, thaw at room temperature only. Do not heat in a microwave oven as this may damage the product and cause leakage.
- If delivery of an infusion is affected by a time or date change, an alarm message appears and must be confirmed.
- Do not immerse the pump in cleaning fluid or water. Do not allow solution to soak into the pump, accumulate on the keypad, or enter the battery compartment, USB port, remote dose cord jack, or power jack areas. Moisture buildup inside the pump may damage the pump.
- Do not clean the pump with acetone, other plastic solvents, or abrasive cleaners, as damage to the pump may occur.
- The pump should *not* be directly irradiated by therapeutic levels of ionizing radiation due to the risk of permanent damage to the electronic circuitry. The best procedure to follow is to remove the pump from the patient during therapeutic radiation sessions or diagnostic levels of radiographic and fluoroscopic radiation. If the pump must remain in the vicinity during a diagnostic or therapy session, it should be shielded, and its ability to function properly should be confirmed following treatment.
- Do not expose the pump directly to ultrasound, as permanent damage to the electronic circuitry may occur.
- Magnetic fields produced by magnetic resonance imaging (MRI) equipment may adversely affect the operation of the pump. Remove the pump from the patient during MRI procedures and keep it at a safe distance from magnetic energy.
- Use of this pump on patients monitored by electronic equipment may cause artifactual interference. As with all electronic equipment, electrical artifacts which affect the performance of other equipment, such as ECG monitors, can occur. The user should check the correct function of the equipment prior to use.
- Do not use the pump in hyperbaric chambers as they affect how the pump works and may also cause damage to the pump.

Table of Contents

Technical Assistance 2
 Contraindications 3
 Warnings 3
 Cautions 6

General Description

Introduction 9
 Indications 9
 Contraindications 9
 Epidural/Subarachnoid Administration 9
 Analgesics 9
 Anesthetics 9
 Symbols 10
 Pump Diagram 12
 Pump Components, Ports, and Connectors 13
 Indicator Lights 13
 Display with Backlighting 13
 Keypad 13
 Power Switch 14
 Power Jack 14
 USB Port 14
 Remote Dose Cord Jack 14
 Battery Compartment 14
 Cassette Latch 14
 Cassette/Keypad Lock 14
 Air Detector 14
 Downstream Occlusion Sensor 14
 Upstream Occlusion Sensor 14
 Accessories 15
 CADD®-Solis Medication
 Safety Software 15
 AC Adapter 15
 Rechargeable Battery Pack 15
 Remote Dose Cord 15
 Polemount Bracket 15
 Polemount Bracket Adapter 15
 Polemount Swivel 15
 Pump Key 15
 Pouches and Backpacks 15
 LockBox 15
 Pump Screens and Display 16
 Pump Screens 16
 Using Color 17
 Security Settings 18
 Security Level Tables 18
 Autolock 21
 Entering Security Codes 21
 Customized Security Settings 22
 Keypad Security 23
 Keypad Code 23
 Clinician Code 24
 Administrator Code 24
 Manual Programming 25

Pump Programming

Protocols and Protocol Libraries 26
 New Protocol Selection 26
 Single Protocol Download 26
 Protocol Selection from a Library 27
 Start New Patient 27
 Start New Protocol Same Patient 27
 Manually Programmed Protocols 28
 Delayed Start/Next Dose Start Time 31
 Delivery Hard and Soft Limits 33
 Delivery Settings 34
Editing Delivery Settings 35
 Review Pump Settings 36
 PCA Delivery Mode 37
 Home Screen 37
 Programming Screens 38
 Programming Example 39
 Give a Clinician Bolus 43
 Start a PCA Dose 45
 Stopping a Clinician Bolus or PCA Dose 46
 Continuous Delivery Mode 47
 Home Screen 47
 Programming Screens 47
 Programming Example 48
 Intermittent Delivery Mode 51
 Home Screen 51
 Programming Screens 52
 Programming Example 53
 Stopping the Pump During an Infusion 58
 Step Delivery Mode 59
 Home Screen 59
 Programming Screens 60
 Programming Example 62
 Starting Each infusion 67
 Stopping and Restarting the Pump
 During an Infusion 67
 Step Up 67
 Step Down 68
 Taper Delivery Mode 69
 Home Screen 69
 Programming Screens 70
 Programming Example 72
 Starting Daily Infusion 76
 Stopping and Restarting the Pump
 During an Infusion 77
 Taper Down Now 77

Operating the Pump

Batteries	78
Install the Batteries or Battery Pack	79
Replace the Battery Door	79
Power Up	80
Power Down	80
Cassettes	81
Remove a Cassette	81
Attach a Cassette	82
Prime Tubing	83
Prime Tubing After Changing a Cassette	83
Prime Tubing, No Cassette Change	84
Start the Pump	85
Stop the Pump	86
Reset Reservoir Volume	87

Tasks and Advanced Tasks

Tasks Menu Overview	88
Display and Sound Settings	89
Backlight Intensity	89
Alarm Volume	89
Sound Theme	90
Key Beep On/Off	91
Numeric Format	91
Time and Date	92
Current Time	92
Current Date	93
Time Format	93
Date Format	94
Daylight Savings Time	94
Reports	95
Total Given	95
Given and PCA Dose Counters	96
PCA Dose Graph	96
Delivery History and Pie Chart	97
Delivery Log	97
Event Log	98
Protocol Library Summary	98
Device Information	98
Advanced Tasks Menu Overview	99
Patient Permissions	100
Priming Security On/Off	100
Delayed Start Security On/Off	101
Air and Occlusion Settings	102
Air Detector On/Off	102
Air Detector Sensitivity	103
Upstream Sensor On/Off	103
Downstream Sensor Sensitivity	104

Alarm Settings	105
Infusion Alerts On/Off	105
Pump Stopped Alarm Type	105
Reservoir Low Trip Point	106
Reservoir Low Alarm Type	106
Reservoir Empty Alarm Type	107
PM Reminder On/Off	107
PM Reminder	108
PM Interval	108
Reset to Factory Settings	109

References and Troubleshooting

Alarms and Messages	110
Types of Alarms	110
System Fault Alarm	110
High Priority Alarm	110
Medium Priority Alarm	110
Low Priority Alarm	111
Informational Message	111
Alarm Help Screens	111
Troubleshooting	112
Alarms and Messages,	
Alphabetical List	112
Cleaning the Pump and Accessories	119
Radiation and Magnetic Resonance	
Imaging (MRI)	120
Standards used in development	
of the pump	121
PCA Delivery Mode Scroll Ranges	123
Military Time	124
Specifications (Nominal)	125
General Pump Specifications	125
Delivery Specifications	130
Common Delivery Specifications	130
PCA Delivery Specifications	131
Continuous Delivery	
Specifications	132
Intermittent Delivery	
Specifications	132
Step Delivery Specifications	133
Taper Delivery Specifications	133
Electromagnetic Emissions and	
Immunity Declarations	134
Collect Separately	136
Menu Maps	137
Tasks Menus	137
Advanced Tasks Menus	139
Default Factory Settings	140
Manual Mode Initial Settings	141
Accuracy Test Results	143
Limited Warranty	145

Index	146
--------------------	-----

General Description

Introduction

The CADD[®]-Solis VIP ambulatory infusion pump provides measured drug therapy to patients in hospital or outpatient settings. Therapy should always be overseen by a physician or a certified, licensed healthcare professional. As appropriate, the caregiver and patient should be instructed in using the pump.

Indications

The CADD[®]-Solis VIP ambulatory infusion pump is indicated for intravenous, intra-arterial, subcutaneous, intraperitoneal, perineural, surgical site, epidural space, or subarachnoid space infusion.

- **PCA** (patient-controlled analgesia) delivery is used for therapies that require a continuous rate of infusion, patient-controlled demand doses, or both, such as patient-controlled analgesia.
- **Continuous** delivery allows the infusion of drug/fluid at a constant, programmed rate.
- **Intermittent** delivery allows the infusion of a specific volume of drug/fluid at a regular, programmed interval.
- **Step** delivery allows an incremental increase in infusion rate to a specified maximum infusion rate for a specified total infusion volume.
- **Taper** delivery allows a plateau rate of infusion with the option of tapering at the beginning and/or end and has a programmable KVO rate at the end of the infusion.

Contraindications

- The pump is not to be used for delivery of blood or cellular blood products, as blood and cellular blood products will be damaged by the pumping mechanism.
- This pump is not to be used in any intra-articular space infusion.

Epidural/Subarachnoid Administration

The selected drug must be used in accordance with the indications included in the package insert accompanying the drug. Administration of any drug by this pump is limited by any warnings, precautions, or contraindications in the drug labeling.

Analgesics

Administration of analgesics to the epidural space is limited to use with indwelling catheters specifically indicated for either short-term or long-term drug delivery.

Anesthetics




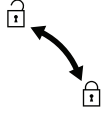
























Administration of anesthetics to the epidural space is limited to use with indwelling catheters specifically indicated for short-term drug delivery.

WARNING:







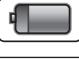



























- **Do not administer drugs to the epidural space or subarachnoid space unless the drug is indicated for administration to those spaces. Drugs not intended for epidural or subarachnoid space infusion could result in serious patient injury or death.**
 - **To prevent the infusion of drugs that are not indicated for epidural space or subarachnoid space infusion, DO NOT use administration sets that incorporate injection sites. The inadvertent use of injection sites for infusion of such drugs may cause serious patient injury or death.**
 - **If a CADD™ medication cassette reservoir, CADD[®] extension set, or CADD[®] administration set is used for epidural space or subarachnoid space drug delivery, it is strongly recommended that it be clearly differentiated from those used for other routes of infusion (for example, by color coding or other means of identification). Drugs not intended for epidural or subarachnoid space infusion could result in serious patient injury or death.**
-

Symbols

Symbols on the Pump and Labels


	Caution		Power on/off switch
	Consult instructions for use		Cassette unlock/lock
	Catalog number		Soft key
	Serial number		Up button
	Date of manufacture		Down button
	Authorized representative in the European Community		Select button
	Australian representative		Stop/Start button
	Temperature limitation		PCA Dose button
	Humidity limitation		AA battery location, positive terminal faces up
	Atmospheric pressure limitation		AA battery location, negative terminal faces up
	Type CF equipment		AA battery location, positive terminal faces up
	Class II equipment		AA battery location, negative terminal faces up
IPX4	Splashproof—water splashed against the pump housing will have no harmful effects.		USB port
	Collect separately		Direct current (power jack)
			Remote dose cord jack

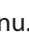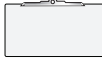Symbols on the Pump Display


	Reservoir volume		PCA delivery mode
	Charge level of rechargeable battery pack		Continuous delivery mode
	Charge level of rechargeable battery pack. AC adapter connected.		Intermittent delivery mode
	Charge level of AA batteries		Step delivery mode
	Charge level of AA batteries. AC adapter connected.		Taper delivery mode
	No battery installed, AC power only		Keypad is locked
	Incompatible battery		Keypad is unlocked
	Incompatible battery. AC Adapter connected.		Appears next to a parameter that has been reviewed and accepted.
	Pump status is Started or Running.		On edit screens. Indicates current parameter value. Press  or  to scroll up or down to edit the value.
	Pump status is Paused.		On edit screens with a menu of options. Indicates which setting is being selected.
	Pump status is Stopped.		The requested action can not be performed.
	There are more items to see in the work area. Press  or  to scroll up or down.		More data is available. On PCA dose report only.
	Item highlighted in the work area is at the top of the menu. Press  to scroll down.		Review screen
	Item highlighted in the work area is at the bottom of the menu. Press  to scroll up.		Saving

Pump Diagram

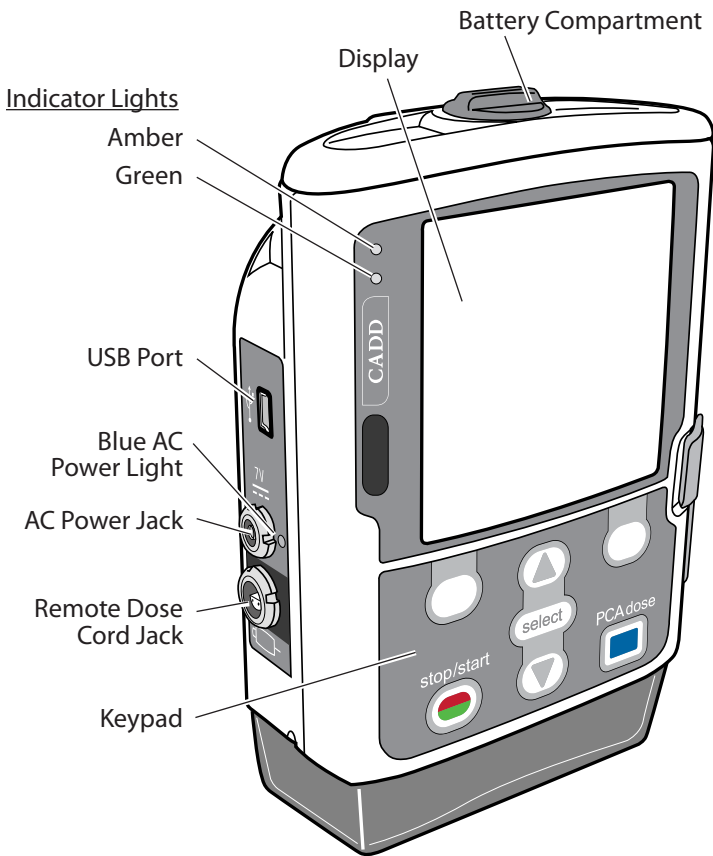Front View



Rear View

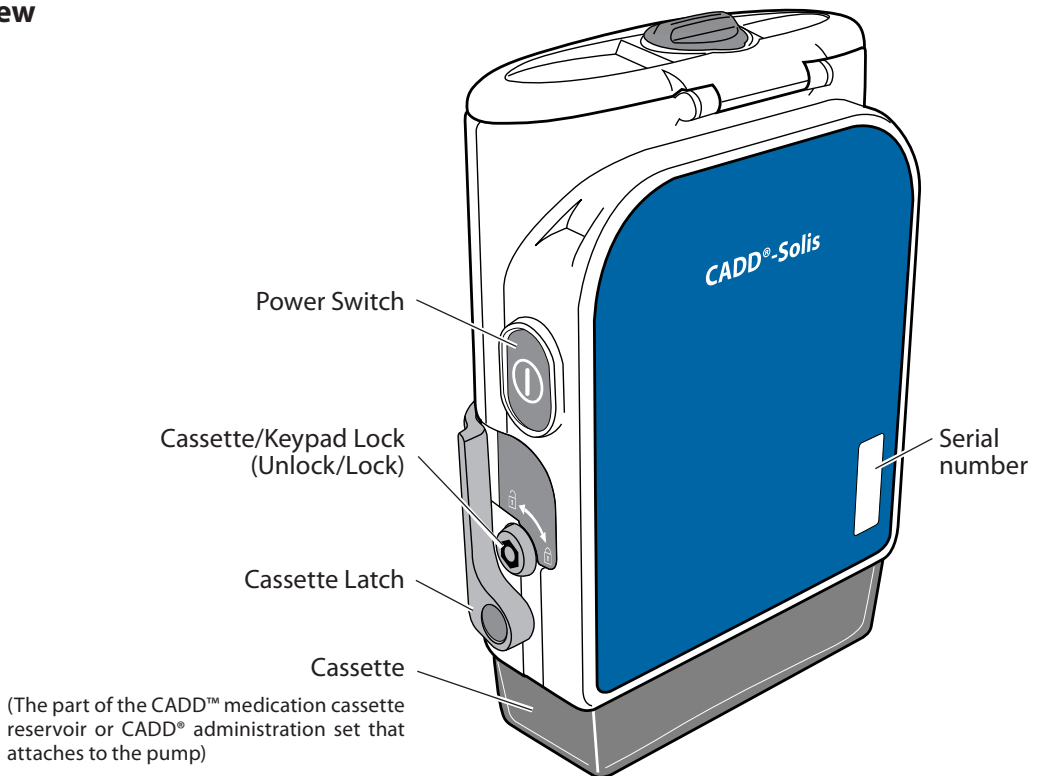

Pump Components, Ports, and Connectors

Indicator Lights

When the pump is powered, one or both of the indicator lights flash.

Green: The green light flashes to indicate that the pump is running and delivering fluid as programmed.

Amber: The amber light flashes when the pump is stopped, an alarm condition exists, or the battery or the reservoir volume is low. It stays on continuously when the pump is inoperable. The display briefly describes the alarm condition when the amber light is flashing.



NOTE: At times both lights may flash. This indicates that the pump is running, but there is a condition that needs attention (for example, low battery or low reservoir volume).

Display with Backlighting

The liquid crystal display (LCD) shows programming information and messages. Backlighting helps keep the display visible in low light. In this manual, “display” is synonymous with display panel or LCD.

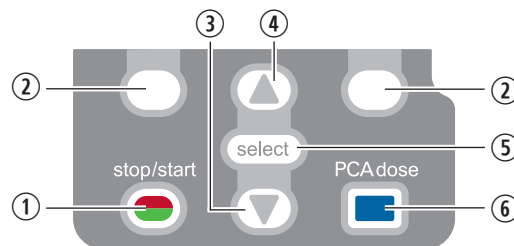
After a period of time in which no keys are pressed, the backlighting turns off and the display goes blank to save battery power (except during an alarm or when the AC adapter is in use). To turn the display back on, press any key except the PCA dose key when in PCA mode.

NOTE:

- When the display is blank, you can determine that the pump is powered by observing either the green or amber (or both) indicator lights periodically flashing.
- If you press stop/start , the display reappears with a message asking if you wish to start or stop the pump.
- If you press PCA dose  while the pump screen is blank, the pump delivers a PCA dose, if that option is available.

Keypad

The keys on the keypad are described below. A key beeps when pressed if it is operable in the current state of the pump. However, the keys do not beep if the key beep function has been turned off in the protocol or the administrator settings.



- ① Starts and stops pump delivery.
- ② Referred to as “soft keys.” Allow you to answer a question on the pump’s display. For example, the screen above this key may display “Yes,” in which case pressing this key gives the question displayed on the screen an answer of “Yes.” Also allows you to navigate through some of the pump screens (for example, canceling an action, or backing out of an open screen).
- ③ Allows navigation through the menus on the pump, scrolling down, or decreasing a value.
- ④ Allows navigation through the menus on the pump, scrolling up, or increasing a value.
- ⑤ Used to select a menu item.
- ⑥ Allows the patient to request a PCA dose if the remote dose cord is not connected, and if the PCA dose option is available. If the remote dose cord is connected, this key is inactive.

Power Switch

Turns the pump on or off. Press and hold the switch to turn the pump on. Press the switch to turn the pump off, and confirm that you want to power down by selecting **Yes**.

Power Jack

You may plug the AC adapter into the power jack. When the AC adapter is plugged in, the blue power light turns on. This light is on regardless of the pump's on or off status. See *AC Adapter* on page 15 for more information.

USB Port

A mini-B USB cord may be attached to the USB port for communication with the CADD[®]-Solis Medication Safety Software.

Remote Dose Cord Jack

Used for attaching the remote dose cord. See *Remote Dose Cord* on page 15 for more information.

Battery Compartment

Four AA batteries or the rechargeable battery pack fit into this compartment. The batteries serve as the primary source of power, or as a backup when the AC adapter is in use. See *Batteries* on page 78 for additional information.

Cassette Latch

Used to attach the cassette to the pump. When the pump is turned on, it detects whether the cassette is latched properly. Delivery stops and an alarm occurs if the cassette becomes unlatched. See *Attach a Cassette* on page 82 and *Remove a Cassette* on page 81.

Cassette/Keypad Lock

Secures the cassette to the pump using the pump key provided. The cassette latch must be latched before it can be locked. The cassette/keypad lock can be configured to unlock only the cassette latch or to unlock the cassette latch as well as the keypad. This is configured by the CADD[®]-Solis system administrator. See *Security Settings* on page 18.

Air Detector

The air detector may be turned on or off, depending on facility or therapy requirements (see *Air Detector On/Off* on page 102). If air is detected in the part of the tubing that passes the air detector sensor, an alarm sounds and delivery stops. If an air detector is not required, it may be turned off.

WARNING: If the air detector is turned off, the pump does not detect air in the fluid path. It is recommended that you periodically inspect the fluid path and remove any air to prevent air embolism, which could result in serious patient injury or death.

Downstream Occlusion Sensor

The pump contains a downstream occlusion sensor. When a downstream occlusion (between the pump and the patient) is detected, an alarm sounds, delivery stops, and the display shows “Downstream occlusion. Clear occlusion between pump and patient”.

Upstream Occlusion Sensor

The pump contains an upstream occlusion sensor. This feature may be turned on or off (see *Upstream Sensor On/Off*, on page 103). When an upstream occlusion (between the pump and the reservoir) is detected, an alarm sounds, delivery stops, and the display shows “Upstream occlusion. Clear occlusion between pump and reservoir”.

WARNING: When the upstream occlusion sensor is turned off, the pump does not detect occlusions upstream between the pump and the reservoir. Periodic inspection of the fluid path for kinks, a closed clamp, or other upstream obstructions is recommended. Upstream occlusions may result in underdelivery of medication. If undetected, these occlusions could lead to death or serious injury to the patient.

Accessories

CAUTION: Use only Smiths Medical accessories that are specified for use with the CADD[®]-Solis ambulatory infusion pump, as other brands may adversely affect pump operations.

All CADD[®]-Solis accessories may be obtained through the Customer Service department at Smiths Medical. For detailed instructions and warnings applicable to each accessory, refer to the instructions for use provided with the product.

CADD[®]-Solis Medication Safety Software

The CADD[®]-Solis Medication Safety Software allows you to create and manage protocol libraries and then send them to the pump. See the installation guide and online help for more information.

AC Adapter

The AC adapter can be used as an alternate source of power for the pump and/or to recharge the rechargeable battery pack. The pump requires 4 AA batteries or the rechargeable battery pack to be installed as a backup while using the AC adapter.



Rechargeable Battery Pack

The rechargeable battery pack is an alternative to using 4 AA batteries. The rechargeable battery pack can be recharged with the AC adapter, either inside or outside of the pump.



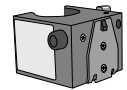
Remote Dose Cord

The remote dose cord can be attached to the pump and provided to the patient as an alternative to pushing the PCA dose  key to request a PCA dose.




Polemount Bracket

The polemount bracket is used along with the polemount bracket adapter to attach the pump to an IV pole.



Polemount Bracket Adapter

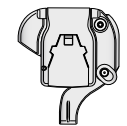The polemount bracket adapter attaches to the back of the pump so that it can be hung from an IV pole or in a lockbox.



Polemount Swivel

The polemount swivel allows the pump to rotate while it is attached to the CADD[®]-Solis polemount bracket. With the swivel attached, the pump can be tilted or rotated 230° for easy viewing of the display. The swivel can not be used in a LockBox.



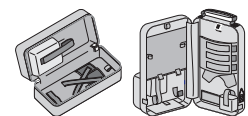
Pump Key

The pump key is used to securely lock a cassette to the pump. It can also be used to lock and unlock the keypad, if allowed by the protocol setting.




Pouches and Backpacks

Several styles of pump pouches and backpacks are available. They can carry a CADD[®]-Solis ambulatory infusion pump with either a CADD[™] medication cassette reservoir (50 mL to 250 mL) or a flexible, plastic IV bag up to 3 liters.



LockBox

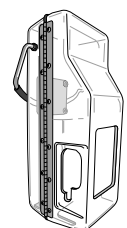
Clear and yellow lockboxes are available to hold both the CADD[®]-Solis ambulatory infusion pump and a CADD[™] medication reservoir. The pump attaches to the lockbox with a polemount bracket adapter and the lockbox can be attached to an IV pole. The lockbox allows for viewing of its contents and access to the keypad, remote dose cord jack, and AC power jack.



Pump Screens and Display

Pump Screens

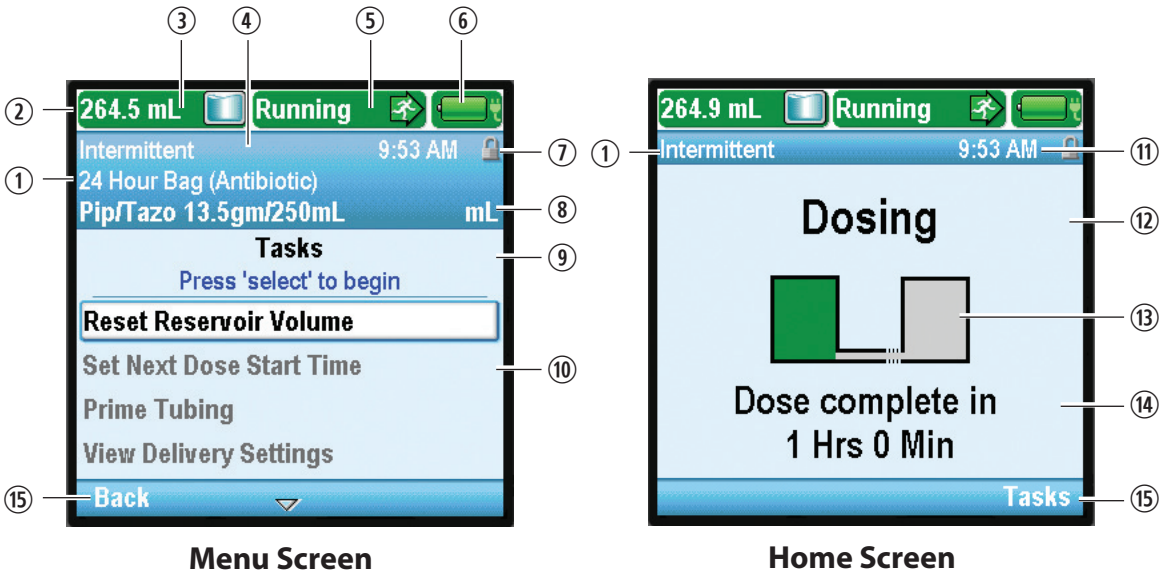The sample screens in this manual are only examples of what might be displayed. The protocols, consisting of therapies, qualifiers, drugs and concentrations, and all associated pump settings in the pump library will be established by your facility.



- ① The therapy, qualifier, and drug in the current protocol. On a home screen, only the therapy is listed.
- ② The status bar shows the status of the pump. It may also display messages and alerts.
- ③ Current reservoir volume.
- ④ The screen color is unique for each therapy.
- ⑤ Delivery status of the pump—stopped or running.
- ⑥ The type of battery in use and the approximate amount of life remaining. Also indicates if an AC adapter is in use.



AA batteries



AA batteries with AC adapter



Rechargeable battery pack



Rechargeable battery with AC adapter

- ⑦ Keypad lock status—locked or unlocked.
- ⑧ Units of measurement and concentration (if applicable) for the drug or solution used in the current protocol.
- ⑨ Screen name and help text, if there is any.
- ⑩ The work area/contents for the displayed screen.
- ⑪ The current time.
- ⑫ The current status of the programmed infusion.
- ⑬ Graphic that identifies which therapy is programmed, and the status of the infusion. See explanations of the different home screens and graphics in *Delivery Settings* beginning on page 34.

⑭ Information indicating upcoming events that are important to the user. The message is mode-specific, and *if the pump is running*, indicates one of the following:

- When a delayed start will begin
- When the reservoir will be empty
- The current rate of delivery
- When the next dose will begin
- When the dose or infusion will be complete

If a new patient or new infusion is programmed but *the pump has not been started*, the message states that the infusion (or dose) will begin when the pump starts.

⑮ Options for navigating the pump. These options change depending on the screen and the functions being performed with the pump.

Using Color

Therapy Screen Colors

The pump display uses color to help you recognize critical information quickly and easily. The unique screen color for each type of therapy is factory-set:

Therapy	Screen Color
PCA	Purple*
Continuous	Gray
Intermittent	Blue
Step	Olive
Taper	Green

* May be customized to yellow.

For PCA therapy only, the screen color may be changed to yellow by using the CADD[®]-Solis Medication Safety Software. Depending on the needs of your pain management program, your facility may choose to relate a purple or yellow screen to a specific protocol, including:

- Route of administration (for example, all epidural protocols), or
- Patient type (for example, all pediatric protocols), or
- Any other hierarchy which fits the needs of your institution.

Refer to your facility's policies and procedures to understand how colors are used to identify your pain management protocols.

Pump Status Colors

The colors green, amber, red, and blue are used to help clinicians and patients recognize pump status. Similar to a traffic control light, green means go, amber indicates caution, and red means stop:

- **Green:** Pump conditions are satisfactory.
- **Amber:** There is a condition to watch, but the current pump conditions are satisfactory. Medium priority alarms and values outside soft limits display in amber.
- **Red:** There is a warning condition that requires immediate attention and infusion has stopped. All high priority and system fault alarms display in red.
- **Blue:** Low priority alarms and informational messages display in blue.

For more information on alarm screens, refer to *Alarms and Messages* on page 110.

Security Settings

Security settings are used to limit patient and unauthorized clinician access to certain programming and operating functions of the pump. Pump functions are protected by 3 different security codes. The security level tables list the functions that are available under each security code.

WARNING:

- **Do not disclose pump security codes or any other information that would allow the patient or unauthorized clinician complete access to programming and operating functions. Improper programming could result in serious patient injury or death.**
 - **Do not leave the pump unattended while unlocked. All programming functions are accessible while the pump is unlocked, and improper programming could result in serious patient injury or death.**
-
-

Security Level Tables

- The keypad code is for clinicians who need to modify and review patient-specific parameters, and manage air and occlusion parameter settings.
- The clinician code allows access to all the functions the keypad code allows. In addition, it allows users to modify most advanced tasks parameters, format the time and date, change delayed start and priming security options, and select new protocols.
- The administrator code allows access to all pump functions and allows users to change protocol ranges for patient-specific parameters, reset the pump to factory defaults, and select manual mode protocols. Its use should be restricted to the CADD[®]-Solis system administrator and key individuals.

Security Levels, Pump Operations and Tasks

Pump Operations and Tasks	Available without a security code		Available with keypad code		Available with clinician code		Available with administrator code	
	Running	Stopped	Running	Stopped	Running	Stopped	Running	Stopped
Stop/Start ☹	✓	✓	✓	✓	✓	✓	✓	✓
PCA Dose ☐ (PCA only)	*		*		*		*	
Reset Reservoir Volume		✓		✓		✓		✓
Taper Down Now (Taper only)	✓		✓		✓		✓	
Taper Down Period of Taper Down Now (Taper only)			✓		✓		✓	
Set Delayed Start		*		✓		✓		✓
Next Dose Start Time (Intermittent only)		*		✓		✓		✓
Prime Tubing		*		✓		✓		✓
View Delivery Settings	✓	✓	✓	✓	✓	✓	✓	✓
Edit Delivery Settings (Therapy type determines which submenus appear)				✓		✓		✓
Display and Sound Settings menu	✓	✓	✓	✓	✓	✓	✓	✓
Backlight Intensity	✓	✓	✓	✓	✓	✓	✓	✓
Alarm Volume	✓	✓	✓	✓	✓	✓	✓	✓
Sound Theme	✓	✓	✓	✓	✓	✓	✓	✓
Key Beep On/Off	✓	✓	✓	✓	✓	✓	✓	✓
Numeric Format						✓		✓
Change Time and Date menu		✓		✓		✓		✓
Current Time		✓		✓		✓		✓
Current Date		✓		✓		✓		✓
Time Format						✓		✓
Date Format						✓		✓
View Reports	✓	✓	✓	✓	✓	✓	✓	✓
View Advanced Tasks (see <i>Security Levels, Advanced Tasks</i> on page 20)	✓	✓	✓	✓	✓	✓	✓	✓

Table Key: Yes ✓ No Availability based on facility protocol *

General Description

General Description

Security Levels, Advanced Tasks

Advanced Tasks	Available without a security code		Available with keypad code		Available with clinician code		Available with administrator code	
	Running	Stopped	Running	Stopped	Running	Stopped	Running	Stopped
Step Down (Step only)	✓	✓	✓	✓	✓	✓	✓	✓
Step Up (Step only)			✓	✓	✓	✓	✓	✓
Give Clinician Bolus (PCA only)					✓		✓	
Patient Permissions						✓		✓
Air and Occlusion Settings				✓		✓		✓
Alarm Settings						✓		✓
Security Settings menu						✓		✓
Keypad Security						✓		✓
Keypad Code						✓		✓
Clinician Code						✓		✓
Admin Code								✓
Manual Programming Security								✓
Start New Patient						✓		✓
Select a manual mode protocol						*		✓
Start New Protocol Same Patient						✓		✓
Select a manual mode protocol						*		✓
Delivery Hard and Soft Limits menu								✓
(Therapy type determines which submenus appear)						✓		✓
Reset to Factory Settings								✓

Table Key: Yes ✓ No Availability based on facility protocol *

Autolock

The CADD[®]-Solis VIP ambulatory infusion pump is designed to meet both safety and usability needs. The autolock feature reduces the chance of unauthorized pump programming. When the keypad is unlocked with a security code and left unlocked, the software automatically locks the keypad.

If the pump is on the home screen, the autolock feature takes effect 30 seconds after the last key press.

Autolock takes longer on programming or task screens when the user typically needs more time to perform an action. Depending on which screen was on the pump last, and if the pump is not alarming, it can take up to 4 minutes after the last key press before the pump reverts to the home screen and autolocks immediately. When the pump is alarming, autolock does not take effect.

NOTE: The keypad can be re-locked by pushing the right soft key twice from the home screen or once from the Tasks or Advanced Tasks menu. As a recommended safety precaution, always manually lock the pump using this feature.





Entering Security Codes

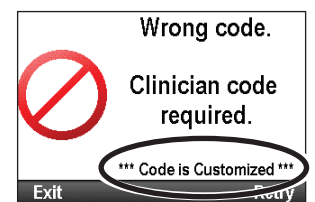
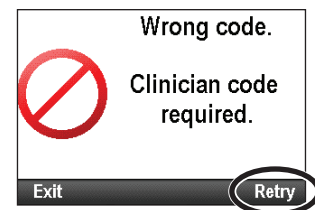
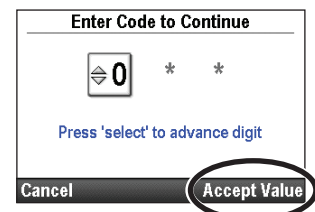
WARNING:

- Do not disclose pump security codes or any other information that would allow the patient or unauthorized clinician complete access to programming and operating functions. Improper programming could result in serious patient injury or death.
- Do not leave the pump unattended while unlocked. All programming functions are accessible while the pump is unlocked, and improper programming could result in serious patient injury or death.

NOTE: Many of the instructions in this manual include steps to unlock the keypad. This step is required only if the keypad is locked.


To enter a security code when prompted:

1. Press  or  to scroll up or down to the correct digit. Press  or **Accept Value** to advance to the next digit. Once the complete code has been entered, press  or **Accept Value**.
2. If an incorrect code is entered, a wrong code error appears. It includes the security level required to access the function. Select **Retry** to enter the code again.




NOTE: If you enter a code that you believe is correct and receive a wrong code message, check the screen to see if the code has been customized. If the code is not customized, it was entered incorrectly and you should retry. If the code is customized, enter the customized code. If you do not know the code, contact your CADD[®]-Solis system administrator.

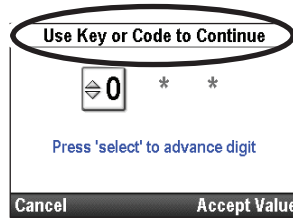- Depending on the level of security required and how the CADD®-Solis system administrator has programmed the security settings, you may be able to use the pump key to unlock the keypad. When accessing a menu or function that requires a security code, 1 of 2 prompts appears, as shown in the sample screens to the right.



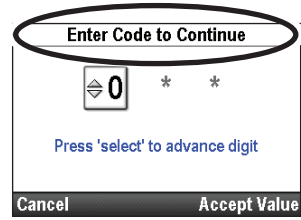
Either the key or a code unlocks the keypad.



Only a code unlocks the keypad.

To unlock the keypad with the pump key, turn the key counterclockwise.

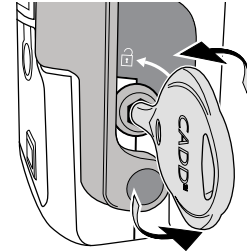NOTE: It is possible for the keypad to be locked while the cassette/keypad lock is unlocked. To use the key to unlock the keypad, first lock the cassette/keypad lock, and then unlock it.




Customized Security Settings

The keypad code, clinician code, and administrator code may be customized by the CADD®-Solis system administrator when setting up protocols with the CADD®-Solis Medication Safety Software, or in the security settings under Advanced Tasks. See your CADD®-Solis system administrator for more information, or to learn the security code you should use if the codes have been customized.


The security settings option in the Advanced Tasks menu allows you to customize the pump security codes. The keypad code, clinician code, and administrator code can all be customized to a unique three-digit number from 001 to 899. If desired, all the codes may be identical. In that case however, all users would have the highest permission level assigned to that code. For example, if the keypad code and clinician codes are the same, each time the code is entered, the user has clinician code permissions. See the security level tables beginning on page 19 for the code permissions.

To access Security Settings:

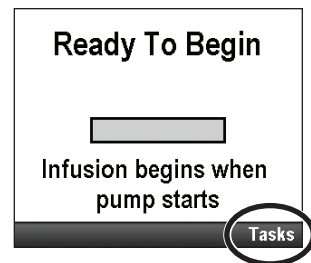
- From the home screen, select **Tasks**.



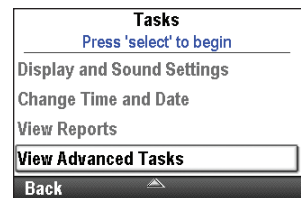- Press **▲** or **▼** to highlight **View Advanced Tasks** and press **select**.



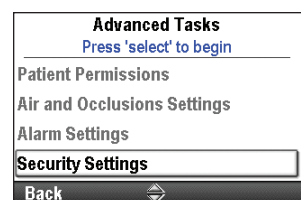- Press **▲** or **▼** to highlight **Security Settings** and press **select**.






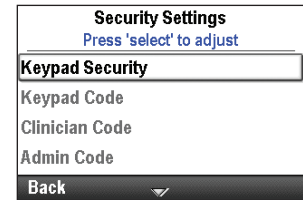
Keypad Security

Under the Keypad Security setting, you may allow the use of the cassette/keypad lock to unlock the keypad. This is done by selecting the Code or Key setting. The use of the key to unlock the keypad allows the same permissions as the keypad code. If the Code Only option is selected, the key will not unlock the keypad.

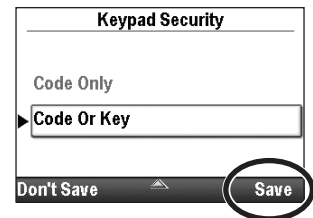
Note: If the Code Or Key setting is selected, using the key to unlock the keypad is equivalent to entering the keypad code. Use of the key allows clinician or administrator code level access *only if* the keypad code is customized to be the same number as the clinician or administrator code.

To set the keypad security:

1. From the Security Settings menu, press  or  to highlight **Keypad Security** and press .



2. Unlock the keypad.
3. Press  or  to highlight **Code Only** or **Code Or Key** and select **Save**.






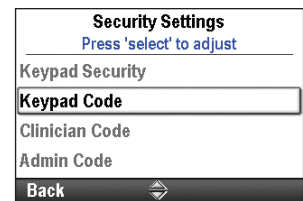
Keypad Code




This screen allows you to create a keypad code that is unique to your facility.

Note: The keypad code changes only for the current protocol. The change remains in effect until another protocol is chosen or the same protocol is chosen by using either the Start New Patient task or the Start New Protocol, Same Patient task.

To set the keypad code:

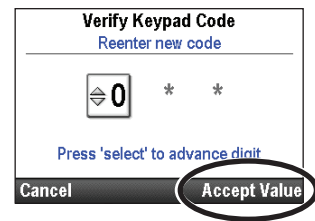
1. From the Security Settings menu, press  or  to highlight **Keypad Code** and press .



2. Unlock the keypad.
3. Choose any 3-digit combination of numbers from 001 to 899. Press  or  to enter the new code. Press  to advance to the next digit. Once the code has been entered, select **Accept Value**.



- Verify the new keypad code by reentering it on the next screen. Once the code has been entered, press **Accept Value**.




Clinician Code

This screen allows you to create a clinician code that is unique to your facility.


Note: The clinician code changes only for the current protocol. The change remains in effect until another protocol is chosen or the same protocol is chosen by using either the Start New Patient task or the Start New Protocol, Same Patient task.

To set the clinician code:

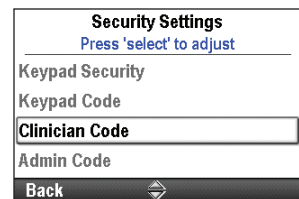
- From the Security Settings menu, press **▲** or **▼** to highlight **Clinician Code** and press **select**.



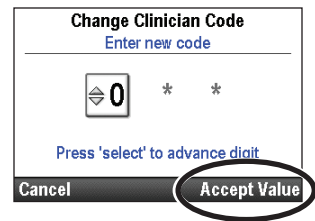- Unlock the keypad.
- Choose any 3-digit combination of numbers from 001 to 899. Press **▲** or **▼** to enter the new code. Press **select** to advance to the next digit. Once the code has been entered, select **Accept Value**.



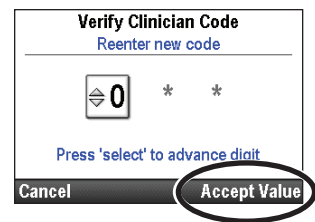- Verify the new clinician code by reentering it on the next screen. Once the code has been entered, press **Accept Value**.




Administrator Code

This screen allows you to create an administrator code that is unique to your facility.

Note: The administrator code changes only for the current protocol. The change remains in effect until another protocol is chosen or the same protocol is chosen by using either the Start New Patient task or the Start New Protocol, Same Patient task.

To set the administrator code:

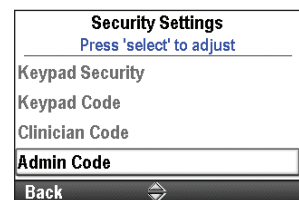
- From the Security Settings menu, press **▲** or **▼** to highlight **Admin Code** and press **select**.



2. Unlock the keypad.
3. Choose any 3-digit combination of numbers from 001 to 899. Press ▲ or ▼ to enter the new code. Press select to advance to the next digit. Once the code has been entered, select **Accept Value**.



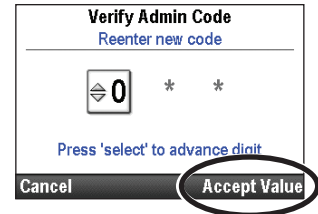4. Verify the new administrator code by reentering it on the next screen. Once the code has been entered, press **Accept Value**.




Manual Programming


The manual programming security option allows you to choose either the clinician or administrator security code to allow access to the manual program protocol.

To select the code for manual programming:

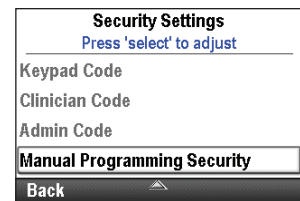1. From the Security Settings menu, press ▲ or ▼ to highlight **Manual Programming Security** and press select.



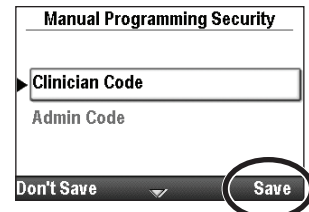2. Unlock the keypad.
3. Press ▲ or ▼ to highlight the desired code and select **Save**.



Pump Programming

Protocols and Protocol Libraries

A **protocol** is a collection of medication delivery and other pump settings, including administrator settings for the CADD®-Solis VIP ambulatory infusion pump, a therapy, qualifier, and drug combination, and other delivery parameters. A collection of protocols in CADD®-Solis Medication Safety Software form a **protocol library**. The CADD®-Solis system administrator creates, maintains, and loads the libraries onto the CADD®-Solis VIP ambulatory infusion pump.

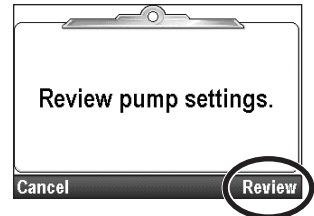
Each protocol includes a therapy, qualifier, and drug.

- **Therapies** are pre-defined and correspond to the delivery modes available on the pump: PCA, continuous, intermittent, step, and taper. No additional therapies can be created.
- **Qualifiers** are specific to a selected therapy and further identify the type of infusion to be delivered. Examples of qualifiers include pediatric, TPN, and epidural.
- **Drugs** are designated by a name/concentration/units/drug ID combination.

New Protocol Selection


Single Protocol Download

1. Make sure the pump is stopped and on the home screen.
2. Connect the pump to the PC using the USB cable. Use the CADD®-Solis Medication Safety Software to load a single protocol on the pump.
3. Acknowledge the low priority alarm.
4. The “Review pump settings” screen appears. Select **Review**.




5. Press **▲** or **▼** to scroll through and review all the settings. If the highlighted setting is correct, press **Accept Value**. A green check mark appears next to the setting.

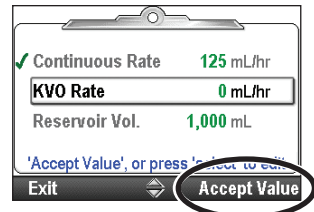
If the highlighted setting is incorrect, press **select** to adjust the setting. Press **▲** or **▼** to modify the setting and then press **Save**.



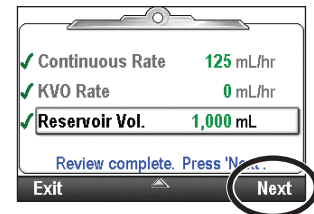6. When your review is complete, select **Next**.



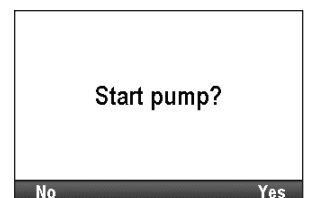7. If a cassette is attached, latched, and locked (if required), the “Start pump?” screen appears. Press **Yes** to start the pump immediately or **No** to start the pump at a later time.



Protocol Selection from a Library

Start New Patient













Each time a new patient is started, it is recorded in the event log. All other reports are cleared. Protocol libraries are created with the CADD[®]-Solis Medication Safety Software and are loaded onto the pump by the CADD[®]-Solis system administrator. You may select a standard therapy protocol from a library created by the CADD[®]-Solis administrator, or manually program the pump.

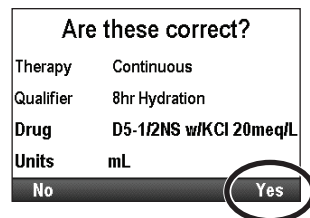
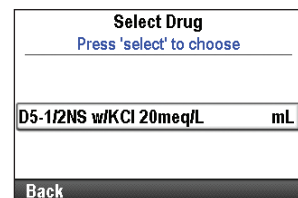
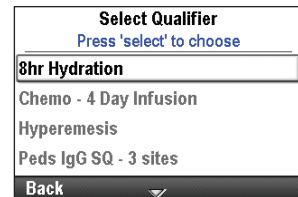
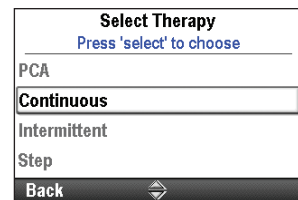
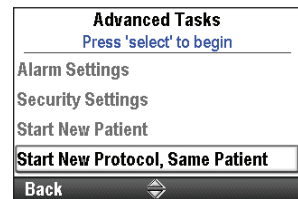
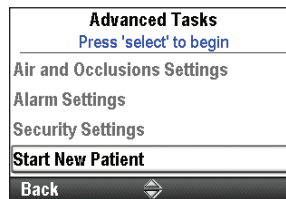
Start New Protocol Same Patient

The process for starting a new protocol for the same patient is the same as starting a new patient. However, the event log does not insert a new patient marker. All other reports are cleared except the delivery log.


NOTE: When starting a new protocol, attach a *new* reservoir with the proper medication.

To start a new patient and run an existing protocol from a library:

1. Make sure the pump is stopped.
2. In the Advanced Tasks menu (see page 99), press  or  to highlight **Start New Patient** or **Start New Protocol Same Patient** and press .
3. Unlock the keypad.
4. Press  or  to highlight a therapy and press .
5. Press  or  to highlight a qualifier and press .
6. Press  or  to highlight a drug and press .
7. A screen appears with a summary of your selections. Review them to confirm that you have entered the correct therapy, qualifier, and drug. Select **Yes**.




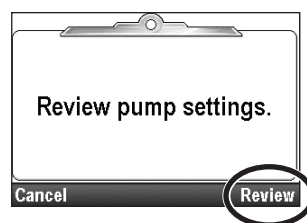8. The “Review pump settings” screen appears. Select **Review**.




9. Press **▲** or **▼** to scroll through and review all the settings. If the highlighted setting is correct, press **Accept Value**. A green check mark appears next to the setting.

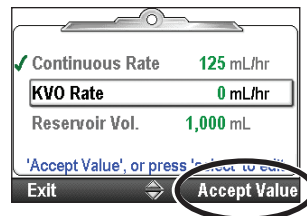If the highlighted setting is incorrect, press **select** to adjust the setting. Press **▲** or **▼** to modify the setting and then press **Save**.




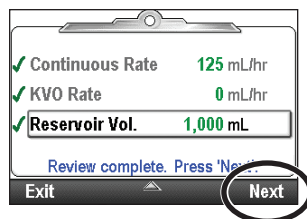10. When the review is complete, select **Next**.



11. When a cassette is attached, latched, and locked (if required by PCA therapies), the "Prime tubing?" screen *may* appear. For more information on priming, see *Prime Tubing* on page 83.

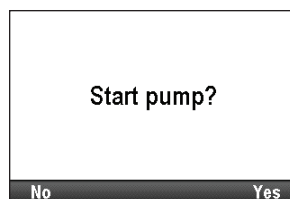The “Start pump?” screen appears. Select **Yes** to start the pump immediately or **No** to start the pump at a later time.



Manually Programmed Protocols

NOTE: Manual programming is for unusual circumstances when an order is received which is different than any of your standard protocols. Programming manually does not use all of the medication safety features built into the CADD®-Solis pump. If you find that you frequently program manually, consider contacting your CADD®-Solis system administrator to discuss if additional protocols should be created.

Manually programming a protocol or making changes to an individual protocol for a single use may occur when the doctor’s orders do not match any of the protocols in the library, or when a protocol library is not available.

Unlike protocols that are created and downloaded into the pump by the CADD®-Solis Medication Safety Software, manual programming allows you to choose units and concentration (in PCA mode), and does not contain any programming limits. Many of the mode-independent parameters may remain the same as the parameters from the previously used protocol. For example, if the previously used protocol had the pump stopped alarm type set as informational, the manual program also has the pump stopped alarm type set as informational. However, the actual programming limits of the previous protocol are cleared. The hard and soft limits are all set to the factory default, which means the delivery ranges are not limited. See the following chart for details.


NOTE: If the pump was set to the factory default, or if you are using the pump for the first time, the previously used protocol is the factory default settings. Refer to the tables beginning on page 140.

The following table lists the initial pump settings when programming manually. *Attributes not listed in the table remain the same as the previously used protocol.*


PCA Mode, Specific Attributes	Value Defaults to:
Clinician bolus amount	0 mL
Concentration	Selected by user
Continuous rate	0 mL/hr (or mg or mcg equivalent)
Delivery limit method	Not in use
KVO rate	0 mL/hr (or mg or mcg equivalent)
Maximum delivery rate	125 mL/hr
PCA dose	0 mL (or mg or mcg equivalent)
PCA dose lockout	1 hr
Programming units	Selected by user
Continuous Mode, Specific Attributes	Value Defaults to:
Continuous rate	0.1 mL/hr
KVO rate	0 mL/hr
Intermittent Mode, Specific Attributes	Value Defaults to:
Dose cycle	24 hr
Dose duration	2 hr
Dose volume	0.1 mL/hr
KVO rate	0 mL/hr
Step Mode, Specific Attributes	Value Defaults to:
Infusion alerts	Off
Initial rate	0.4 mL/hr
KVO rate	0 mL/hr
Plateau rate	0.4 mL/hr
Step duration	30 min
Step rate increment	0.4 mL/hr
Total infusion volume	1 mL
Taper Mode, Specific Attributes	Value Defaults to:
Infusion duration	24 hr
Infusion volume	24 mL
KVO rate	1 mL/hr
Plateau rate upper limit	500 mL/hr
Taper down	0 hr
Taper up	0 hr
Mode-Independent Attributes	Value Defaults to:
Air detector on/off	On
Delayed start requires security	On
Priming requires security	On
Protocol identifiers	Qualifier: [Program Manually] Drug: Non-Library Protocol
Reservoir volume	0 mL
Reservoir volume low trip point	5 mL
Reservoir volume reset value	100 mL
Start time	Now
Unlock with key	False

To start a new patient and manually program the protocol:

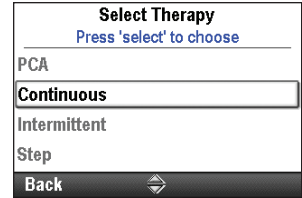
1. Make sure the pump is stopped.
2. In the Advanced Tasks menu (see page 99), press **▲** or **▼** to highlight **Start New Patient** and press **select**.
3. Unlock the keypad.
4. Press **▲** or **▼** to highlight a therapy and press **select**.



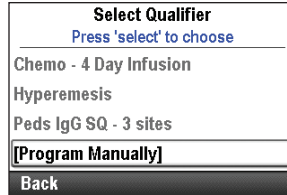5. Press **▲** or **▼** to highlight **[Program Manually]** and press **select**.




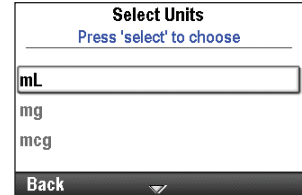If you are programming a PCA therapy, you must also press **▲** or **▼** to highlight the desired units and press **select**. If the units are mg or mcg, press **▲** or **▼** to highlight the desired concentration and select **Save**.



6. A screen appears with a summary of your selections. Confirm that the correct information appears, and select **Yes**.

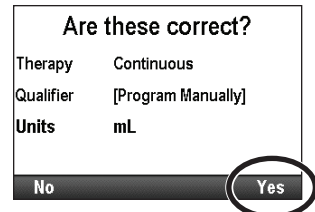NOTE: If the information is incorrect, select the left soft key as many times as required to back out of the screen, allowing you to start over.




7. The “Review pump settings” screen appears. Select **Review**.



8. Press **▲** or **▼** to scroll through all the settings. If the highlighted setting is correct, select **Accept Value**. A green check mark appears next to the setting.


If the setting is incorrect, press **select** to edit. Press **▲** or **▼** to modify the setting and then select **Save** or **Confirm** (see note).

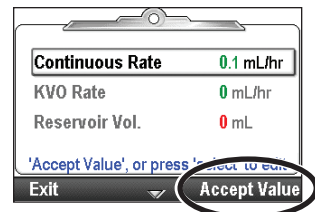NOTE: If you edit any of the parameters, you may exceed the soft limit range and must **Confirm** the soft limit override on each screen. The parameter and value will be displayed in amber.



WARNING: The manual program mode does not contain programming limits. Carefully review each parameter to ensure it accurately matches the prescription. Failure to set the manual program to the correct values could result in serious patient injury or death.

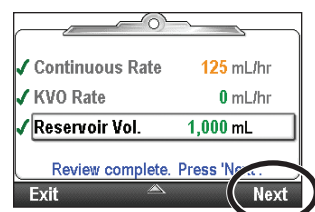9. When the review is complete, select **Next**.



10. Follow the instructions on the pump to attach the cassette, prime, and start the pump.

Delayed Start/Next Dose Start Time

Delayed Start






NOTE: This task is available for all therapies *except* intermittent. If an intermittent therapy is programmed, this option is not visible under the Tasks menu.

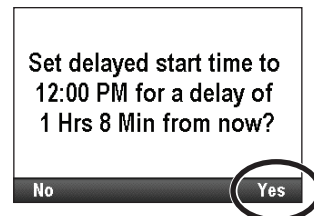
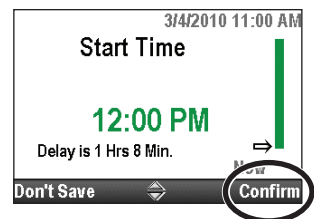
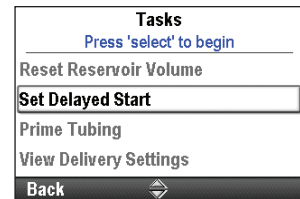
The Set Delayed Start task allows you to delay the start of an infusion by up to 96 hours by selecting the date and time at which the infusion should begin. The pump must be running on the selected date and time in order for delivery to begin.

If Set Delayed Start is programmed, the keep vein open (KVO) rate is active until the infusion begins.

NOTE: In step and taper therapies, you may not set a delayed start after an infusion has been started.

To Set Delayed Start:

1. From the Tasks menu, press  or  until **Set Delayed Start** is highlighted, and then press .
2. Unlock the keypad, if required.
3. Press  or  until the desired start time appears on the screen, and select **Confirm**.
4. Confirm the new start time by selecting **Yes**.



Next Dose Start Time

NOTE: This task is available for intermittent therapies only. If another therapy type is programmed, this option is not visible under the Tasks menu.

A Next Dose Start Time may be programmed to delay the start of the next dose of an infusion. However if a dose is in progress, this cancels the remainder of the current dose. To avoid interrupting the dose in progress, adjust the next dose start time when the pump is in KVO.

NOTE: During dosing, the home screen shows “Dosing”. If a dose is stopped in progress, the home screen shows “Interrupted.” If the infusion was in KVO, “Between Doses” appears on the home screen.

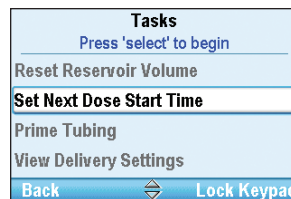
The Next Dose Start Time may also be set to delay delivery of the first dose of the infusion by up to 96 hours. The pump must be running on the selected date and time in order for delivery to begin.



To set the next dose start time:


1. Stop the pump if it is running. If a dose was delivering when the pump was stopped, the “Interrupted” screen appears.

NOTE: If a dose is in progress, this cancels the remainder of the dose.


2. From the Tasks menu, press  or  until **Set Next Dose Start Time** is highlighted, and then press .



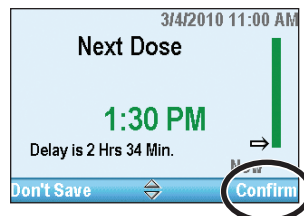3. Unlock the keypad, if required.
4. Press  or  until you reach the time that you want the next dose start time to begin. The length of the delay is listed below the time. Select **Confirm**.



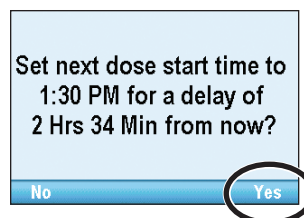5. A screen appears asking you to confirm the new start time and length of the delay. Select **Yes**.



Delivery Hard and Soft Limits














Pump protocol libraries consisting of therapies, qualifiers, and drugs are created and loaded on the pump using the CADD®-Solis Medication Safety Software. Many of the delivery settings have hard and soft limits assigned to them by the CADD®-Solis system administrator. In the unlikely event that any of these values need to be changed manually on the pump, it can be done under the Delivery Hard and Soft Limits setting by a user with the administrator code.

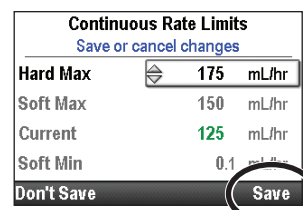
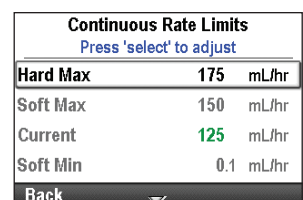
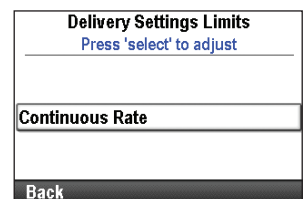
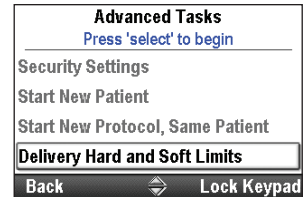
The delivery hard and soft limits define the ranges that are available when programming delivery settings in the View Delivery Settings menu. Not all of the settings have hard and soft maximum and minimum values. However, the settings still effect the limits that appear in the View Delivery Settings menu. For example, in PCA mode, the Delivery Limit Period, Delivery Limit Amount, and Max Doses/Hr may not appear, based on the value of the Delivery Limit Method.

Default delivery hard and soft limits are listed in the table *Manual Mode Initial Delivery Hard and Soft Limits* on page 142. The settings are defined below:

- **Hard Maximum:** The maximum value for the parameter being edited. It is not possible to program the setting above the hard maximum value.
- **Soft Maximum:** The soft maximum value should be the highest common value for the setting you are programming. The setting may be programmed above this value, although extra confirmation is required.
- **Current Value:** Displays the current value for the setting in the therapy. Changing this setting changes the programmed value in the current therapy.
- **Soft Minimum:** The soft minimum value should be the lowest common value for the setting you are programming. The setting may be programmed below this value, although extra confirmation is required.
- **Hard Minimum:** The minimum value for the parameter being edited. It is not possible to program the setting below the hard minimum value.

To set the delivery hard and soft limits:

1. Make sure the pump is stopped.
2. In the Advanced Tasks menu (see page 99), press  or  to highlight **Delivery Hard and Soft Limits** and press .
3. Press  or  to scroll to and highlight the setting you want to adjust and press . In this example, only 1 setting is adjustable.
4. Unlock the keypad.
5. Press  or  to scroll to and highlight the setting you want to adjust and press .
6. Press  or  to increase or decrease the limit for the setting and select **Save**.
7. Continue to scroll through the settings until you have set limits for all the values that you want to set. Select **Back** to go back to the Delivery Settings Limits menu. Press  or  to scroll to any other setting you want to adjust, or press **Back** to finish.



Delivery Settings

The delivery settings are patient-specific parameters of a therapy that are directly related to the drug being infused and can be edited within limits established in the protocol. The CADD[®]-Solis system administrator establishes the initial values of the parameters, any applicable programming units and drug concentration, other parameter limits, and which parameters can be viewed and/or edited. A security code (or the pump key, if enabled) are required to edit the parameters.

The delivery settings are programmed by the CADD[®]-Solis system administrator through the CADD[®]-Solis Medication Safety Software. They may also be manually programmed through the Advanced Tasks menu when setting up a new protocol or new patient.

WARNING: The manual program mode does not contain programming limits. Carefully review each parameter to ensure it accurately matches the prescription. Failure to set the correct values could result in serious patient injury or death.

When programming the initial values for a protocol, the CADD[®]-Solis system administrator also sets the hard and soft limits, which may allow you to modify the parameters as necessary. Any programmed value within and including the soft limits is displayed in green, and is the range most commonly used for the protocol. Any programmed value outside the soft limits but within the hard limits is displayed in amber. The hard limits are the highest and lowest values that can be programmed for each protocol. Changes to any of these settings have a direct impact on the amount of drug delivered to the patient, and any changes made in the Delivery Settings menu affect only the therapy that is currently programmed. If the prescribed orders do not match the initial values, the values should be edited to match. Editing the parameters above or below the soft limits results in a screen that requires you to confirm the soft limit override.

Delivery settings specific to each therapy are shown in the table below.

Delivery Settings

PCA (See page 37)	Continuous (See page 47)	Intermittent (See page 51)	Step (See page 59)	Taper (See page 69)
Continuous Rate	Continuous Rate	Dose Volume	Infusion Volume	Infusion Volume
PCA Dose		Dose Duration	Initial Rate	Taper Up
PCA Lockout *		Dose Cycle	Rate Increment	Taper Down
Delivery Limit †		Dose Rate	Plateau Rate	Infusion Duration
Max Doses/Hr †		Next Dose	Step Duration	Plateau Rate
			Infusion Duration	
KVO Rate	KVO Rate	KVO Rate	KVO Rate	KVO Rate
Reservoir Vol.	Reservoir Vol.	Reservoir Vol.	Reservoir Vol.	Reservoir Vol.
Start Time ‡	Start Time ‡		Start Time ‡	Start Time ‡









* This setting does not appear if the PCA dose is zero.

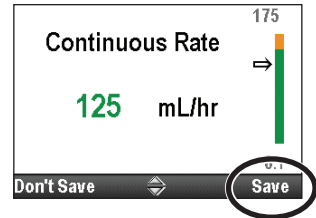
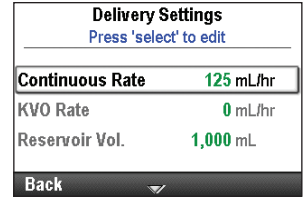
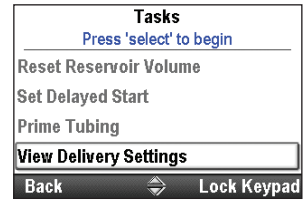
† This setting may not appear, depending on the delivery limit method selected.

‡ This setting appears only if a delayed start is in effect.



Editing Delivery Settings

To view and edit delivery settings:


1. Stop the pump if it is running.
2. In the Tasks menu (see page 88), press  or  until **View Delivery Settings** is highlighted, and then press .
3. Press  or  until the desired setting is highlighted, and then press .
4. If requested, unlock the keypad.
5. Press  or  until the desired value appears on the screen, then select **Save**.



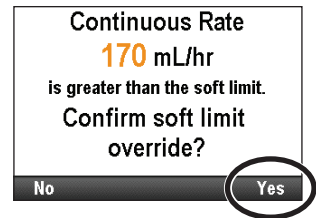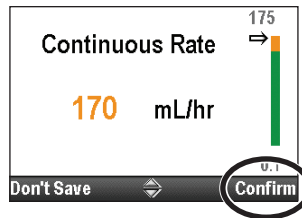NOTE: If the new value is above the maximum soft limit or below the minimum soft limit, it appears in amber on the screen. To select a value above the maximum soft limit or below the minimum soft limit, select **Confirm** and then select **Yes**.



Change any additional settings by scrolling through the remaining delivery settings until all values are correct.

NOTE: Changes to any of the delivery settings must be confirmed and accepted in a pump settings review. The confirmation requires a security code and must take place before the pump will start. Press stop/start  to display the "Review pump settings" screen (see page 36).

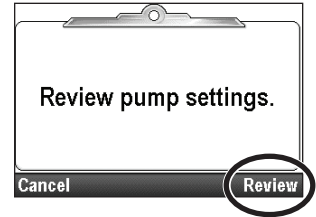
Review Pump Settings

When starting a new patient or a new protocol for the same patient (see page 27), or after adjusting any delivery settings (see page 35), the pump settings must be reviewed and accepted as accurate.


WARNING: Always carefully review the program on the pump after it has been programmed to verify that the pump is programmed correctly. Failure to do so could result in undesired programming of the pump and depending on the drug being administered, death or serious injury to the patient.

To review pump settings:



1. When the “Review pump settings” screen appears, select **Review**.



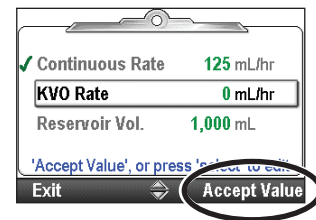2. Press **▲** or **▼** to scroll through and review all the settings. If the highlighted setting is correct, press **Accept Value**. A green check mark appears next to the setting and the next setting is highlighted.




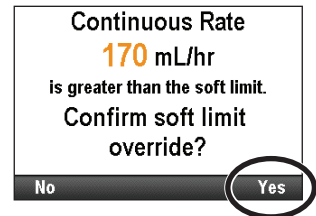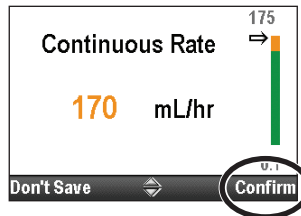If the highlighted setting is incorrect, press **select** to edit the setting. Press **▲** or **▼** to modify the setting and then press **Save** or **Confirm**.



NOTE: If the new value is above the maximum soft limit or below the minimum soft limit, it appears in amber on the screen. To select a value above the maximum soft limit or below the minimum soft limit, select **Confirm**, and then select **Yes** to confirm the soft limit override.

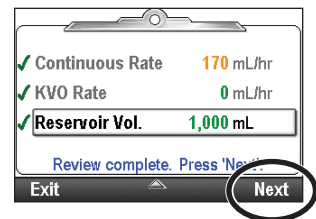3. When the review is complete, select **Next**.




PCA Delivery Mode

PCA (patient-controlled analgesia) delivery is used for therapies that require a continuous rate of infusion, patient-controlled demand doses, and/or a clinician-activated bolus. The delivery methods available are:

- Continuous rate
- PCA dose, a demand dose activated by the patient
- Clinician bolus, a dose activated by the clinician



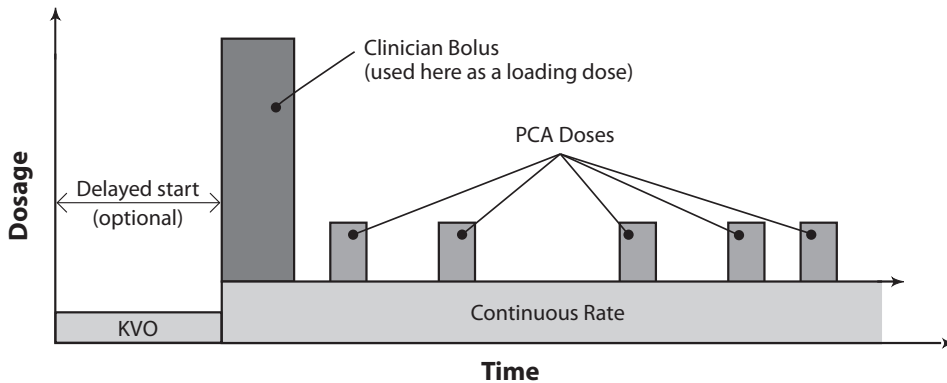
Each of the methods may be programmed individually or in combination with each other. The graph below illustrates the combined delivery methods. Ranges and programming increments are listed in *Specifications (Nominal)* on page 125.



Home Screen

The delivery mode and status of the infusion are indicated on the home screen with messages appearing in text, a graphic in a shape unique to the therapy, and with the colors green, red, and gray. The colors indicate that the pump is running (green), or stopped (red). If a new patient or new infusion is programmed but *the pump has not been started*, or if a delayed start time was programmed and the start time passed while the pump was stopped and not yet started, the graphic is entirely gray. The location of the color on the graphic also indicates the status of the infusion. The examples below are specific to the PCA delivery mode.

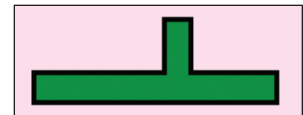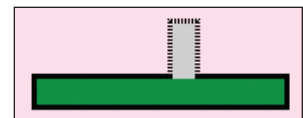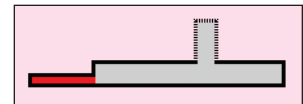- The pump is stopped with a programmed delayed start.
- The pump is running and delivering at a continuous rate.
- The pump is running and delivering a PCA bolus.



PCA graphics do **not** depict the following:

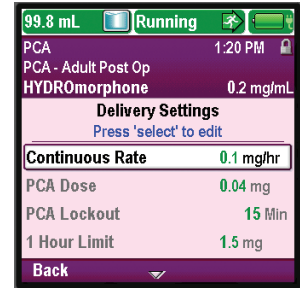
- The amount of the PCA dose compared to the continuous rate.
- The number of PCA doses attempted or delivered.

Programming Screens


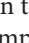
Continuous Rate

This is the desired continuous rate of medication delivery. If the prescription does not call for a continuous rate during a PCA therapy, enter 0.

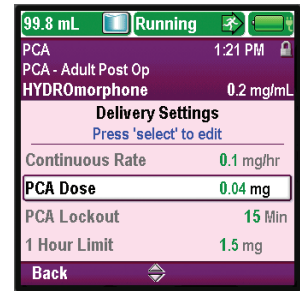
The continuous rate value may be changed by scrolling up or down between the hard maximum and hard minimum settings.



PCA Dose


This is the amount of drug delivered when the patient requests a PCA dose by either pressing the PCA dose  key or using the remote dose cord button. The PCA dose  key on the pump is inactive when the remote dose cord is connected to the pump. If the prescription does not call for a PCA dose, enter 0.

WARNING: The remote dose cord is for patient use only. Operation by anyone other than the patient may cause overdosing of medication that could result in serious injury or death of the patient.



PCA Lockout

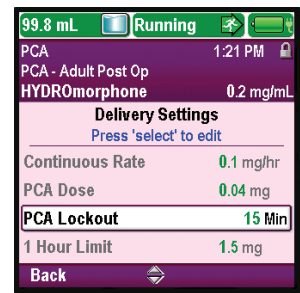
This is the minimum time that must elapse between the time one PCA dose starts and the time that the next PCA dose is available.



Delivery Limit


This is used to restrict the amount of drug delivered to the patient within a specified time frame. The programmable range is 1 hour to 12 hours. The limit includes the continuous rate and PCA doses, but does not include clinician boluses.

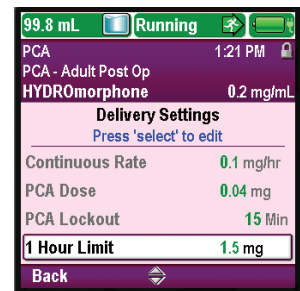The frequency of available doses are limited by the programmed PCA lockout time. If the delivery limit is reached while a PCA dose is in progress, the PCA dose will not be completed. If the delivery limit is not visible in the menu, your CADD[®]-Solis system administrator has elected not to require it.



Max Doses/Hr

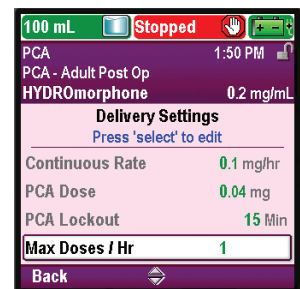This is the maximum number of PCA doses allowed in a 1 hour period. The frequency of available doses are limited by the programmed PCA lockout time. If the PCA Lockout is 1 hour or greater, this setting will not appear in the menu. If the setting is not visible in the menu, your CADD[®]-Solis system administrator has elected not to require it.



KVO Rate

The KVO or “keep vein open” rate is optional. It allows the delivery of a minimal amount of drug to help maintain catheter patency. If a continuous rate is programmed, the KVO rate automatically displays 0.1 mL/mg (or the mg or mcg equivalent).

The KVO is active under 2 conditions: If a delayed start is programmed, the KVO rate is active until the infusion begins. If a delivery limit is programmed, the KVO rate becomes active when the delivery limit is reached.

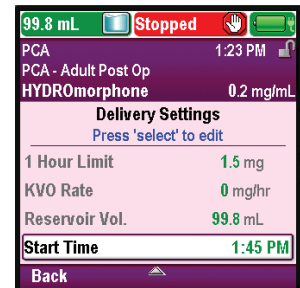
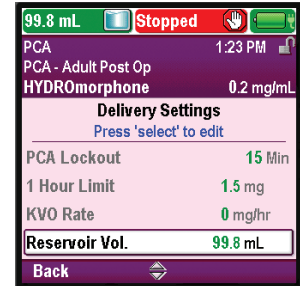
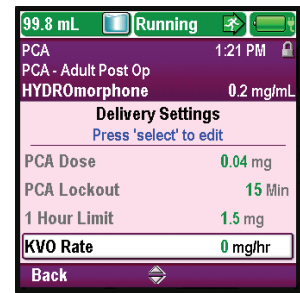
Reservoir Vol.

This is the volume of fluid contained in the reservoir. The reservoir volume value decreases as the pump is primed or delivers fluid.

The CADD[®]-Solis system administrator may choose a standard reservoir volume which allows the reservoir volume to be quickly reset to that value. The clinician can edit the reservoir volume to another amount. This change is for the current protocol only and remains in effect until the Start New Patient or Same Patient New Protocol task is selected.

Start Time

This is the time that the next infusion will begin. It is displayed only if a delayed start is programmed.




Programming Example

WARNING: Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is used to deliver critical or life-sustaining medication, the interruption in the delivery of medication could result in patient injury or death.


Example: Medication is provided in a 100 mL CADD™ medication cassette reservoir at a concentration of 1 mg/mL. The patient must receive medication at a continuous rate of 5 mg/hr. PCA doses of 2 mg are allowed with a 15 minute lockout time between doses and a maximum of 2 doses per hour.

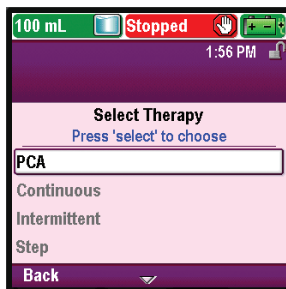1. Select the Therapy




- Stop the pump if it is running.
- Select **Tasks**, then **View Advanced Tasks**, and then **Start New Patient**.
- Unlock the keypad.
- Press **▲** or **▼** to highlight the **PCA** therapy and press **select**. See *Start New Patient* on page 27 if you need more information on selecting the therapy.

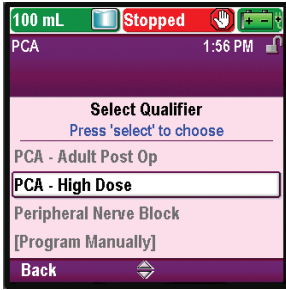2. Select the Qualifier




- Press **▲** or **▼** to highlight the desired qualifier and press **select**.

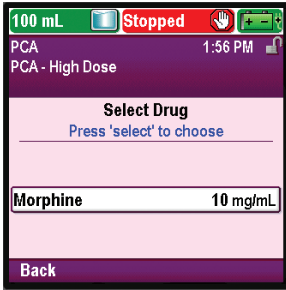3. Select the Drug




- Press **▲** or **▼** to highlight the desired drug and concentration, and press **select**.

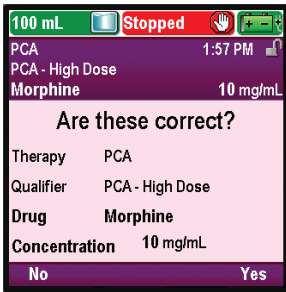4. Confirm and Review the Settings




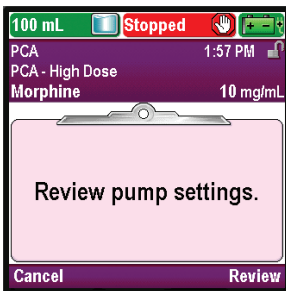- Review the information on the screen and select **Yes** if they are correct, or **No** to edit the settings.




- The pump displays the **Review pump settings** screen. Select **Review** to continue.

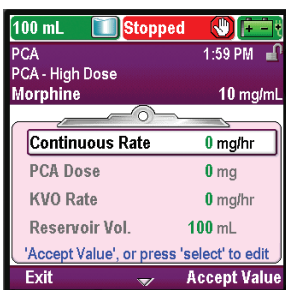5. Enter the Continuous Rate




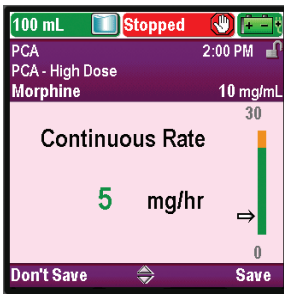- Press **▲** or **▼** to highlight **Continuous Rate**. If the value shown is set at the desired rate, select **Accept Value** and go to step 6. If the value shown does not show the desired rate, press **select**.




- Press **▲** or **▼** until the value reads **5 mg/hr** and select **Save**.
- The pump returns to the Review screen.

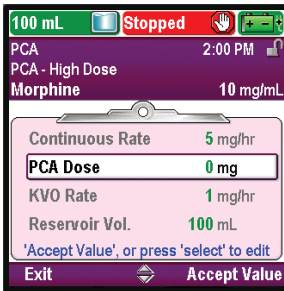6. Enter the PCA Dose




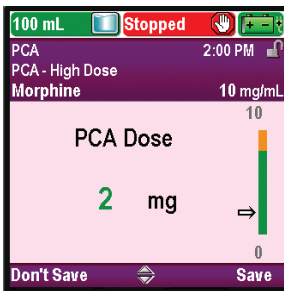- Press **▲** or **▼** to highlight **PCA Dose**. If the value shown is set at the desired dose, select **Accept Value** and go to step 7. If the value shown does not show the desired dose, press **select**.




- Press **▲** or **▼** until the value reads **2 mg** and select **Save**.
- The pump then returns to the Review screen.

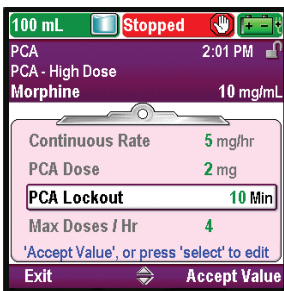7. Enter the PCA Lockout




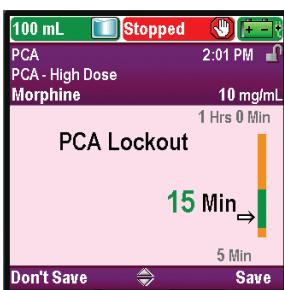- Press **▲** or **▼** to highlight **PCA Lockout**. If the time shown is set at the desired value, select **Accept Value** and go to step 8. If the value shown does not show the desired time, press **select**.




- Press **▲** or **▼** until the value reads **15 Min** and select **Save**.
- The pump then returns to the Review screen.

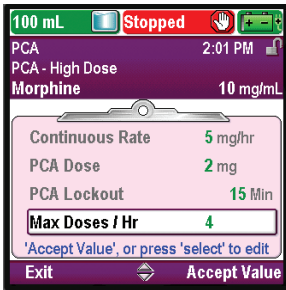8. Enter the Maximum Doses per Hour




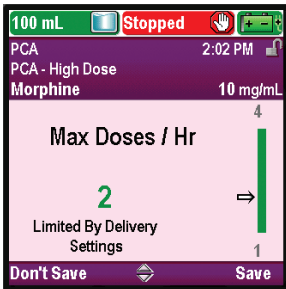- Press **▲** or **▼** to highlight **Max Doses/Hr**. If the value shown is correct, select **Accept Value** and go to step 9. If the value is not correct, press **select**.




- Press **▲** or **▼** until the value reads **2** and select **Save**.
- The pump then returns to the Review screen.

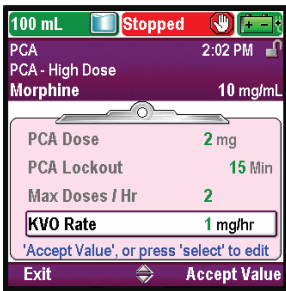9. Review the KVO Rate




- In this example, the KVO Rate is a calculated value and cannot be edited.

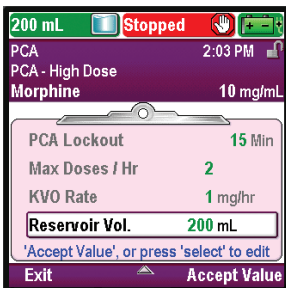10. Enter the Reservoir Volume




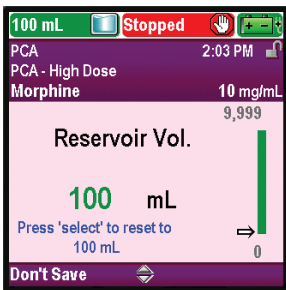- Press **▲** or **▼** to highlight **Reservoir Vol.**. If the value shown is correct, select **Accept Value** and go to step 11. If the value is not correct, press **select**.



- Press **▲** or **▼** until the value reads **100 mL** and select **Save**.
- The pump then returns to the Review screen.

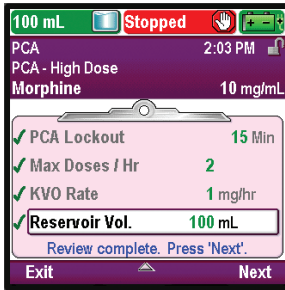11. Accept Values



- If any values have not already been accepted, press or to highlight each setting and select **Accept Value** if the setting is correct. A green check mark appears next to each accepted value.
- After your review is complete, select **Next**.

12. Prepare the Pump for the Patient

- Follow the instructions for attaching a cassette, priming, and attaching the pump to the patient in *Operating the Pump* on page 78).

Give a Clinician Bolus

WARNING: Exercise care when using the clinician bolus function. Because there are no limits to the frequency of delivering a bolus, and because the bolus amount may be set as high as 50 mL (or the mg or mcg equivalent), you should not permit the patient or unauthorized clinician to become familiar with the procedure for giving a clinician bolus. Improper programming could result in serious patient injury or death.

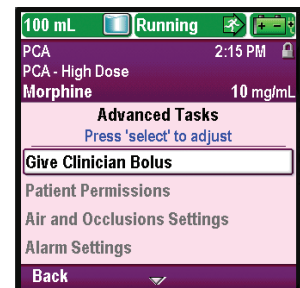
A clinician bolus may be delivered only when the pump is running. It allows delivery of a specified amount of drug, for example, as a loading dose. A clinician bolus cannot be started while a PCA dose is in progress. The amount delivered decreases the reservoir volume and increases the given amount, but does not add to the dose counters or to the delivery limit. A clinician bolus may be stopped in progress.

NOTE: If a clinician bolus is manually stopped by a clinician, or automatically stopped by an alarm, power failure, or other condition that stops delivery, the pump remembers at what point the bolus was when it stopped. The next time a clinician bolus is given, you will have the option to restart the clinician bolus where it left off, or to start with a new clinician bolus.

NOTE: The maximum clinician bolus may be limited by the settings in the protocol, which is determined by the CADD®-Solis system administrator.



To start a clinician bolus:

1. Make sure the pump is running.
2. From the Advanced Tasks menu (see page 99), press or to highlight **Give Clinician Bolus** and press .



3. Unlock the keypad.

WARNING: To prevent the patient from accessing the clinician bolus function, do not let the patient know the clinician or administrator security codes. Improper programming could result in serious patient injury or death.

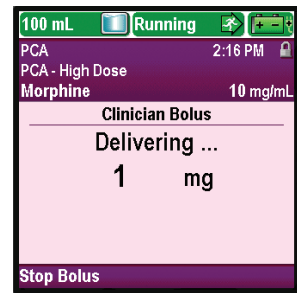
- Press  or  to change the value of the clinician bolus. Select **Deliver**.




- If you enter a value outside the soft limit, a screen appears asking you to confirm the soft limit override. Select **Confirm** to deliver the bolus.

WARNING: Never leave the pump unattended while on the clinician bolus edit screen. You must deliver the programmed value, or cancel and leave the screen. Failure to do so could result in serious patient injury or death.

- The screen shows the decreasing amount as the bolus is delivered. You may stop the bolus at any time by selecting **Stop Bolus**.



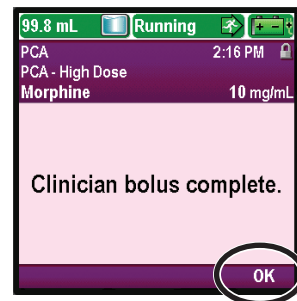
- When the bolus has been delivered, the screen says, "Clinician bolus complete." Select **OK**.



Start a PCA Dose

If a PCA dose is programmed, the patient may start a PCA dose while the pump is running. The amount delivered is added to the amount provided by the continuous rate. Each time the patient requests a PCA dose, the pump automatically adds it to the Given and PCA Dose Counters report. If no PCA dose is programmed, the pump displays the message “PCA dose not available because no dose programmed”.

If the patient attempts to deliver a PCA dose during the lockout time, “PCA dose not available. Currently locked out.” appears on the display and the pump will not deliver the dose. The lockout time is determined by the PCA Lockout or the Max Doses/Hr (if in use), **whichever limits the dose frequency more**. The attempt is added to the Attempted counter on the Given and PCA Dose Counters and the PCA Dose Graph reports.

The remote dose cord can be attached to the pump and provided to the patient as an alternative to pushing the PCA dose  key when requesting a PCA dose. The LED on the remote dose cord indicates PCA dose status:



Off: A PCA dose is not available.

Flashing: A PCA dose is available.


On: A PCA dose has been requested and delivery has started.

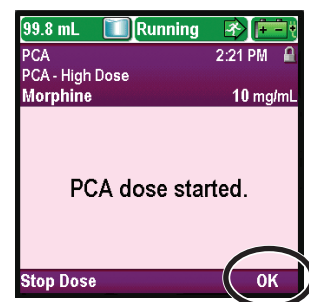
WARNING: The remote dose cord is for *patient use only*. Operation by anyone other than the patient may cause overdosing of medication that could result in serious injury or death of the patient.

NOTES:

- A PCA dose cannot be started while another PCA dose or a clinician bolus is in progress.
- The PCA dose  key on the pump is inactive when the remote dose cord is connected to the pump
- Even if the display has automatically blanked, pressing the PCA dose  key turns the display back on and delivers a PCA dose, if available.

To start a PCA dose:


Press the PCA dose  key on the pump or the remote dose cord. Select **OK**. (If OK is not selected, the screen returns to the previous screen after 5 seconds.)



Stopping a Clinician Bolus or PCA Dose

To stop a clinician bolus:

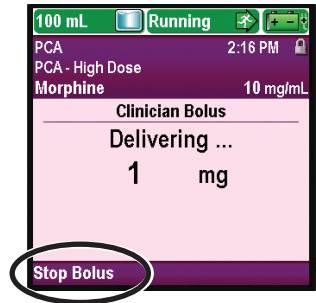
While the clinician bolus is delivering, the Clinician Bolus screen appears on the pump and shows the decreasing amount of the bolus remaining to be delivered. You may stop the clinician bolus from this screen by selecting **Stop Bolus**.




To stop a PCA dose:

1. From the "PCA dose started" screen you may stop the PCA dose while it is delivering by selecting **Stop Dose**.

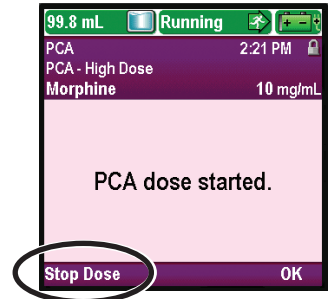
OR



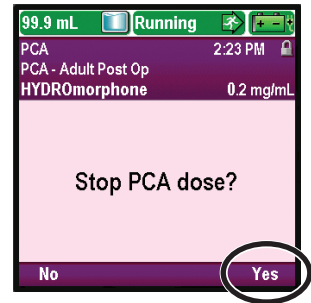2. If the pump is on the home screen, press the stop/start  button.




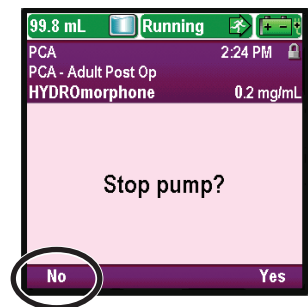You will be asked if you want to stop the pump as well as the PCA dose. Select **No** to stop only the PCA dose.



Continuous Delivery Mode

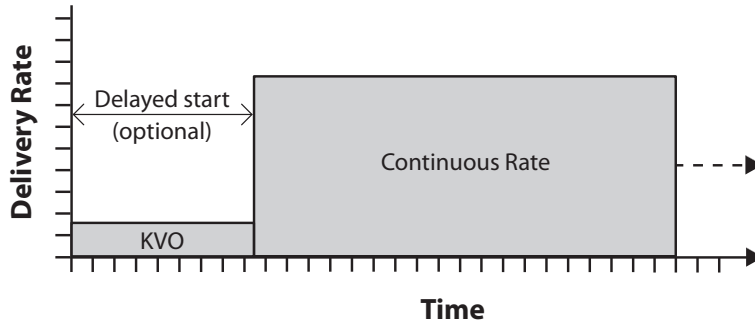
Continuous delivery allows the infusion of drug at a constant, programmed rate. The following graph illustrates an example continuous delivery.



Home Screen

The delivery mode and status of the infusion are indicated on the home screen with messages appearing in text, a graphic in a shape unique to the therapy, and with the colors green, red, and gray. The colors indicate that the pump is running (green), or stopped (red). If a new patient or new infusion is programmed but *the pump has not been started*, or if a delayed start time was programmed and the start time passed while the pump was stopped and not yet started, the graphic is entirely gray. The location of the color on the graphic also indicates the status of the infusion. The examples below are specific to the continuous delivery mode.

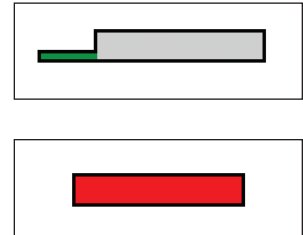- The pump is running with a programmed delayed start.
- The pump is stopped. It will deliver at a continuous rate once it is started.




Continuous delivery mode graphics do *not* depict the relative rate of delivery.

Programming Screens


Continuous Rate

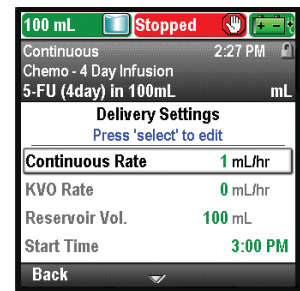This is the desired continuous rate of medication delivery. Rates above 250 mL/hr require a CADD® high volume administration set. The continuous rate value may be changed by scrolling up or down between the hard maximum and hard minimum settings.



KVO Rate

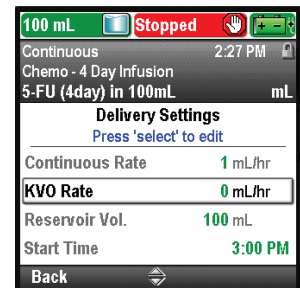The KVO or “keep vein open” rate is optional. It allows the delivery of a minimal amount of drug to help maintain catheter patency. The KVO rate is active only if a delayed start is programmed. The pump infuses at the KVO rate until the programmed infusion begins.





Reservoir Vol.

This is the volume of fluid contained in the reservoir. The reservoir volume value decreases as the pump is primed or delivers fluid.

The CADD®-Solis system administrator may choose a standard reservoir volume which allows the reservoir volume to be quickly reset to that value. The clinician can edit the reservoir volume to another amount. This change is for the current protocol only and remains in effect until the Start New Patient or Same Patient New Protocol task is selected.

Start Time

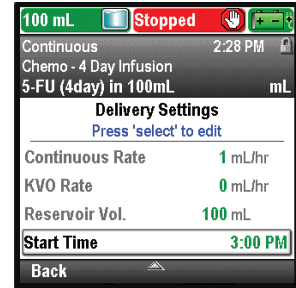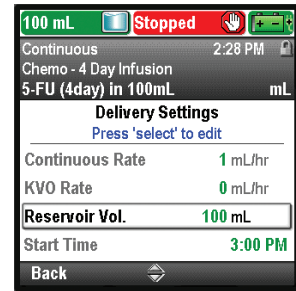This is the time that the next infusion will begin. It is displayed only if a delayed start is programmed.




Programming Example

WARNING: Ensure that the ± 6% system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is used to deliver critical or life-sustaining medication, the interruption in the delivery of medication could result in patient injury or death.


Example: Medication is provided in a 1000 mL standard IV bag. The patient must receive medication at a continuous rate of 125 mL/hr. The KVO rate should be set to 0 mL/hr.

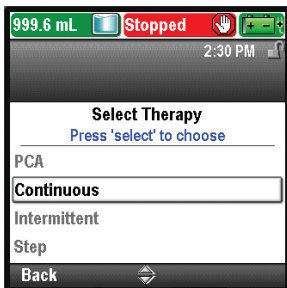1. Select the Therapy




- Stop the pump if it is running.
- Select **Tasks**, then **View Advanced Tasks**, and then **Start New Patient**.
- Unlock the keypad.
- Press **▲** or **▼** to highlight the **Continuous** therapy and press **select**. See *Start New Patient* on page 27 if you need more information on selecting the therapy.

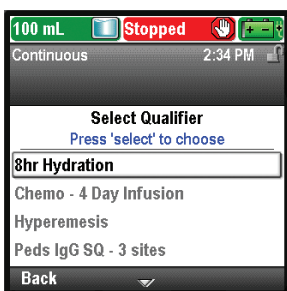2. Select the Qualifier




- Press **▲** or **▼** to highlight the desired qualifier and press **select**.

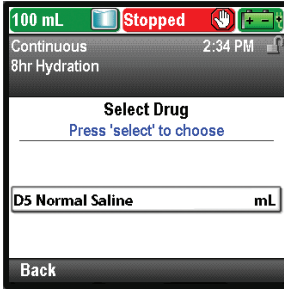3. Select the Drug



- Press **▲** or **▼** to highlight the desired drug, and press **select**.

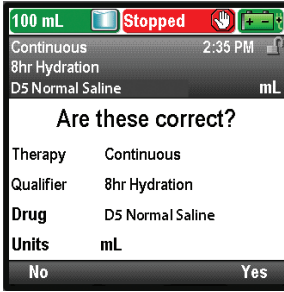4. Confirm and Review the Settings




- Review the information on the screen and select **Yes** if they are correct, or **No** to edit the settings.




- The pump displays the **Review pump settings** screen. Select **Review** to continue.

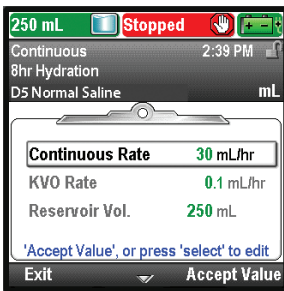5. Enter the Continuous Rate




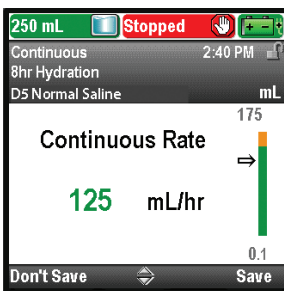- Press **▲** or **▼** to highlight **Continuous Rate**. If the value shown is set at the desired rate, select **Accept Value** and go to step 6. If the value shown does not show the desired rate, press **select**.




- Press **▲** or **▼** until the value reads **125 mL/hr** and select **Save**.
- The pump returns to the Review screen.

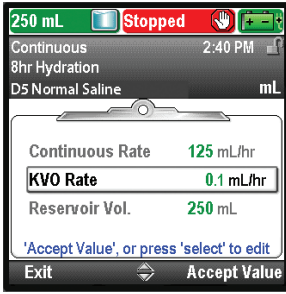6. Enter the KVO rate




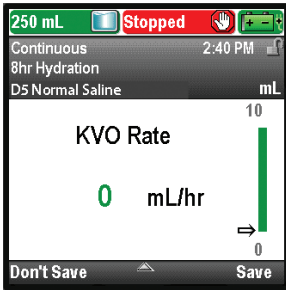- Press **▲** or **▼** to highlight **KVO Rate**. If the value shown is set at the desired rate, select **Accept Value** and go to step 7. If the value shown does not show the desired rate, press **select**.




- Press **▲** or **▼** until the value reads **0 mL/hr** and select **Save**. The pump then returns to the Review screen.

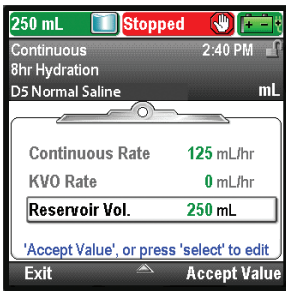7. Enter the Reservoir Volume




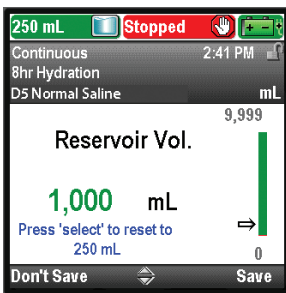- Press **▲** or **▼** to highlight **Reservoir Vol.**. If the value shown is correct, select **Accept Value** and go to step 8. If the value is not correct, press **select**.



- Press **▲** or **▼** until the value reads **1000 mL** and select **Save**.
- The pump then returns to the Review screen.

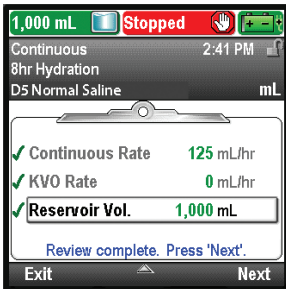8. Accept Values



- If any values have not already been accepted, press **▲** or **▼** to highlight each setting and select **Accept Value** if the setting is correct. A green check mark appears next to each accepted value.
- After your review is complete, select **Next**.


9. Prepare the Pump for the Patient

- Follow the instructions for attaching a cassette, priming, and attaching the pump to the patient in *Operating the Pump* on page 78).

Intermittent Delivery Mode

Intermittent delivery allows the infusion of a specific volume of drug at a regular, programmed interval and over a specified duration. A dose may be repeated in a cycle up to 96 hours. The KVO (keep vein open) feature allows delivery of a minimal amount of drug between doses to maintain catheter patency. You may also delay the start of the first or next dose by using the Next Dose Start Time task.




The following graph illustrates an example intermittent delivery:

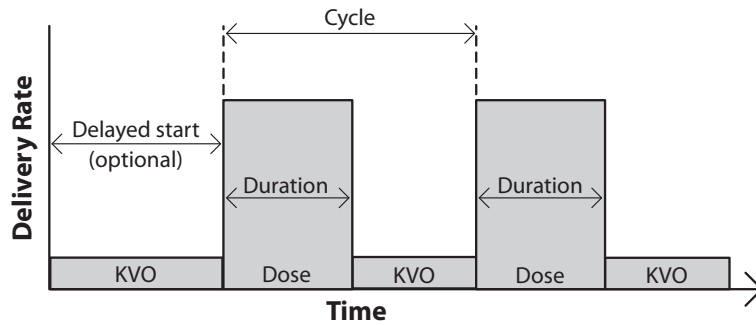

Home Screen

The delivery mode and status of the infusion are indicated on the home screen with messages appearing in text, a graphic in a shape unique to the therapy, and with the colors green, red, and gray. The colors indicate that the pump is running (green), or stopped (red). If a new patient or new infusion is programmed but *the pump has not been started*, or if a next dose start time was programmed and the next dose start time passed while the pump was stopped and not yet started, the graphic is entirely gray. The location of the color on the graphic also indicates the status of the infusion. The examples below are specific to the intermittent delivery mode.

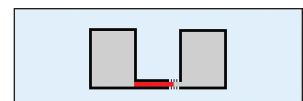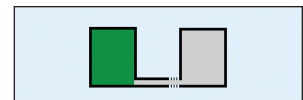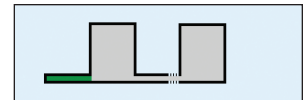- The pump is running with a manually programmed next dose start time.
- The pump is running and delivering a dose.
- The pump is stopped between doses while at the KVO rate.




Intermittent delivery mode graphics do **not** depict the following:

- The relative ratio of dose period to dose cycle.
- The number of doses remaining in the current reservoir volume.

Programming Screens


Dose Volume

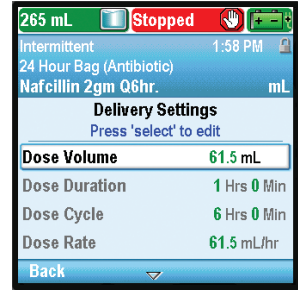This is the volume of the dose. The dose volume cannot be programmed to a value that would cause the rate of delivery for the programmed dose duration to exceed the dose rate maximum of 500 mL/hr. Therefore, programming the dose volume may result in the automatic lengthening of the dose duration to accommodate the maximum dose rate. Additionally, there must be at least 5 minutes between the end of one dose and the start of the next dose (the dose cycle must be 5 minutes longer than the dose duration). Therefore, any automatic lengthening of the dose duration may result in an automatic lengthening of the dose cycle in order to maintain the 5 minute difference. Entering a dose volume automatically resets the dose cycle and sets the next dose start time to begin immediately.




Dose Duration

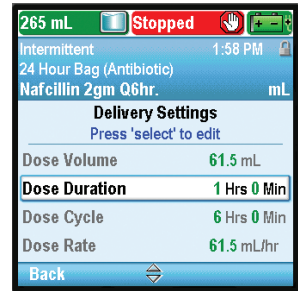This is the length of time required to deliver the dose. The dose duration cannot be programmed to a value that would cause the rate of delivery for the programmed dose volume to exceed the maximum dose rate of 500 mL/hr. Additionally, there must be at least 5 minutes between the end of one dose and the start of the next dose (the dose cycle must be 5 minutes longer than the dose duration). Therefore, programming the dose duration may result in an automatic lengthening of the dose cycle in order to maintain the 5 minute difference. Entering a dose duration automatically resets the dose cycle and sets the next dose start time to begin immediately.




Dose Cycle

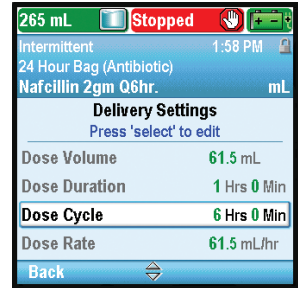This is the time from the start of one dose to the start of the next dose. The programmable values for the dose cycle are dependent on the dose duration. There must be at least 5 minutes between the end of one dose and the start of the next dose (the dose cycle must be 5 minutes longer than the dose duration). Therefore, the minimum programmable cycle is the dose duration plus 5 minutes. The maximum dose cycle is 96 hours. Entering a dose cycle automatically sets the next dose start time to begin immediately.




Dose Rate

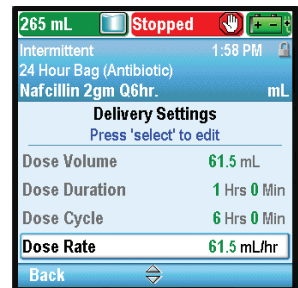This setting is for review only and shows the rate at which the dose is delivered based on the programmed dose volume and dose duration. The maximum allowable rate is 500 mL/hr. Rates above 250 mL/hr require a CADD® high volume administration set.



Next Dose

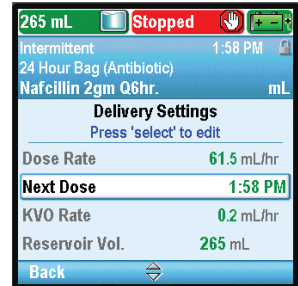
This is the time the next dose is scheduled to begin. The pump must be running on the selected date and time in order for the dose to begin. The next dose time can be programmed at any time to delay the start of the next dose, but if a dose is in progress, this cancels the remainder of the dose. While dosing, the home screen shows "Dosing". If a dose is stopped in progress, the home screen shows "Interrupted".



KVO Rate

The KVO or “keep vein open” rate is optional. It allows the delivery of a minimal amount of drug to help maintain catheter patency.

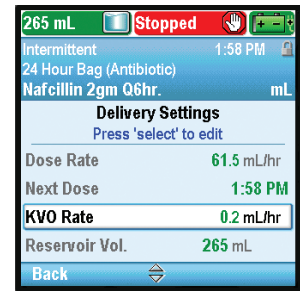
If a next dose start time is programmed, the KVO rate is active until the dose begins. It is also active between doses.



Reservoir Vol.

This is the volume of fluid contained in the reservoir. The reservoir volume value decreases as the pump is primed or delivers fluid.

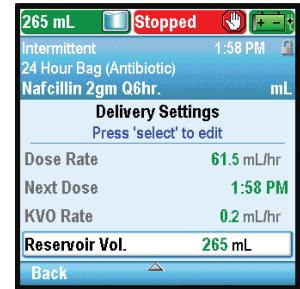The CADD®-Solis system administrator may choose a standard reservoir volume which allows the reservoir volume to be quickly reset to that value. The clinician can edit the reservoir volume to another amount. This change is for the current protocol only and remains in effect until the Start New Patient or Same Patient New Protocol task is selected.




Programming Example

WARNING: Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is used to deliver critical or life-sustaining medication, the interruption in the delivery of medication could result in patient injury or death.


Example: Medication is provided in a 100 mL CADD™ medication cassette reservoir. The patient must receive a 23.5 mL dose over 1 hour. The dose must be given every 6 hours with a 0.2 mL/hr KVO rate between doses. The first dose must begin at 8:00 PM tonight.

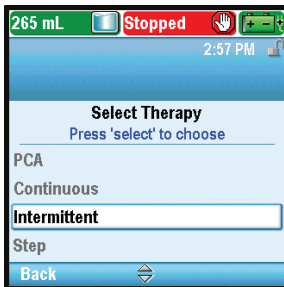1. Select the Therapy




- Stop the pump if it is running.
- Select **Tasks**, then **View Advanced Tasks**, and then **Start New Patient**.
- Unlock the keypad.
- Press **▲** or **▼** to highlight the **Intermittent** therapy and press **select**. See *Start New Patient* on page 27 if you need more information on selecting the therapy.

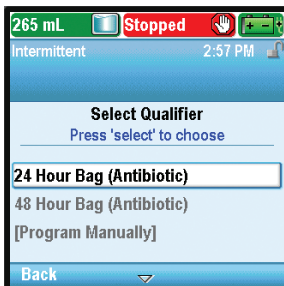2. Select the Qualifier




- Press **▲** or **▼** to highlight the desired qualifier and press **select**.

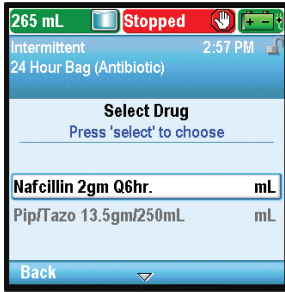3. Select the Drug



- Press **▲** or **▼** to highlight the desired drug, and press **select**.

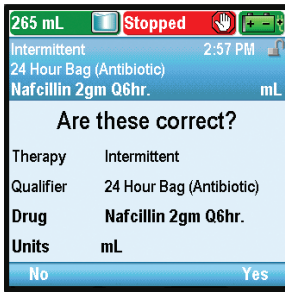4. Confirm and Review the Settings




- Review the information on the screen and select **Yes** if they are correct, or **No** to edit the settings.



- The pump displays the **Review pump settings** screen. Select **Review** to continue.

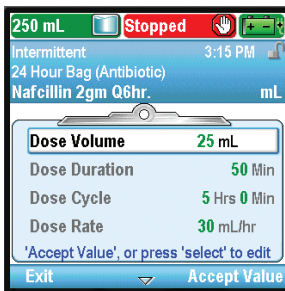5. Enter the Dose Volume




- Press **▲** or **▼** to highlight **Dose Volume**. If the value shown is set at the desired volume, select **Accept Value** and go to step 6. If the value shown does not show the desired volume, press **select**.




- Press **▲** or **▼** until the value reads **23.5 mL** and select **Save**.
- The pump returns to the Review screen.

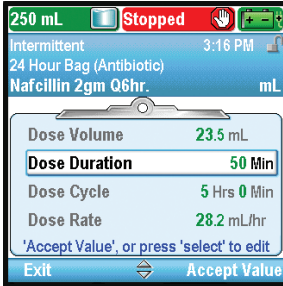6. Enter the Dose Duration




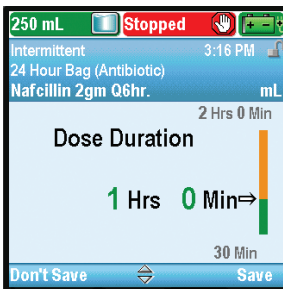- Press **▲** or **▼** to highlight **Dose Duration**. If the value shown is set at the desired duration, select **Accept Value** and go to step 7. If the value shown does not show the desired duration, press **select**.




- Press **▲** or **▼** until the value reads **1 Hrs 0 Min** and select **Save**. The pump then returns to the Review screen.

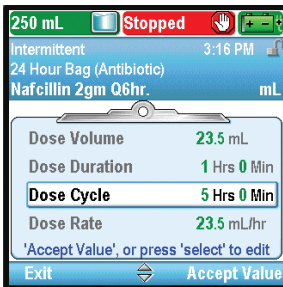7. Enter the Dose Cycle




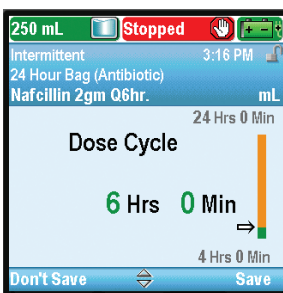- Press **▲** or **▼** to highlight **Dose Cycle**. If the value shown is set at the desired cycle, select **Accept Value** and go to step 8. If the value shown does not show the desired cycle, press **select**.



- Press **▲** or **▼** until the value reads **6 Hrs 0 Min** and select **Save**.
- The pump then returns to the Review screen.

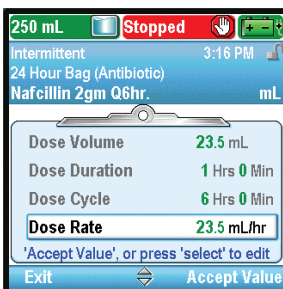8. Enter the Dose Rate




- In this example, Dose Rate is a calculated value and cannot be edited.


9. Enter the Next Dose



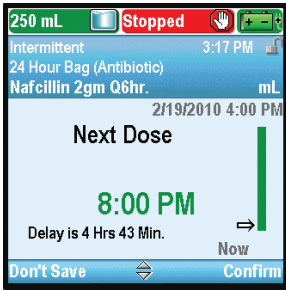- Press **▲** or **▼** to highlight **Next Dose**. If the time shown is correct, select **Accept Value** and go to step 9. If the value is not correct, press **select**.




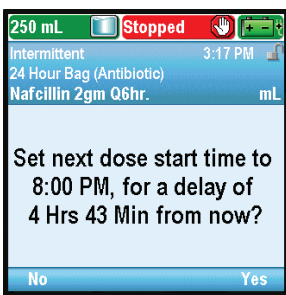- Press **▲** or **▼** until the value reads **8:00 PM** and select **Confirm**.




- Confirm the next dose start time by selecting **Yes**.
- The pump then returns to the Review screen.

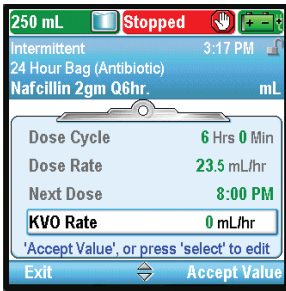10. Enter the KVO Rate



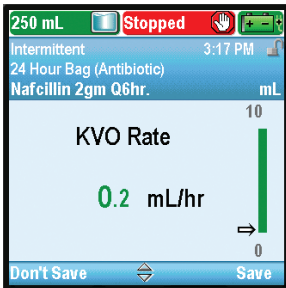- Press **▲** or **▼** to highlight **KVO Rate**. If the value shown is correct, select **Accept Value** and go to step 10. If the value is not correct, press **select**.




- Press **▲** or **▼** until the value reads **0.2 mL/hr** and select **Save**.
- The pump then returns to the Review screen.

11. Enter the Reservoir Volume




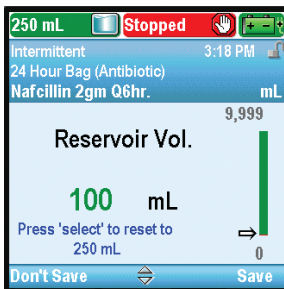- Press **▲** or **▼** to highlight **Reservoir Vol.** If the value shown is correct, select **Accept Value** and go to step 11. If the value is not correct, press **select**.



- Press **▲** or **▼** until the value reads **100 mL** and select **Save**.
- The pump then returns to the Review screen.

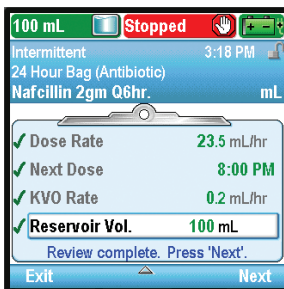12. Accept Values



- If any values have not already been accepted, press **▲** or **▼** to highlight each setting and select **Accept Value** if the setting is correct. A green check mark appears next to each accepted value.
- After your review is complete, select **Next**.

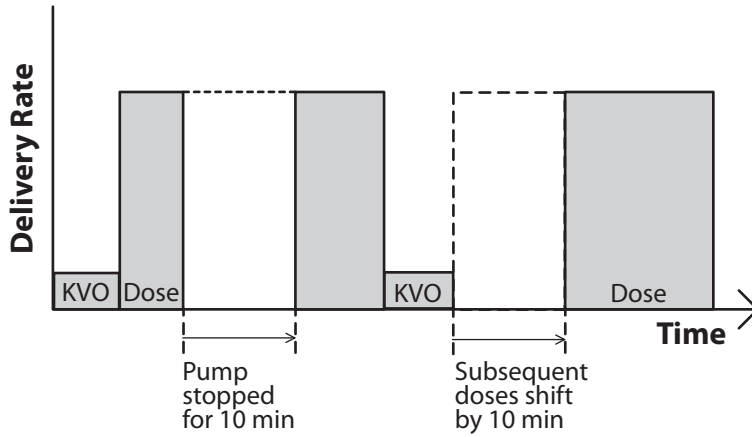
13. Prepare the Pump for the Patient

- Follow the instructions for attaching a cassette, priming, and attaching the pump to the patient in *Operating the Pump* on page 78).

Stopping the Pump During an Infusion


Stopping the pump *between doses* does not affect the start time of subsequent doses.

Stopping the pump *while a dose is in progress* shifts all subsequent doses by the amount of time the pump is stopped (see illustration).



If you want to make up for the time lost in the cycle because the pump was stopped during a dose, *wait until the current dose is completed*. Then, stop the pump and reprogram the Next Dose Start Time. For additional information, see *Next Dose Start Time* on page 32.

Resuming the Dose

To resume delivery of a dose stopped in progress, simply restart the pump by pressing stop/start .

Resetting the Cycle

To reset the cycle, stop the pump between doses, and reprogram the Next Dose Start Time.


When a dose is stopped in progress, the “Interrupted” screen appears. To cancel the current dose and reset the cycle, reprogram the Next Dose Start Time through the Tasks menu. The next cycle will start at the time you select.

See *Next Dose Start Time* on page 32 for more information on setting the Next Dose Start Time.

Step Delivery Mode

Step delivery allows the infusion of a specified volume of drug at an initial rate with linear step increases up to a plateau rate. Multiple steps down may also be applied, reducing the continuous rate as low as the initial rate. An optional KVO (keep vein open) rate may be delivered at the end of the infusion, depending on protocol settings.



The following graph illustrates an example step delivery:

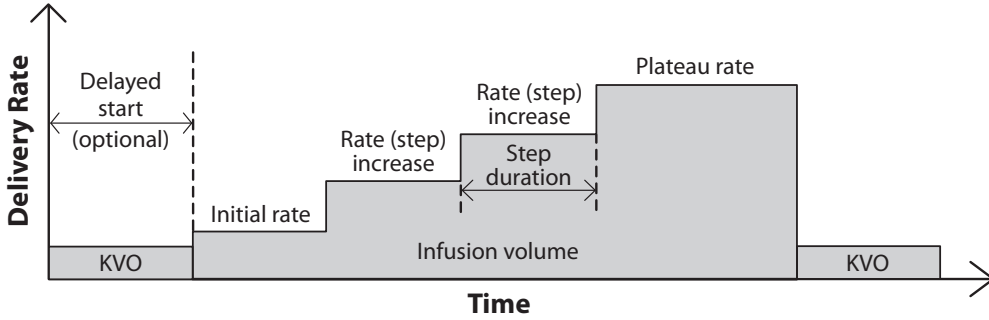

Home Screen

The delivery mode and status of the infusion are indicated on the home screen with messages appearing in text, a graphic in a shape unique to the therapy, and with the colors green, red, and gray. The colors indicate that the pump is running (green), or stopped (red). If a new patient or new infusion is programmed but *the pump has not been started*, or if a delayed start time was programmed and the start time passed while the pump was stopped and not yet started, the graphic is entirely gray. The location of the color on the graphic also indicates the status of the infusion. The examples below are specific to the step delivery mode.

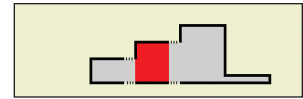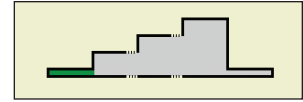- The pump is running with a programmed delayed start.
- The pump is stopped and has not yet reached the plateau.
- The pump is running at the plateau rate.
- The pump is running at continuous hold after a manual step down.
- The pump is running at the KVO rate, nearing the end of the infusion.




Step delivery mode graphics do **not** depict how many steps will be made to reach the plateau. Nor do they depict the relative step rate increases in comparison to the initial and plateau rates.

Programming Screens

Infusion Volume


This is the total volume of fluid to be delivered. The reservoir volume must be large enough to support the entire infusion volume. Therefore, programming the infusion volume may result in an automatic increase in the reservoir volume.

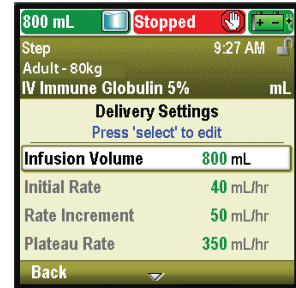Entering a new infusion volume resets the infusion so that delivery starts at the beginning of the infusion duration.




Initial Rate

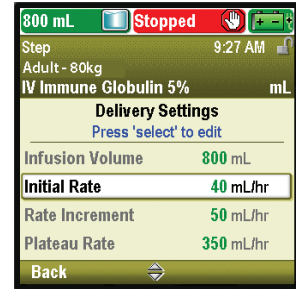This is the initial rate of medication delivery. Entering a new initial rate resets the infusion so that delivery starts at the beginning of the infusion duration.



Rate Increment

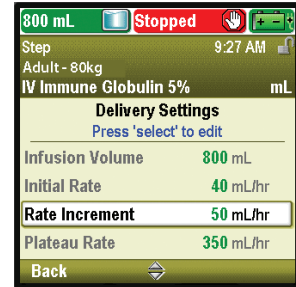This is the amount that you want medication delivery to increase by for each step. Entering a new rate increment resets the infusion so that delivery starts at the beginning of the infusion duration.




Plateau Rate

This is the maximum rate at which you want the medication to be delivered. The pump begins at the initial rate, and increases by the rate increment until the plateau rate is reached. The pump runs at the plateau rate until the infusion volume is fully delivered.

Entering a new plateau rate resets the infusion so that delivery starts at the beginning of the infusion duration.

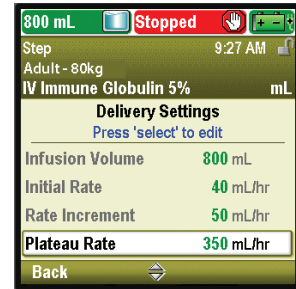
Rates above 250 mL/hr require a CADD® high volume administration set.



Step Duration

This is the length of time programmed for each step during medication delivery. After infusion occurs for the step duration, the rate increases by the rate increment unless the current rate is equal to the plateau rate.

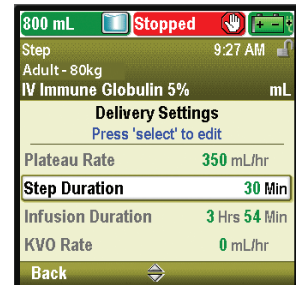Entering a new step duration resets the infusion so that delivery starts at the beginning of the infusion duration.



Infusion Duration

This is the time required to deliver the infusion volume. It is calculated by the pump based on the initial values programmed for infusion volume, initial rate, rate increment, step duration, and plateau rate. It appears on the screen for review only. Using Step Up or Step Down from the Advanced Tasks menu will change the actual infusion duration, although the displayed value will not change.

800 mL Stopped 9:28 AM
Step
Adult - 80kg
IV Immune Globulin 5% mL
Delivery Settings
Press 'select' to edit
Plateau Rate 350 mL/hr
Step Duration 30 Min
Infusion Duration 3 Hrs 54 Min
KVO Rate 0 mL/hr
Back

KVO Rate

The KVO or “keep vein open” rate is optional. It allows the delivery of a minimal amount of drug to help maintain catheter patency.

If a delayed start is programmed, the KVO rate is active until the infusion begins. It is also active after the infusion volume is delivered if the reservoir volume is programmed greater than the infusion volume.

800 mL Stopped 9:28 AM
Step
Adult - 80kg
IV Immune Globulin 5% mL
Delivery Settings
Press 'select' to edit
Plateau Rate 350 mL/hr
Step Duration 30 Min
Infusion Duration 3 Hrs 54 Min
KVO Rate 0 mL/hr
Back

Reservoir Vol.

This is the volume of fluid contained in the reservoir. The reservoir volume value decreases as the pump is primed or delivers fluid. Entering a new reservoir volume resets the infusion so that delivery starts at the beginning of the infusion duration.

The CADD[®]-Solis system administrator may choose a standard reservoir volume which allows the reservoir volume to be quickly reset to that value. The clinician can edit the reservoir volume to another amount. The reservoir volume should not be edited to less than the programmed infusion volume. This change is for the current protocol only and remains in effect until the Start New Patient or Same Patient New Protocol task is selected.

If KVO delivery is intended at the end of the infusion, the reservoir volume must be larger than the infusion volume to compensate for fluid used for priming, and so that automatic KVO delivery may occur. KVO delivery continues until the reservoir volume reaches 0 mL or until the pump is stopped. If the reservoir volume is the same as the infusion volume, KVO delivery will not occur at the completion of the infusion.

800 mL Stopped 9:28 AM
Step
Adult - 80kg
IV Immune Globulin 5% mL
Delivery Settings
Press 'select' to edit
Infusion Duration 3 Hrs 54 Min
KVO Rate 0 mL/hr
Reservoir Vol. 800 mL
Start Time 10:00 AM
Back

Start Time

This is the time that the next infusion will begin. It is displayed only if a delayed start is programmed.

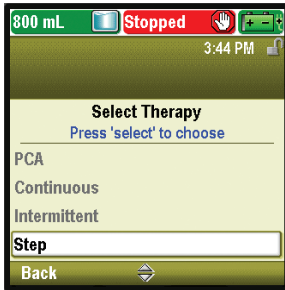
800 mL Stopped 9:28 AM
Step
Adult - 80kg
IV Immune Globulin 5% mL
Delivery Settings
Press 'select' to edit
Infusion Duration 3 Hrs 54 Min
KVO Rate 0 mL/hr
Reservoir Vol. 800 mL
Start Time 10:00 AM
Back

Programming Example

WARNING: Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is used to deliver critical or life-sustaining medication, the interruption in the delivery of medication could result in patient injury or death.


Example: A total of 700 mL of medication must be delivered to the patient. The initial dose must be given at 30 mL/hr, stepping up by 50 mL/hr in 30 minute increments until a plateau rate of 350 mL/hr is reached.

1. Select the Therapy




- Stop the pump if it is running.
- Select **Tasks**, then **View Advanced Tasks**, and then **Start New Patient**.
- Unlock the keypad.
- Press \uparrow or \downarrow to highlight the **Step** therapy and press **select**. See *Start New Patient* on page 27 if you need more information on selecting the therapy.

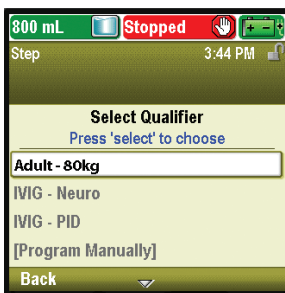2. Select the Qualifier




- Press \uparrow or \downarrow to highlight the desired qualifier and press **select**.

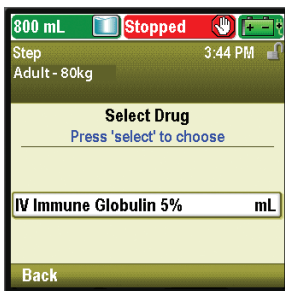3. Select the Drug




- Press \uparrow or \downarrow to highlight the desired drug, and press **select**.

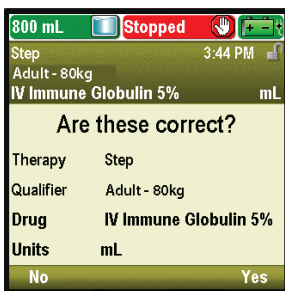4. Confirm and Review the Settings




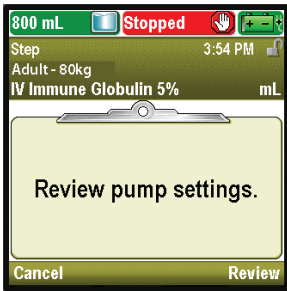- Review the information on the screen and select **Yes** if they are correct, or **No** to edit the settings.




- The pump displays the **Review pump settings** screen. Select **Review** to continue.

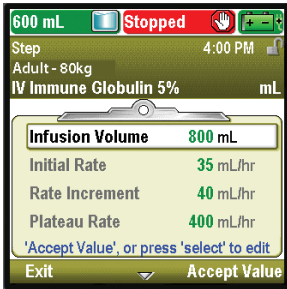
5. Enter the Infusion Volume



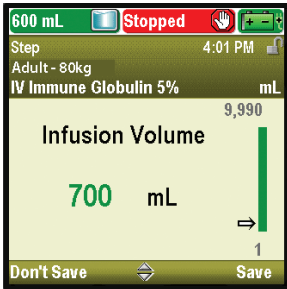- Press **▲** or **▼** to highlight **Infusion Volume**. If the value shown is set at the desired volume, select **Accept Value** and go to step 6. If the value shown does not show the desired volume, press **select**.




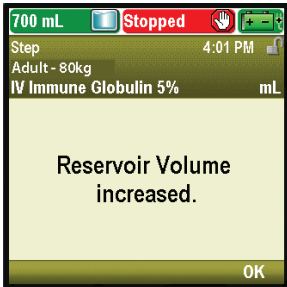- Press **▲** or **▼** until the value reads **700 mL** and select **Save**.




- If the infusion volume is greater than the displayed reservoir volume, you will be asked to confirm the change in the reservoir volume by selecting **OK**.
- The pump returns to the Review screen.

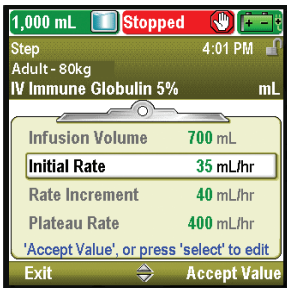6. Enter the Initial Rate




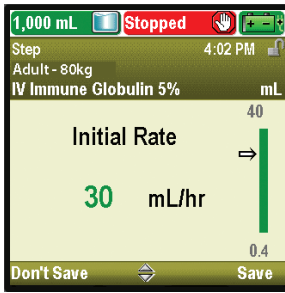- Press **▲** or **▼** to highlight **Initial Rate**. If the value shown is set at the desired rate, select **Accept Value** and go to step 7. If the value shown does not show the desired rate, press **select**.




- Press **▲** or **▼** until the value reads **30 mL/hr** and select **Save**. The pump then returns to the Review screen.

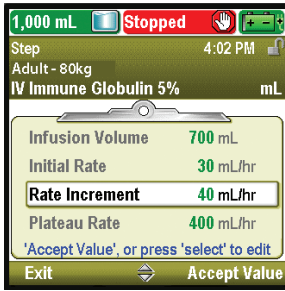7. Enter the Rate Increment




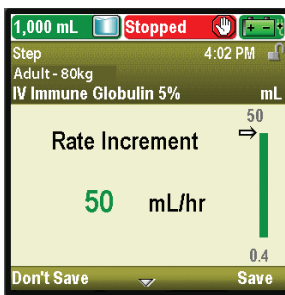- Press **▲** or **▼** to highlight **Rate Increment**. If the value shown is set at the desired increment, select **Accept Value** and go to step 8. If the value shown does not show the desired increment, press **select**.




- Press **▲** or **▼** until the value reads **50 mL/hr** and select **Save**.
- The pump then returns to the Review screen.

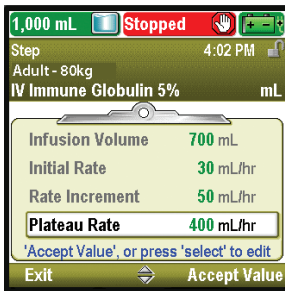8. Enter the Plateau Rate




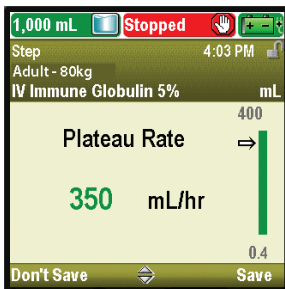- Press **▲** or **▼** to highlight **Plateau Rate**. If the rate shown is correct, select **Accept Value** and go to step 9. If the value is not correct, press **select**.




- Press **▲** or **▼** until the value reads **350 mL/hr** and select **Save**.
- The pump then returns to the Review screen.

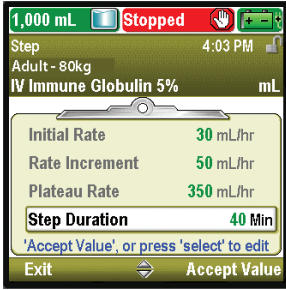9. Enter the Step Duration




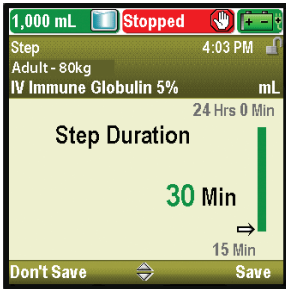- Press **▲** or **▼** to highlight **Step Duration**. If the value shown is correct, select **Accept Value** and go to step 10. If the value is not correct, press **select**.




- Press **▲** or **▼** until the value reads **30 Min** and select **Save**.
- The pump then returns to the Review screen.

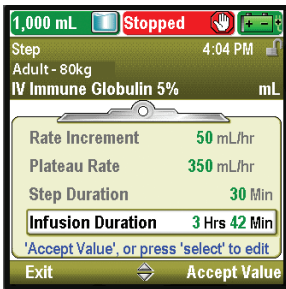10. Review the Infusion Duration




- Infusion Duration is a calculated value and cannot be edited.

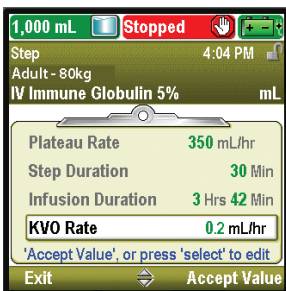11. Enter the KVO Rate




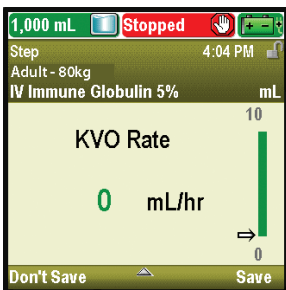- Press **▲** or **▼** to highlight **KVO Rate**. If the value shown is correct, select **Accept Value** and go to step 12. If the value is not correct, press **select**.




- Press **▲** or **▼** until the value reads **0 mL/hr** and select **Save**.
- The pump then returns to the Review screen.

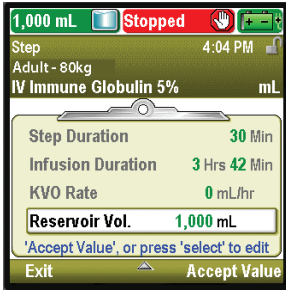12. Enter the Reservoir Volume




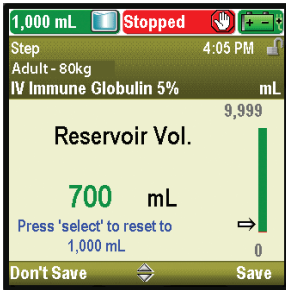- Press **▲** or **▼** to highlight **Reservoir Vol.** If the value shown is correct, select **Accept Value** and go to step 13. If the value is not correct, press **select**.



- Press **▲** or **▼** until the value reads **700 mL** and select **Save**.
- The pump then returns to the Review screen.

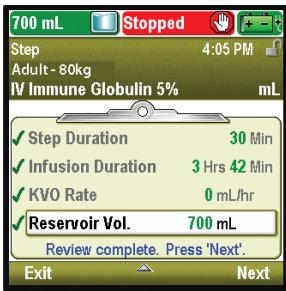13. Accept Values




- If any values have not already been accepted, press **▲** or **▼** to highlight each setting and select **Accept Value** if the setting is correct. A green check mark appears next to each accepted value.
- After your review is complete, select **Next**.

14. Prepare the Pump for the Patient

- Follow the instructions for attaching a cassette, priming, and attaching the pump to the patient in *Operating the Pump* on page 78).


Starting Each infusion

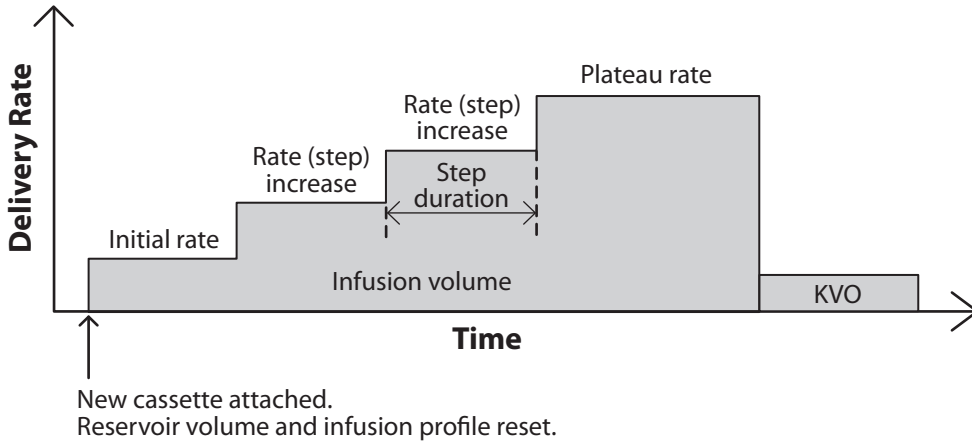When a new cassette is attached to the pump at the beginning of infusion, the reservoir volume should be reset. **This also resets the infusion.** Follow the instructions in *Operating the Pump* beginning on page 78, for attaching a cassette and resetting the reservoir. When the pump is started, delivery will start at the beginning of the infusion.



Stopping and Restarting the Pump During an Infusion

If you stop the pump before a step infusion was fully delivered, you may choose to either restart the infusion from where delivery left off, or you may reset the infusion.

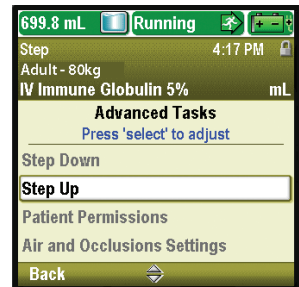
When you press stop/start , a screen appears with the message, “Infusion was interrupted before completion. Finish interrupted infusion or start from beginning with new bag.” Select **New Bag** to reset the reservoir volume and reset the infusion. Or, select **Finish** to resume the infusion from the point it was at when it was last running.

Step Up


This option allows you to increase the rate of infusion at a time other than what is programmed. You cannot step up to a higher rate than the plateau rate.

To step up from the current rate:

1. In the Advanced Tasks menu, press  or  to highlight **Step Up** and press .



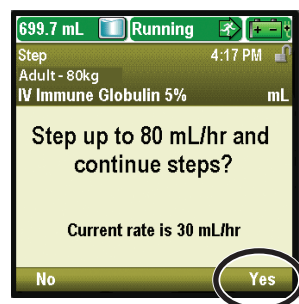2. Unlock the keypad.
3. The pump displays a screen that tells you the current rate, and asks you if you want to step up to the next step rate. Select **Yes**.




Step Down


This option allows you to reduce the rate of infusion at any time. You cannot step down lower than the initial rate. Once the step down task has been performed successfully, the automatic step increases after each step duration are cancelled. Delivery continues at the step down (hold) rate. The step up task can be used to restart the automatic step increases.

To step down from the current rate:

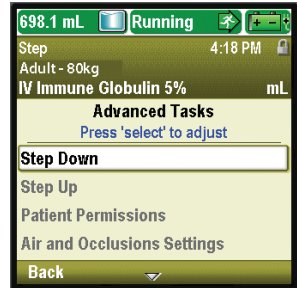1. In the Advanced Tasks menu, press  or  to highlight **Step Down** and press .



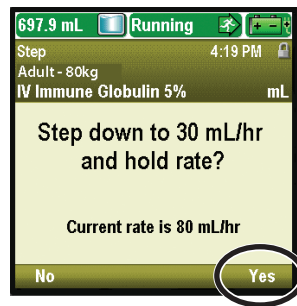
2. The pump displays a screen that tells you the current rate, and asks you if you want to step down to the previous step rate. Select **Yes**.



Taper Delivery Mode

Taper delivery allows the infusion of nutritional solutions (TPN, total parenteral nutrition) or other fluids, with optional tapering. Delivery can be gradually increased, or tapered up, at the beginning of the infusion, and it can be gradually decreased, or tapered down, at the end of the infusion. An optional KVO (keep vein open) rate may be delivered at the end of the infusion, depending on protocol settings.


The following graph illustrates an example taper delivery:

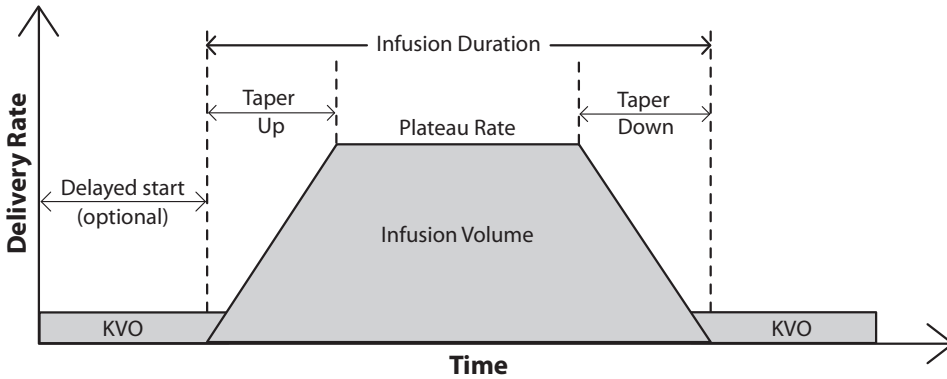

Home Screen

The delivery mode and status of the infusion are indicated on the home screen with messages appearing in text, a graphic in a shape unique to the therapy, and with the colors green, red, and gray. The colors indicate that the pump is running (green), or stopped (red). If a new patient or new infusion is programmed but *the pump has not been started*, or if a delayed start time was programmed and the start time passed while the pump was stopped and not yet started, the graphic is entirely gray. The location of the color on the graphic also indicates the status of the infusion. The graphic used for taper delivery also identifies if taper up and/or taper down are programmed as part of the delivery setting. The examples below are specific to the taper delivery mode.

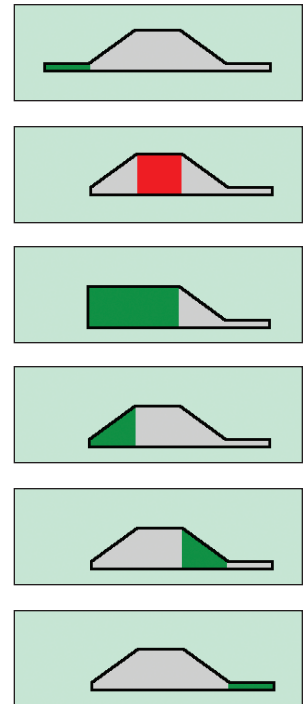- The pump is running with a programmed delayed start.
- The pump is stopped at the plateau rate.
- The pump is running at the plateau rate. No taper up was programmed.
- The pump is running at taper up.
- The pump is running at taper down.
- The pump is running at the KVO rate, nearing the end of the infusion.




Taper delivery mode graphics do **not** depict the relative slope of the taper or the relative duration of the tapers in comparison to the plateau.

Programming Screens


Infusion Volume

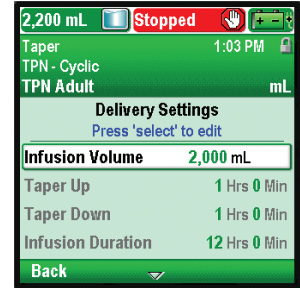This is the total volume of fluid to be delivered. If you enter an infusion volume that causes the plateau rate to exceed the maximum allowable rate, the pump automatically lengthens the infusion duration to accommodate the new volume. Additionally, the reservoir volume must be large enough to support the entire infusion volume. Therefore, programming the infusion volume may result in an automatic increase of the reservoir volume. Entering a new infusion volume resets the infusion so that delivery starts at the beginning of the infusion duration.




Taper Up

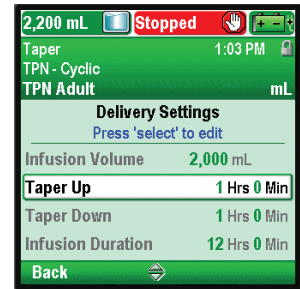This is the duration for the taper up period. The maximum programmable taper up is limited by a combination of the plateau rate, infusion duration, infusion volume, and taper down. You cannot select a taper up duration that would cause the plateau rate to exceed the maximum allowable rate. In addition, the infusion duration must be at least 10 minutes more than the combined total of the taper up and taper down periods. Therefore, changing taper up may result in the pump automatically lengthening the infusion duration to accommodate the new taper up duration. Entering a new taper up duration resets the infusion so that delivery starts at the beginning of the infusion duration.




Taper Down

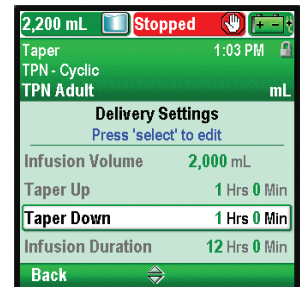This is the duration for the taper down period. The maximum programmable taper down is limited by a combination of the plateau rate, infusion duration, infusion volume, and taper up. You cannot select a taper down duration that would cause the plateau rate to exceed the maximum allowable rate. In addition, the infusion duration must be at least 10 minutes more than the combined total of the taper up and taper down periods. Therefore, changing taper down may result in the pump automatically lengthening the infusion duration to accommodate the new taper down duration. Entering a new taper down duration resets the infusion so that delivery starts at the beginning of the infusion duration.




Infusion Duration

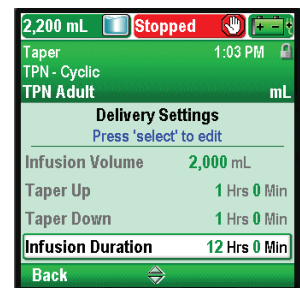This is the time required to deliver the infusion volume. The pump automatically calculates the rate necessary to deliver the infusion volume. The minimum programmable infusion duration is limited by a combination of the plateau rate, infusion volume, taper up, and taper down. You cannot select an infusion duration that would cause the plateau rate to exceed the maximum allowable rate. In addition, the infusion duration must be at least 10 minutes more than the combined total of the taper up and taper down periods. Entering a new infusion duration resets the infusion so that delivery starts at the beginning of the infusion duration.




Plateau Rate

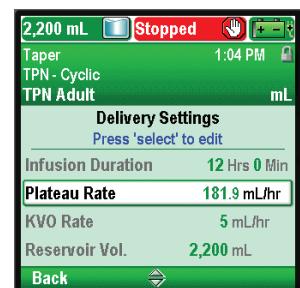This is the maximum rate at which medication is delivered during the plateau portion of the infusion. It appears on the screen for review only and is calculated by the pump. The plateau rate is based on the infusion volume, infusion duration, and any programmed tapering. Rates above 250 mL/hr require a CADD® high volume administration set.



KVO Rate

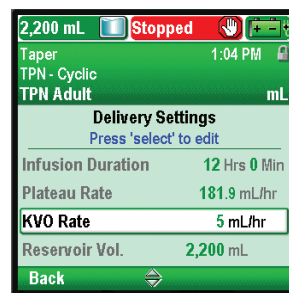
The KVO or “keep vein open” rate is optional. It allows the delivery of a minimal amount of drug to help maintain catheter patency. If a delayed start is programmed, the KVO rate is active until the infusion begins. It is also active after the infusion volume is delivered if the reservoir volume is programmed greater than the infusion volume.



Reservoir Vol.

This is the volume of fluid contained in the reservoir. The reservoir volume value decreases as the pump is primed or delivers fluid. Entering a new reservoir volume resets the infusion so that delivery starts at the beginning of the infusion duration.

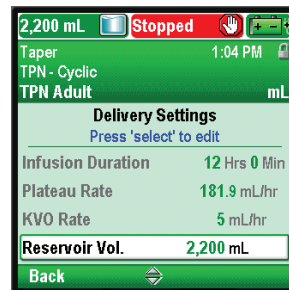The CADD®-Solis system administrator may choose a standard reservoir volume which allows the reservoir volume to be quickly reset to that value. The clinician can edit the reservoir volume to another amount that should not be less than the programmed infusion volume. This change is for the current protocol only and remains in effect until the Start New Patient or Same Patient New Protocol task is selected.



If KVO delivery is intended at the end of the infusion, the reservoir volume must be larger than the infusion volume to compensate for fluid used for priming, and so that automatic KVO delivery may occur. KVO delivery continues until the reservoir volume reaches 0 mL or until the pump is stopped. If the reservoir volume is the same as the infusion volume, KVO delivery will not occur at the completion of the infusion.

Start Time

This is the time that the next infusion begins. It is displayed only if a delayed start is programmed.



Taper Programming

Programming Example

WARNING: Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is used to deliver critical or life-sustaining medication, the interruption in the delivery of medication could result in patient injury or death.

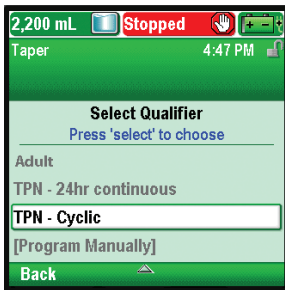
Example: A total of 1 L of TPN solution must be delivered to the patient over 12 hours. The TPN bag contains 1050 mL of solution. Delivery must taper up over 90 minutes at the beginning of the delivery, and taper down over 90 minutes at the end of delivery.

1. Select the Therapy




- Stop the pump if it is running.
- Select **Tasks**, then **View Advanced Tasks**, and then **Start New Patient**.
- Unlock the keypad.
- Press **▲** or **▼** to highlight the **Taper** therapy and press **select**. See *Start New Patient* on page 27 if you need more information on selecting the therapy.

2. Select the Qualifier




- Press **▲** or **▼** to highlight the desired qualifier and press **select**.

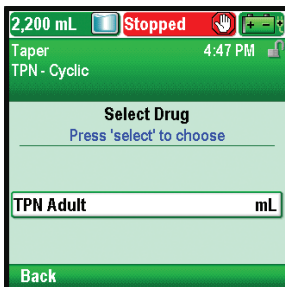3. Select the Drug



- Press **▲** or **▼** to highlight the desired drug, and press **select**.

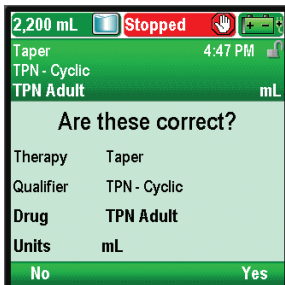4. Confirm and Review the Settings




- Review the information on the screen and select **Yes** if they are correct, or **No** to edit the settings.



- The pump displays the **Review pump settings** screen. Select **Review** to continue.

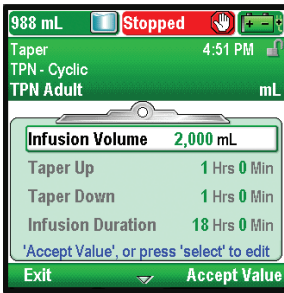
5. Enter the Infusion Volume



- Press **▲** or **▼** to highlight **Infusion Volume**. If the volume shown is set at the desired value, select **Accept Value** and go to step 6. If the value shown does not show the desired volume, press **select**.




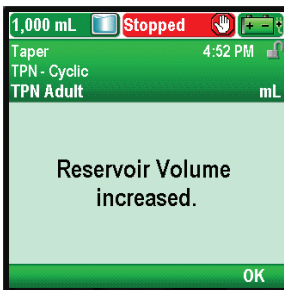- Press **▲** or **▼** until the value reads **1000 mL** and select **Save**.




- If the infusion volume is greater than the displayed reservoir volume, you will be asked to confirm the change in the reservoir volume by selecting **OK**.
- The pump returns to the Review screen.

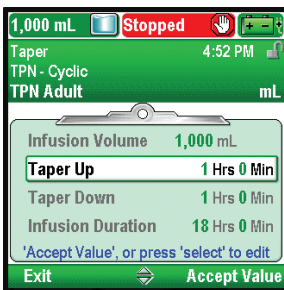6. Enter Taper Up




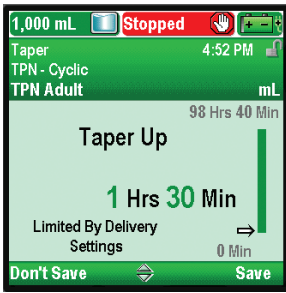- Press **▲** or **▼** to highlight **Taper Up**. If the value shown is correct, select **Accept Value** and go to step 7. If the value shown is not correct, press **select**.




- Press **▲** or **▼** until the value reads **1 Hrs 30 Min** and select **Save**. The pump then returns to the Review screen.

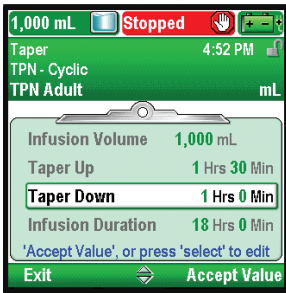7. Enter Taper Down




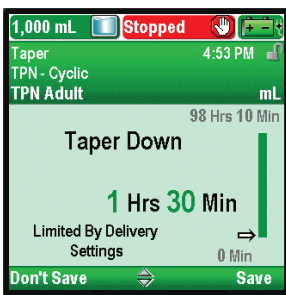- Press **▲** or **▼** to highlight **Taper Down**. If the value shown is correct, select **Accept Value** and go to step 8. If the value shown is not correct, press **select**.




- Press **▲** or **▼** until the value reads **1 Hrs 30 Min** and select **Save**. The pump then returns to the Review screen.

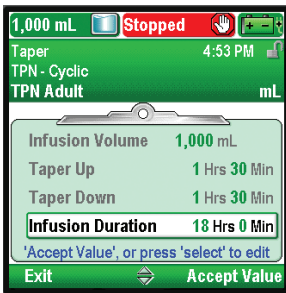8. Enter the Infusion Duration




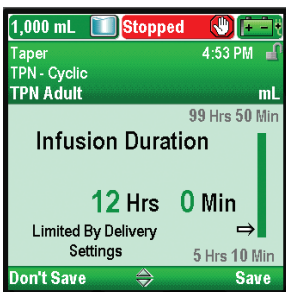- Press **▲** or **▼** to highlight **Infusion Duration**. If the value shown is correct, select **Accept Value** and go to step 9. If the value shown is not correct, press **select**.




- Press **▲** or **▼** until the duration reads **12 Hrs 0 Min** and select **Save**.
- The pump then returns to the Review screen.

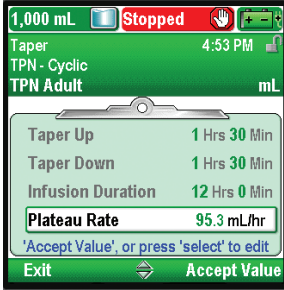9. Review the Plateau Rate




- The plateau rate is a calculated value and cannot be edited.

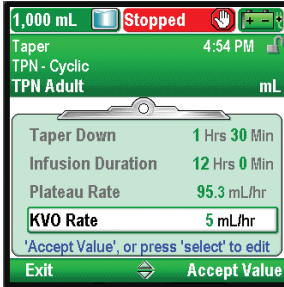10. Enter the KVO Rate




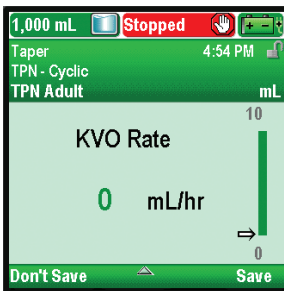- Press **▲** or **▼** to highlight **KVO Rate**. If the value shown is correct, select **Accept Value** and go to step 11. If the value is not correct, press **select**.




- Press **▲** or **▼** until the value reads **0 mL/hr** and select **Save**.
- The pump then returns to the Review screen.

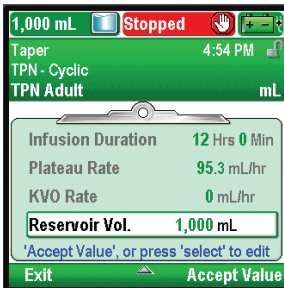11. Enter the Reservoir Volume




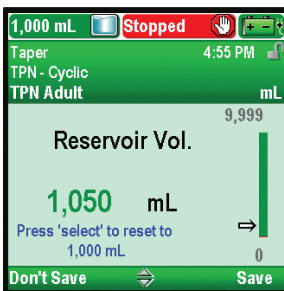- Press **▲** or **▼** to highlight **Reservoir Vol.** If the value shown is correct, select **Accept Value** and go to step 12. If the value is not correct, press **select**.



- Press **▲** or **▼** until the value reads **1050 mL** and select **Save**.
- The pump then returns to the Review screen.

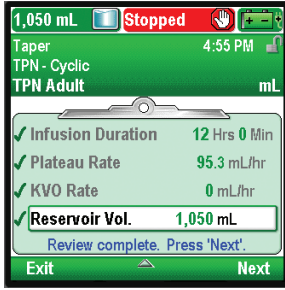12. Accept Values




- If any values have not already been accepted, press \blacktriangle or \blacktriangledown to highlight each setting and select **Accept Value** if the setting is correct. A green check mark appears next to each accepted value.
- After your review is complete, select **Next**.

13. Prepare the Pump for the Patient

- Follow the instructions for attaching a cassette, priming, and attaching the pump to the patient in *Operating the Pump* on page 78).


Starting Daily Infusion

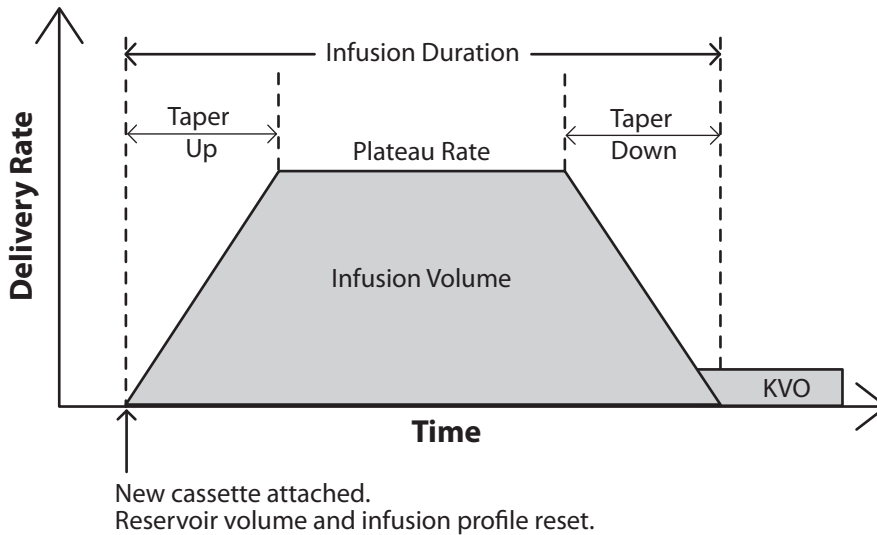When a new cassette is attached to the pump at the beginning of infusion, the reservoir volume should be reset. **This also resets the infusion.** Follow the instructions in *Operating the Pump* beginning on page 78), for attaching a cassette and resetting the reservoir. Then when the pump is started, delivery will start at the beginning of the infusion.



Stopping and Restarting the Pump During an Infusion

If you stop the pump before a taper infusion was fully delivered, you may choose to either restart the infusion from where delivery left off, or you may reset the infusion.

When you press stop/start , a screen appears with the message, “Infusion was interrupted before completion. Finish interrupted infusion or start from beginning with new bag.” Select **New Bag** to reset the reservoir volume and reset the infusion. Or, select **Finish** to resume the infusion from the point it was at when it was last running.



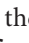
Taper Down Now



This task allows you to stop an infusion early by immediately tapering down. The following conditions are necessary for Taper Down Now:


- The pump must be running.
- The pump must be delivering at the plateau rate or already be in taper down.

To access the Taper Down Now option:


1. From the Tasks menu (see page 88), press  or  until **Taper Down Now** is highlighted, and then press .

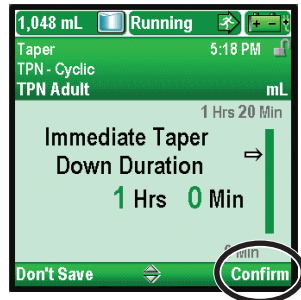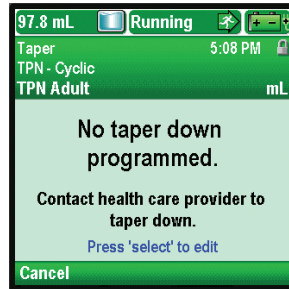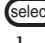NOTE: If a taper down was not programmed, you can edit the setting by pressing . Enter the security code, and then press  or  to set the taper down duration. Press **Confirm** to accept the change.



2. Select **Yes** to continue with taper down now or press  to edit the taper down period. Editing the taper down now period only affects the period for this immediate taper. It has not effect on the taper period of the protocol that is seen in the delivery settings.

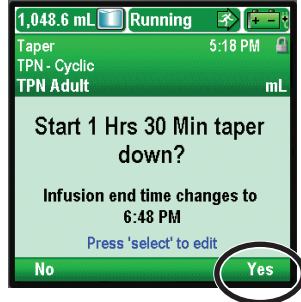NOTE: Once you select **Yes**, you can *not* restart delivery at the plateau rate without resetting the infusion period.



3. When “Beginning Taper Down...” appears, select **OK**.



Operating the Pump

Batteries

Four AA, 1.5 volt primary (non-rechargeable) alkaline batteries or the CADD[®]-Solis rechargeable battery pack are recommended for use in the CADD[®]-Solis VIP ambulatory infusion pump.

NOTE: Smiths Medical does not recommend mixing new and used batteries because it may affect low battery alarm times. Always use 4 new batteries when replacing depleted ones.

CAUTION: Inspect the AA batteries for damage or wear to the metal or plastic insulation prior to use, or after the pump has been dropped or hit. Replace the batteries if any damage is noted.

The pump retains all programmed values when the batteries are removed, but the batteries must be in place during delivery. If the batteries are removed while the pump is delivering, and an AC adapter is connected, delivery stops. If an AC adapter is not connected and the batteries are removed, delivery stops and the pump loses power.

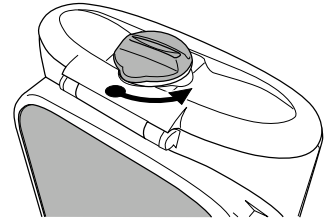
Dispose of used batteries in an environmentally safe manner, and according to any regulations which may apply.

WARNING:


- Do not use rechargeable NiCd or nickel metal hydride (NiMH) batteries. Do not use carbon zinc (“heavy duty”) batteries. They do not provide sufficient power for the pump to operate properly, which could result in death or serious injury to the patient.
 - Always have new batteries available for replacement. If power is lost, nondelivery of drug occurs and depending on the type of drug being administered, could result in death or serious injury to the patient.
 - A rechargeable battery pack that has reached the end of its useful life must be replaced with either another CADD[®]-Solis rechargeable battery pack or with 4 AA batteries. Using rechargeable battery packs from other manufacturers could result in fire or explosion.
 - There is no pump alarm to alert users that a battery has not been properly installed. An improperly installed battery could result in loss of power and nondelivery of drug and depending on the type of drug being administered, could result in death or serious injury to the patient.
 - Always check the battery compartment for fluid or debris before inserting the batteries, and do not allow any fluid or debris to fall into the battery compartment. Fluid or debris in the battery compartment may damage the battery contacts and could result in loss of power and nondelivery of drug and depending on the type of drug being administered, could result in death or serious injury to the patient.
 - If the pump is dropped or hit, the battery door may become broken or damaged. Do not use the pump if the battery door is damaged because the batteries will not be properly secured. This may result in loss of power, nondelivery of drug and depending on the type of drug being administered, could result in death or serious injury to the patient.
-

Install the Batteries or Battery Pack

1. Make sure the pump is stopped and powered off. Using your fingers, the pump key, or a coin, turn the knob on the battery door counterclockwise to open the battery door.




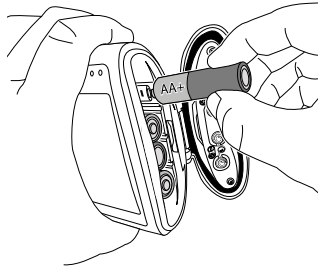2. Hold the pump at an angle and place 4 AA batteries in the pump, from the bottom up (see picture). Match the + and – markings on the new batteries with the markings on the pump.



OR: If using a rechargeable battery pack, insert it into the pump as shown.

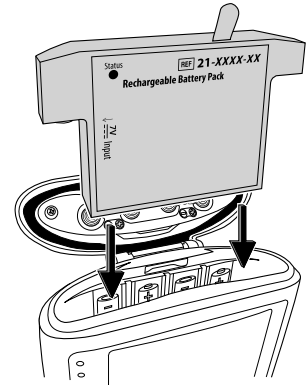OR



3. Close the battery door and using your fingers, the pump key, or a coin, turn the knob on the battery door clockwise to lock.

NOTE: If you put the batteries in backwards, the pump will not power up. Check the batteries, making sure to match the + and – markings.

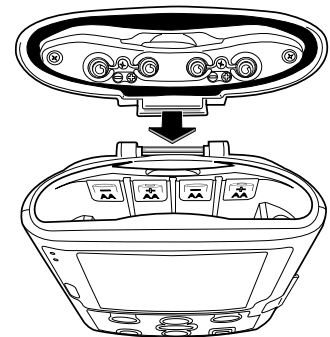
CAUTION: Do not store the pump for prolonged periods with the batteries installed. Battery leakage could damage the pump.

NOTES:

- Battery life is dependent on the amount of medication delivered, delivery rate, battery age, temperature, active display time, and backlight intensity.
- Battery power is quickly depleted at temperatures below 10°C (50°F).

Replace the Battery Door

If the battery door is removed or needs replacing, simply snap the door onto the bar that is located on the pump.



Power Up

Press and hold the power switch to turn the pump on. The pump starts the power up sequence during which it performs various self-tests, tests for alarm conditions. Watch for the following during power-up:

- Both the green and amber indicator lights flash.
- The display quickly flashes gray, then blue. An amber swirl then fills the display, followed by a CADD®-Solis VIP Ambulatory Infusion System display. Look for any stripes or black or white pixels which indicate a faulty display. If you see any indication of a faulty display, remove the pump from service and contact Smiths Medical Customer Service.
- After the power up is completed, listen for the Morse Code “OK” sound (a series of six audible beeps). If you do not hear this sound, there may be a problem with the audible alarms. If you believe there is a problem, remove the pump from service and contact Smiths Medical Customer Service.
- If any issues are found while the pump is performing the self tests, alarms will sound (for example, if the battery is low or a key on the keypad is stuck in the pressed position).

CAUTION: If the power up results in an error message indicating that the protocol library was lost, do not proceed with using the pump. Follow your facility’s procedures for downloading protocol libraries.

Power Down

Press the power switch to turn the pump off. You must acknowledge the “Power down?” prompt by selecting **Yes**. The pump then powers down. Be aware of the following conditions when the pump is powered down:

- The display is blank.
- Keypad presses are not detected.
- The remote dose cord does not function.
- Medication is not infused.
- Alarm conditions are not sensed.
- Audio and visual alarms are not displayed.
- Communication with CADD®-Solis Medication Safety Software is not possible.
- If the AC adapter is plugged in, the blue AC power light is on. No other lights are lit.
- The internal battery used to retain time and date does not charge.
- If the pump is connected to AC power and a rechargeable battery pack is installed, the battery pack continues to charge.

Cassettes

The cassette is the part of the CADD™ medication cassette reservoir or CADD® administration set that attaches to the bottom of the pump. Either a CADD™ medication cassette reservoir used with a CADD® extension set, or a CADD® administration set are compatible with the pump.

WARNING: Follow the instructions for use provided with the CADD™ medication cassette reservoir, CADD® extension set, or CADD® administration set, paying particular attention to all warnings and cautions associated with their use. Incorrect preparation and/or use of these products could result in serious patient injury or death.

NOTE:

- A CADD® set with free flow protection must be used.
- A CADD® high volume administration set is required for rates above 250 mL/hr.
- Appropriate supplies needed to replace cassettes should be available if a cassette is damaged.
- For detailed instructions on preparing the product for use, refer to the instructions for use supplied with the product.

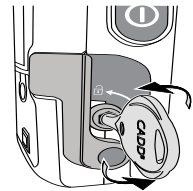
Remove a Cassette

1. Make sure the pump is stopped before removing the cassette.
2. Close the tubing clamp.


WARNING: Per general rules of safe practice, always clamp the tubing before removing the cassette from the pump. Removing the cassette without closing the clamp could potentially cause unregulated gravity infusion, which could result in patient injury or death.

3. If locked, insert the pump key and turn the cassette/keypad lock counter-clockwise into the unlocked position. “Cassette Unlocked” briefly appears in the status bar.

NOTE: Only PCA infusions require the cassette to be locked in order to run the pump.



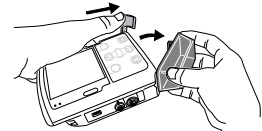4. Push down on the cassette latch until the cassette detaches.



Attach a Cassette

Obtain a new, filled CADD™ medication cassette reservoir, or CADD® administration set attached to a flexible IV bag. Refer to the instructions for use supplied with the product for information on preparing the product for use.

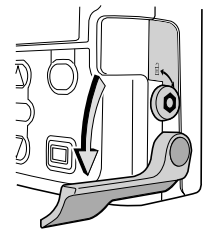
WARNING: Attach the cassette properly. A detached or improperly attached cassette could result in unregulated gravity infusion of medication from the reservoir or a reflux of blood, which could result in death or injury to the patient.

If you are using a CADD® administration set or CADD™ medication cassette reservoir that does not have the Flow Stop feature, you must use a CADD® extension set with anti-siphon valve, or a CADD® administration set with either an integral or add-on anti-siphon valve to protect against unregulated gravity infusion that can result from an improperly attached cassette. Unregulated gravity infusion can result in death or serious injury.




CAUTION: If you are using a CADD™ medication cassette reservoir in which the medication is frozen, thaw at room temperature only. Do not heat in a microwave oven as this may damage the product and cause leakage.

Before attaching a new cassette, make sure the pump is powered on. Once the cassette is attached, the pump automatically displays screens that allow you to verify the cassette type, prime the fluid path, reset the reservoir volume, review pump settings, and/or start the pump.

1. Clamp the tubing.
2. Make sure the cassette latch is unlocked and open the cassette latch.
3. Insert the cassette hooks into the hinge pins on the bottom of the pump.




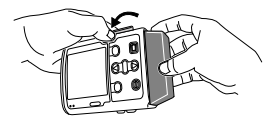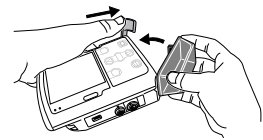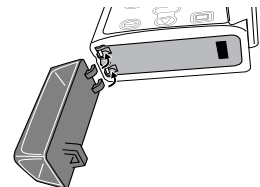4. Push down on the cassette latch, and push up on the cassette until it firmly clicks into place.
5. Lift the cassette latch into the closed position. A message briefly appears in the status bar so you can verify the type of cassette you have attached.



6. If you wish to lock the cassette, insert the pump key into the cassette/keypad lock and turn clockwise into the locked position. "Cassette Locked" appears briefly in the status bar.

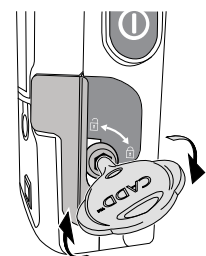NOTE: For PCA infusions, the cassette must be locked in order to start the pump.



Prime Tubing

When priming the fluid path, the tubing downstream of the pump is filled with fluid, removing any air bubbles. Prime the tubing *before* connecting it to the patient's infusion set or indwelling catheter.

Fluid delivered by priming is subtracted from the reservoir volume, but is not added to the amount given because this fluid is not delivered to the patient. Priming is not allowed when the reservoir volume is 0 mL.

NOTE: The air detector is disabled while the pump is priming.

NOTE: If the fluid path contains an air eliminating filter, it is acceptable for air bubbles to be present on the vent side of the filter. See the instructions for use supplied with the disposable sets for more information.

Prime Tubing After Changing a Cassette

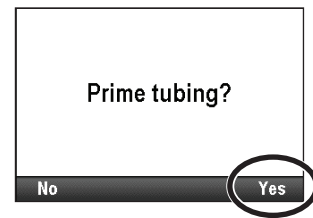
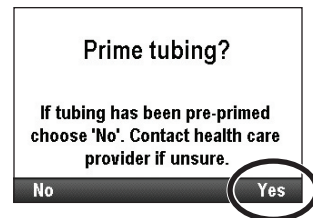
NOTE: If you are not changing the cassette but wish to prime the fluid path, use the *Prime Tubing, No Cassette Change* task described on page 84.

1. When a cassette is attached after the pump is powered on, a **Prime Tubing?** screen always appears.

When priming security is set to *On* under patient permissions, a screen similar to the one on the right appears, and a security code is required to proceed. If you do not have permission to prime the tubing or if the tubing is pre-primed, select **No**. For additional information, see *Patient Permissions* on page 100.


If you have permission to prime the tubing, select **Yes**. Unlock the keypad.

When priming security is set to *Off* under patient permissions, a screen like the one to the right appears. Select **Yes**.




2. If you have not already done so, disconnect the tubing from the patient, open the clamps, and select **Prime**.

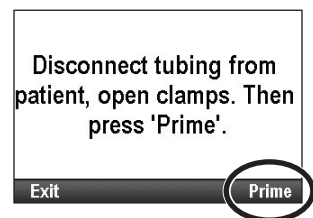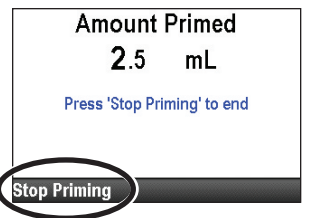WARNING: Do not prime the fluid path with the tubing connected to a patient as this could result in overdelivery of medication or air embolism, which could result in serious patient injury or death.






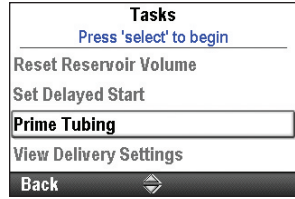
3. Stop priming at any time by selecting **Stop Priming**. Priming automatically stops after 10 mL (or 20 mL if a high volume set is attached) are primed. Continue priming as needed.

WARNING: To prevent air embolism, ensure that the entire fluid path is free of all air bubbles before connecting to the patient. Air embolism could result in serious patient injury or death.



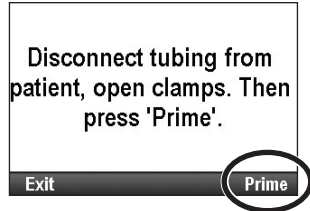
Prime Tubing, No Cassette Change

1. Stop the pump if it is running.
2. In the Tasks menu (see page 88), press  or  until **Prime Tubing** is highlighted, and then press .




3. Unlock the keypad, if required.
4. If you have not already done so, disconnect the tubing from the patient, open the clamps, and select **Prime**.

WARNING: Do not prime the fluid path with the tubing connected to a patient as this could result in overdelivery of medication or air embolism, which could result in serious patient injury or death.



5. Stop priming at any time by selecting **Stop Priming**. Priming automatically stops after 10 mL (or 20 mL if a high volume set is attached) are primed. Continue priming as needed.

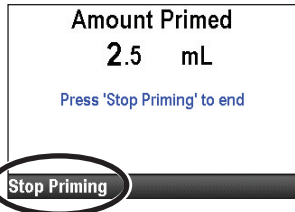WARNING: To prevent air embolism, ensure that the entire fluid path is free of all air bubbles before connecting to the patient. Air embolism could result in serious patient injury or death.



Start the Pump

Infusion begins when the pump starts. When the pump is running, “Running” appears on the status bar, the graphic on the home screen is green, and the green indicator light flashes. If the pump will not start, a message appears on the display. Refer to *Alarms and Messages, Alphabetical List* on page 112.

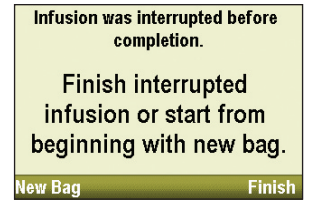
NOTE: Before starting the pump, be sure the tubing is primed and the pump is connected to the patient according to your facility’s standards of practice.

1. Press stop/start .

NOTE: If the delivery settings have not been reviewed and the values have not been accepted, you must do so before the pump will run. Instructions for reviewing settings are available on page 36.




If a step or taper infusion was previously stopped and delivery was not completed (the infusion was not fully delivered), a screen appears with the message, “Infusion was interrupted before completion. Finish interrupted infusion or start from beginning with new bag.” Select **New Bag** to reset the reservoir volume and reset the infusion. Or, select **Finish** to resume the infusion from the point it was at when it was last running.



2. When “Start Pump?” appears, select **Yes**.

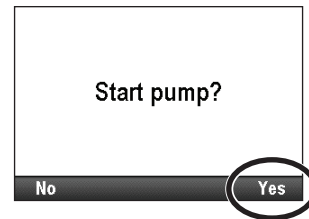
NOTE: To set a delayed start time, refer to *Delayed Start* on page 31 for all therapies except intermittent, or *Next Dose Start Time* on page 32 for intermittent therapies.



3. The pump begins running. The red “Stopped” message in the status bar changes to a green “Running” message, and “Infusion is starting now...” appears briefly on the screen.



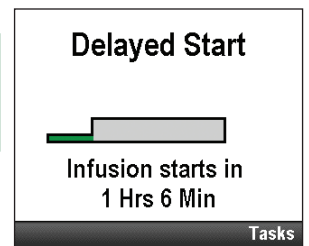If a delayed start was programmed, the display turns green and a message that the infusion is delayed appears briefly on the screen. Then the “Delayed Start” screen appears along with the time remaining until the infusion starts. The pump infuses at the programmed KVO rate.




Stop the Pump



Stopping the pump stops delivery. After the pump is stopped, “Stopped” appears in red on the status bar, the graphic on the home screen is red, the amber indicator light flashes, and the green indicator light is off.

1. Press stop/start .

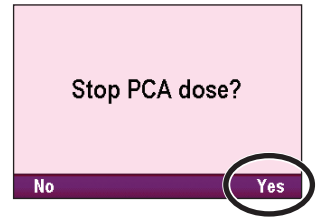In a PCA protocol, if a PCA dose is in progress, “Stop PCA dose?” appears. Select **Yes** to stop the dose.



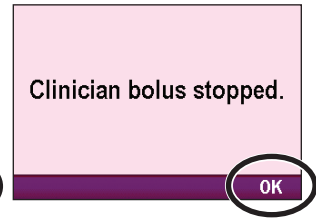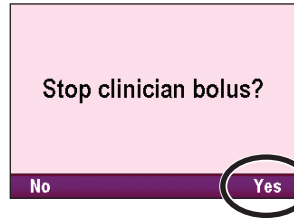
In a PCA protocol, if a clinician bolus is in progress, “Stop clinician bolus?” appears. Select **Yes** to stop the bolus. A confirmation screen appears stating, “Clinician bolus stopped.” Select **OK** to continue.



In a taper protocol, “Taper down instead?” appears. Select **No** to stop the infusion (see *Taper Down Now* on page 77 if you want to taper down rather than stop the pump).



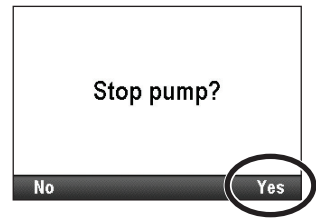2. When “Stop Pump?” appears, select **Yes**.



3. The pump stops running. The green “Running” message in the status bar changes to a red “Stopped” message, and “Pump is stopping...” appears briefly on the screen.



Reset Reservoir Volume

The reservoir volume setting indicates the amount of fluid contained in the reservoir. Once this number is set, the pump keeps track of how much fluid has been delivered and adjusts the reservoir volume setting accordingly.

WARNING: Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is used to deliver critical or life-sustaining medication, the interruption in the delivery of medication could result in patient injury or death.

When you latch and/or lock a cassette onto the pump as described on page 82, a series of messages leads you through resetting the reservoir volume, priming the tubing, and starting the pump. You can, however, reset the reservoir volume without changing the cassette by using the Reset Reservoir Volume task.

The CADD[®]-Solis system administrator may set a standard reservoir volume for each therapy that allows you to quickly reset the reservoir volume to that value. To reset the reservoir volume to a quantity other than the default setting, refer to *Editing Delivery Settings* on page 35.

NOTE: If you are running a step or taper therapy, resetting the reservoir volume also resets the infusion. When you restart the pump, delivery starts at the beginning of the infusion duration.

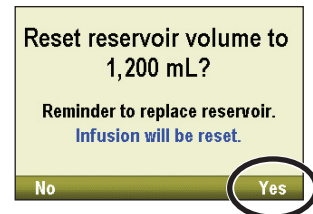
To reset the reservoir volume after attaching a new cassette:

1. The pump displays a question asking if you wish to reset the volume to the default amount. If this screen does not appear, the reservoir volume may already be reset.

In this example, selecting **Yes** resets the volume to 1,000 mL.




2. If you are running a step or taper therapy, you also see, “Infusion will be reset”. Select **Yes** to reset the volume and infusion. Select **No** to keep the reservoir volume at the current setting.



To reset the reservoir volume without changing the cassette:

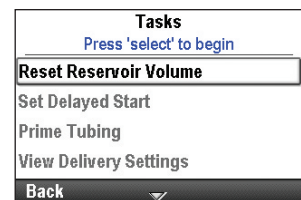1. Stop the pump if it is running
2. In the Tasks menu (see page 88), press or until **Reset Reservoir Volume** is highlighted, and then press .



3. The pump displays a screen asking you to confirm that you want to reset the reservoir volume. Select **Yes**.



Tasks and Advanced Tasks

Tasks Menu Overview


The Tasks menu leads to most of the pump’s operating functions. Some items on the Tasks menu are protected by the various security levels. To learn more about the security codes, refer to *Security Settings* on page 18.

The following functions are on the Tasks menu.

Task Menu Function	Description on Page
Reset Reservoir Volume	87
Taper Down Now (Taper therapy)	77
Set Delayed Start (all therapies except Intermittent)	31
Set Next Dose Start Time (Intermittent therapy)	32
Prime Tubing	83
View Delivery Settings	34
Display and Sound Settings	89
Change Time and Date	92
View Reports	95
View Advanced Tasks	99

To access the Tasks menu:

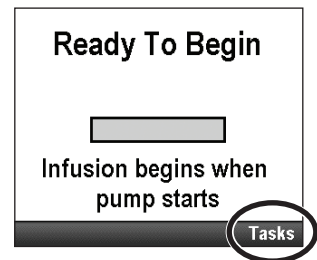
Select **Tasks** from the home screen.



Display and Sound Settings

The Display and Sound Settings menu allows you to adjust various factors such as the intensity of the backlight on the screen, alarm volumes, sound themes, key beeps, and numeric formats. This change is for the current protocol only and remains in effect until the Start New Patient or Same Patient New Protocol task is selected.

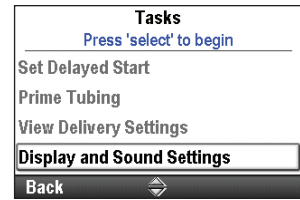In the Tasks menu, press  or  to highlight **Display and Sound Settings** and press .








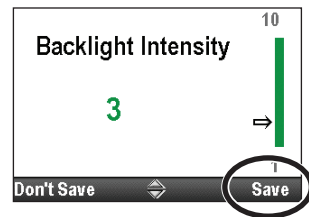
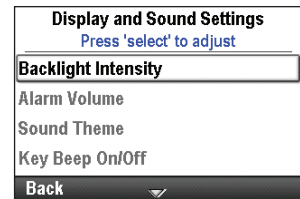
Backlight Intensity

The backlight intensity feature allows you to adjust the screen backlight intensity within the range of 1 to 10, with 1 being low and 10 being high.

NOTE: Increasing the backlight intensity shortens the battery life. For example, an increase of the backlight intensity from 5 to 10, reduces the battery life by 25%.

To adjust the backlight intensity:






- In the Display and Sounds Settings menu, press  or  to highlight **Backlight Intensity** and press .
- Press  or  to scroll from 1 to 10. The pump displays the intensity of each number as it appears. Once you have the desired backlight intensity, select **Save**.

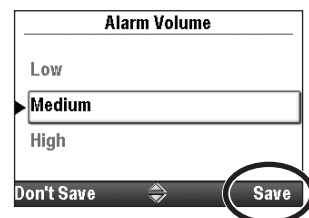
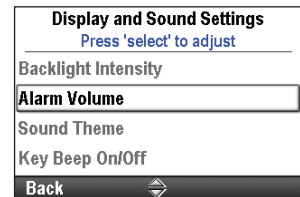


Alarm Volume

The alarm volume feature allows you to set the volume of the alarms in the therapy. You may choose between 3 volumes: low, medium, and high.

To adjust the alarm volume:

- In the Display and Sounds Settings menu, press  or  to highlight **Alarm Volume** and press .
- Press  or  to highlight **Low**, **Medium**, or **High**, and select **Save**.







Sound Theme

The Sound Theme setting allows you to choose between 3 different sound themes for the alarms and beeps that the pump makes. The 3 themes are standard, intense, and distinctive. You may preview the sound themes in the menu. The alarms and beeps that use the sound theme:

- High Priority Alarm
- Medium Priority Alarm
- Low Priority Alarm
- Informational Alarm

To choose a sound theme:

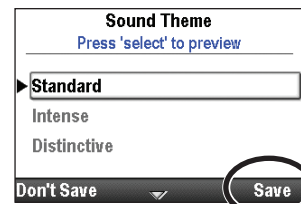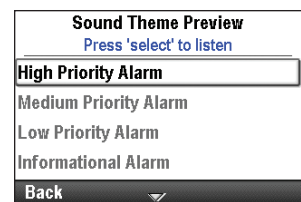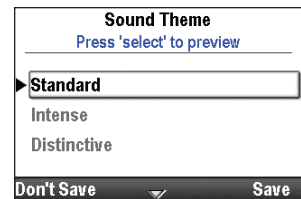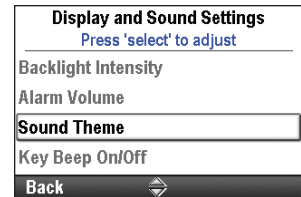1. In the Display and Sounds Settings menu, press **▲** or **▼** to highlight **Sound Theme** and press **select**.
2. Press **▲** or **▼** to highlight **Standard**, **Intense**, or **Distinctive**, and press **select** to preview the theme, or select **Save** to save the theme without a preview.
3. To preview the theme, press **▲** or **▼** to scroll through the alarms and beeps. Press **select** to listen to each alarm or beep. Continue until you have decided which one to apply. To exit, select **Back**.
4. Highlight the desired theme and select **Save**.








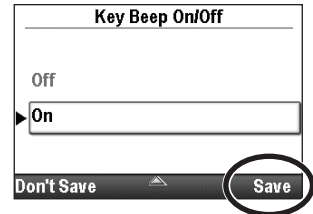
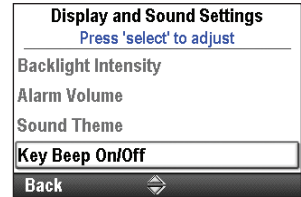
Key Beep On/Off

The key beep on/off feature allows you to turn on or off the audible beep that accompanies each key press. This feature **does not** turn off any audible alarms associated with alarm or alert conditions of the pump.

NOTE: If the key beep feature is turned on, the key beeps will not sound while any alarm screens are being displayed or when entering a security code.

To turn the key beeps on or off:

1. In the Display and Sounds Settings menu, press  or  to highlight **Key Beep On/Off** and press .
2. Press  or  to highlight **Off** or **On** and select **Save**.








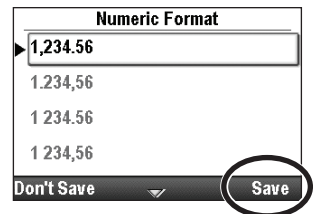
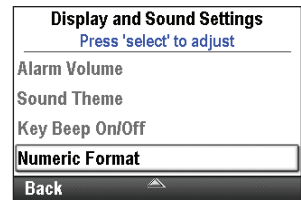
Numeric Format

The numeric format setting allows you to select the character or space used to indicate the decimal point and thousands units separators. The formats available are:

- 1,234.56
- 1.234,56
- 1 234.56
- 1 234,56

To modify the numeric format:

1. In the Display and Sounds Settings menu, press  or  to highlight **Numeric Format** and press .
2. Unlock the keypad.
3. Press  or  to highlight the format you want and select **Save**.



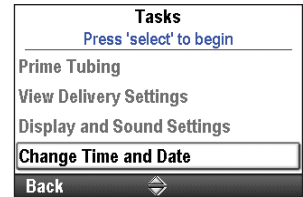
Time and Date

The Change Time and Date setting allows you to edit the time and date, and to choose a time and date format. The time and date options should reflect the current time and date. The clock is powered by a separate, internal battery, which retains the time and date even when the 4 AA batteries or the battery pack is removed. The pump uses this feature to record the time and date of events in the delivery and event logs, and in other reports. It is also used to determine when delayed starts and next dose start times begin, and when preventive maintenance alarms occur.

NOTE: The internal battery is rechargeable and automatically recharges when the pump is powered on. If the pump is off for a long period of time, it may not show the correct time and date when it is powered on. Check the time and date, and edit these settings if necessary.













CAUTION: If delivery of an infusion is affected by a time or date change, an alarm message appears and must be confirmed.

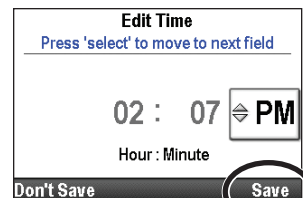
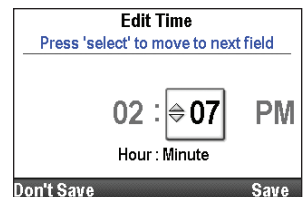
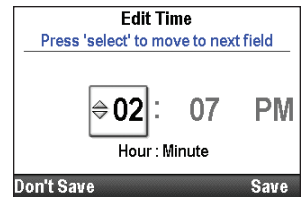
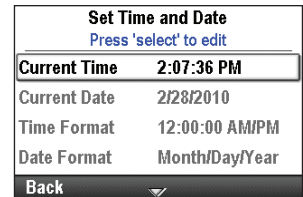
In the Tasks menu, press  or  to highlight **Change Time and Date** and press .



Current Time

To set the current time:


1. In the Change Time and Date menu, press  or  to highlight **Current Time** and press .
2. Press  or  to scroll to the correct hour and press  to navigate to the minutes.
3. Press  or  to scroll to the correct minutes and press  to navigate to the AM or PM field.
4. Press  or  to scroll between AM and PM and press  when the correct setting appears. If you are using a 24-hour clock, the AM/PM field does not appear. Select **Save**.



Current Date


To set the current date:

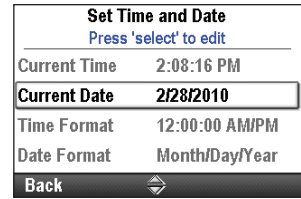1. In the Change Time and Date menu, press or to highlight **Current Date** and press .




2. Press or to scroll to the correct month and press to navigate to the day.

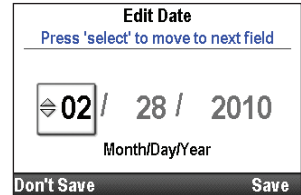NOTE: This example is for the Month/Day/Year date format. The date format is indicated below the date so that you know which field is changing.



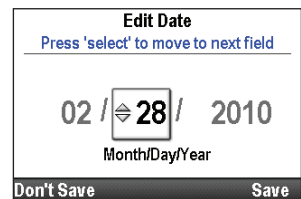3. Press or to scroll to the correct day and press to navigate to the year.



4. Press or to scroll to the correct year and select **Save**.



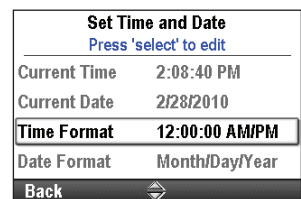
Time Format

The time format setting allows you to select either a 12-hour or 24-hour clock format, and to set the character used for indicating the hours, minutes and seconds separators. The formats available are as follows:

- 12:00:00 AM/PM
- 12.00.00 AM/PM
- 24:00:00
- 24.00.00


To set the time format:

1. In the Change Time and Date menu, press or to highlight **Time Format** and press .



2. Unlock the keypad.

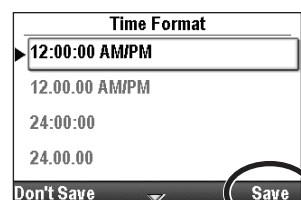3. Press or to highlight the desired time format and select **Save**.






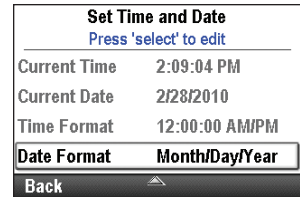
Date Format



The date format option allows you to select 1 of 3 formats for indicating the day, month, and year. The formats available are as follows:

- Month/Day/Year
- Day/Month/Year
- Year/Month/Day

To set the date format:

1. In the Change Time and Date menu, press  or  to highlight **Date Format** and press .



2. Unlock the keypad.
3. Press  or  to highlight the desired date format and select **Save**.



Daylight Savings Time

The pump time and date do not automatically update for daylight savings time. If you live in a geographical area that follows daylight savings time, you must manually update the time and date. You may choose to change the time while the pump is currently being used by a patient, or wait until the patient is finished with the infusion. However, because delivery of some therapies is affected by a time and/or date change, it is recommended that you wait until the infusion is complete.

NOTE: If you update the time while the pump is in use, the timestamps in the event and delivery logs are not updated to reflect daylight savings time prior to the change. All events record the reported time from when the event actually occurred. For your reference, the event log records the time when it was changed.

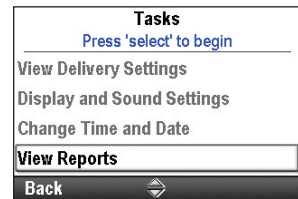
Reports

The reports screen is used to view a variety of reporting and record-keeping functions. The reports available are:

- Total Given (not available for PCA therapies)
- Given and PCA Dose Counters (PCA therapies only)
- PCA Dose Graph (PCA therapies only)
- Delivery History and Pie Chart
- Delivery Log
- Event Log
- Protocol Library Summary
- Device Information

The reports can be viewed at any time, with the pump running or stopped.

In the Tasks menu, press  or  to highlight **View Reports** and press .






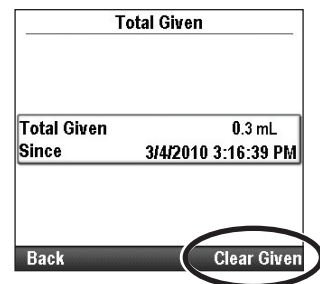
Total Given

NOTE: This report is not available for PCA therapies. See *Given and PCA Dose Counters* on page 96, and *PCA Dose Graph* on page 96 for reports specific to PCA therapies.

The total given report shows the amount of drug given since the date and time indicated, which is the last time the counter was cleared manually, or when a new protocol or new patient was started.

To clear the Given counter:

1. In the View Reports menu, press  or  to highlight **Total Given** and press .
2. Select **Clear Given** to clear the counter. This also updates the time and date on this screen.



Given and PCA Dose Counters

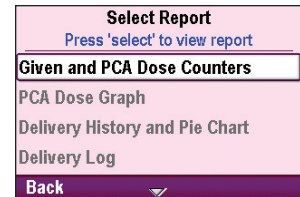
NOTE: This report applies only to PCA therapies. See *Total Given* on page 95 for reports specific to therapies other than PCA.

This report shows the total amount of drug given and the number of PCA doses given and attempted since the date and time indicated, which is the last time the information was cleared manually, or when a new protocol or new patient was started.

- **Total Given** shows the amount of drug (in programmed units) delivered by continuous rate, clinician boluses, and PCA doses.
- **PCA Doses Given** shows the number of PCA doses actually delivered to the patient, including any doses stopped in progress.
- **PCA Doses Attempted** shows the total number of PCA doses attempted by the patient while the pump was running, including those that were delivered, locked out, and stopped in progress.



To clear total given and/or the dose counters:

1. In the View Reports menu, press **▲** or **▼** to highlight **Given and PCA Dose Counters** and press **select**.



2. Press **▲** or **▼** until the counter you wish to clear is highlighted.

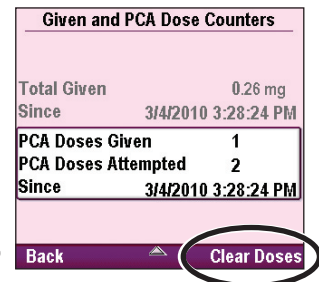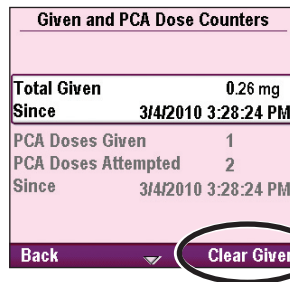3. Select **Clear Given** or **Clear Doses** to clear the counter. This also updates the time and date on this screen.



PCA Dose Graph



NOTE: This report applies only to PCA therapies. For all other therapies, information is available in the report *Delivery History and Pie Chart*.

This graph displays the number of doses attempted and given in 30 minute increments, starting from the current time to 8 hours in the past, or since the start of a new protocol or new patient. Select this report to review the number of attempted doses for a particular time frame.

To view the PCA Dose Graph:

In the View Reports menu, press **▲** or **▼** to highlight **PCA Dose Graph** and press **select**.

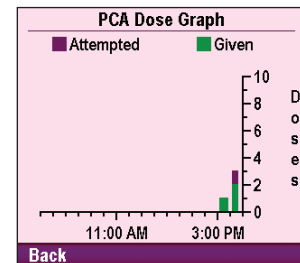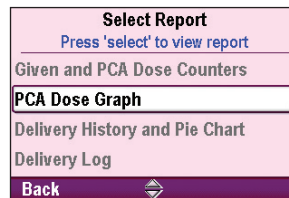In this example, the patient has attempted 3 PCA doses, and 1 was delivered.



Delivery History and Pie Chart

The delivery history and pie chart is a graphical view of the total drug given over a specified time frame, or since the start of a new patient or new protocol. The delivery history is displayed in the units for the current protocol.



The time frame can be adjusted in various intervals from 30 minutes to 7 days. To view the pie chart in different time frames, press ▲ or ▼. This provides a quick review of the methods of delivery over the interval selected.

The information displayed in the pie chart differs for each therapy, as shown in the table below:

PCA	Continuous	Intermittent	Step	Taper
KVO Rate	KVO Rate	KVO Rate	KVO Rate	KVO Rate
Continuous Rate	Continuous Rate	Doses	Step and Plateau	Taper and Plateau
PCA Dose				
Clinician Bolus				

To view the Delivery History and Pie Chart:

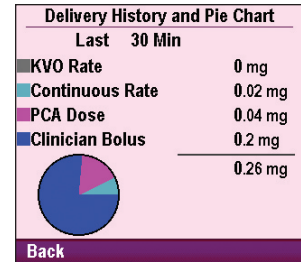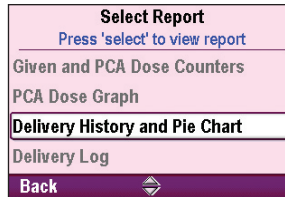In the View Reports menu, press ▲ or ▼ to highlight **Delivery History and Pie Chart** and press **select**.




Delivery Log

The delivery log is a subset of the event log and contains specific information regarding delivery events. The delivery log is maintained by the pump, and displays all entries since the last time a new patient was started. Delivery log information includes but is not limited to:


- Programming changes to the patient-specific parameters
- Changes to infusion status
- Pump started, stopped, powered up, and powered down
- Starting a new protocol

To view the Delivery Log:

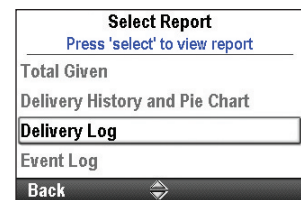1. In the View Reports menu, press ▲ or ▼ to highlight **Delivery Log** and press **select**.



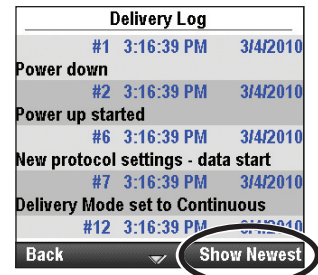2. While viewing the delivery log, you may quickly scroll from the oldest to the newest entries by selecting **Show Oldest** or **Show Newest**.





Event Log

The pump records the time and date of each event, and lists events in order, with the most recent at the bottom of the screen through the last 5000 events. The event log information includes but is not limited to:

- Hourly given totals
- Delivery related events
- Alarms
- Error codes
- Power source changes
- Cassette changes
- Protocol library changes
- Changes to pump programming or settings

To view the Event Log:

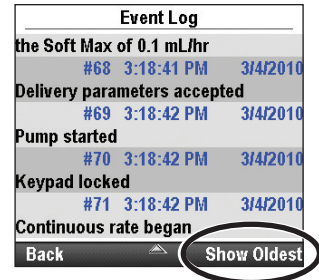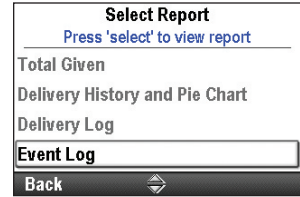- In the View Reports menu, press **▲** or **▼** to highlight **Event Log** and press **(select)**.
- While viewing the event log, you may quickly scroll from the oldest to the newest entries by selecting **Show Oldest** or **Show Newest**.





Protocol Library Summary

The protocol library summary allows you to view information about the protocol library currently installed in the pump. This screen tells you the name of the protocol library, the revision, and the number of protocols in the library.

To view the Protocol Library Summary:

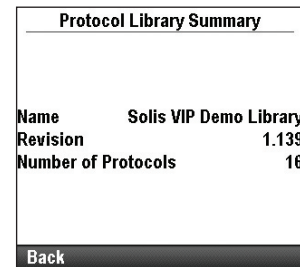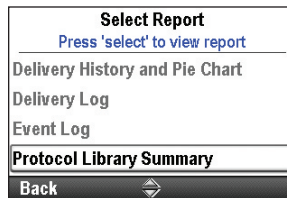In the View Reports menu, press **▲** or **▼** to highlight **Protocol Library Summary** and press **(select)**.



Device Information

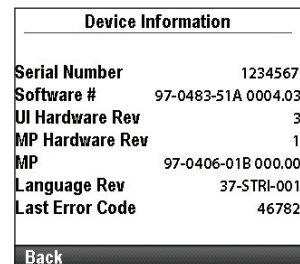
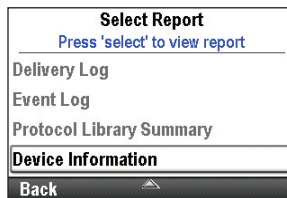
The device information screen allows you to view the information about the pump, including:

- Serial number
- Software version
- User interface (UI) hardware revision
- Motor processor (MP) hardware and software revision
- Language revision
- Last error code (if one exists)

NOTE: Follow your facility's procedure for handling error codes. See *Alarms and Messages* on page 110 for more information.

To view the Device Information:

In the View Reports menu, press **▲** or **▼** to highlight **Device Information** and press **(select)**.



Advanced Tasks Menu Overview

Most pump configurations are set up by a CADD[®]-Solis system administrator. Protocol libraries consisting of therapies, qualifiers, and drugs are created using the CADD[®]-Solis Medication Safety Software. The Advanced Tasks menu allows you to select a standard therapy protocol created by the CADD[®]-Solis administrator, or make changes to an individual protocol for single use, or to manually program a protocol.

NOTE:

- When creating or editing a protocol under Advanced Tasks, it is important to note that the edits made affect only the protocol currently in use.
- The Advanced Tasks screens revert to the home screen if there is a 4 minute delay between button presses. Make certain you are still in the Advanced Tasks menu after any delays.
- Most of the advanced tasks require a security code to unlock the keypad. This step is necessary only if the keypad is not already unlocked.




WARNING:


- **Do not disclose the pump security codes or any other information that would allow the patient or unauthorized clinician complete access to programming and operating functions. Improper programming could result in serious patient injury or death.**
- **Do not leave the pump unattended while unlocked. All programming functions are accessible while the pump is unlocked, and improper programming could result in serious patient injury or death.**

The following functions are on the Advanced Tasks menu.

Advanced Task Menu Function	Description on Page
Give Clinician Bolus (PCA therapy)	43
Step Up (Step therapy)	67
Step Down (Step therapy)	68
Patient Permissions	100
Air and Occlusion Settings	102
Alarm Settings	105
Security Settings	18
Start New Patient	27
Start New Protocol Same Patient	27
Delivery Hard and Soft Limits	33
Reset to Factory Settings	109

To access the Advanced Tasks menu:

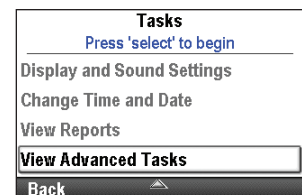1. From the home screen, select **Tasks**.
2. Press  or  to highlight **View Advanced Tasks** and press .




Patient Permissions

Security levels are used to limit unauthorized access to certain programming and operating functions. In some cases, the default security settings may be changed to allow patients additional access to pump functions. Do not reduce security settings without providing proper training to patients.

WARNING:

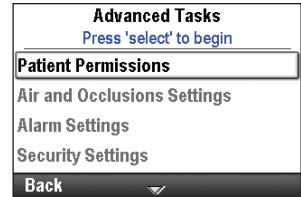- Do not disclose the pump security codes or any other information that would allow the patient or unauthorized clinician complete access to programming and operating functions. Improper programming could result in serious patient injury or death.
- Do not leave the pump unattended while unlocked. All programming functions are accessible while the pump is unlocked, and improper programming could result in serious patient injury or death.






In the Advanced Tasks menu, press  or  to highlight **Patient Permissions** and press .

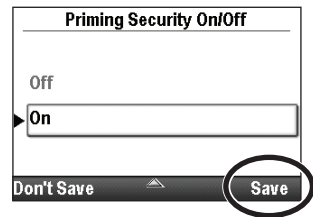
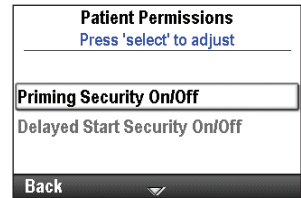


Priming Security On/Off

Setting this value to Off allows patients to prime the tubing without having to enter a security code.

To set patient permissions for priming activity:






1. Stop the pump if it is running.
2. In the Patient Permissions menu, press  or  to highlight **Priming Security On/Off** and press .
3. Unlock the keypad.
4. Press  or  to set the security to **On** (security code required) or **Off** (no security code required) and select **Save**.

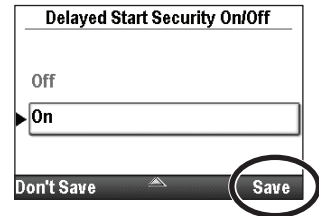
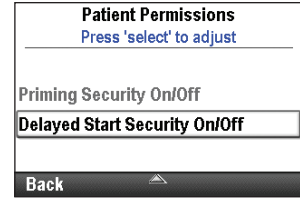


Delayed Start Security On/Off




Setting this value to Off allows patients to set delayed starts (and change next dose start times for Intermittent therapies) without having to enter a security code.


To set patient permissions for delayed start security:

1. Make certain the pump is stopped.
2. In the Patient Permissions menu, press  or  to highlight **Delayed Start Security On/Off** and press .
3. Unlock the keypad.
4. Press  or  to highlight **On** (security code required) or **Off** (no security code required), and select **Save**.



Air and Occlusion Settings

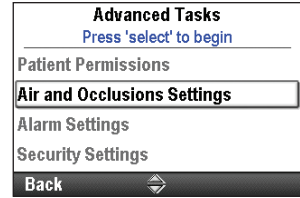In the Advanced Tasks menu, press  or  to highlight **Air and Occlusions Settings** and press .








Air Detector On/Off

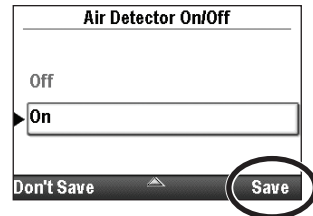
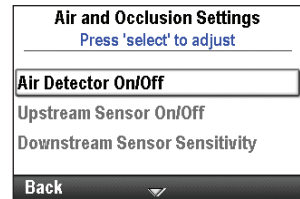
The pump is built with an integrated air detector. The pump may be customized to use the air detector. If air is detected in the part of the tubing that passes through the air detector, an alarm sounds and delivery stops. If an air detector is not required, it may be turned off.

The Air Detector On/Off setting allows you to control whether the air detector is turned on or off. If the air detector is turned on, an alarm sounds when air is detected in the fluid path. For certain therapies, for example, subcutaneous infusions, it may be desirable to turn the air detector off.

WARNING: If the air detector is turned off, the pump does not detect air in the fluid path. It is recommended that you periodically inspect the fluid path and remove any air to prevent air embolism, which could result in serious patient injury or death.

To turn the air detector on or off:

1. Make certain the pump is stopped.
2. In the Air and Occlusions Settings menu, press  or  to highlight **Air Detector On/Off** and press .
3. Unlock the keypad.
4. Press  or  to highlight **On** or **Off**, and select **Save**.








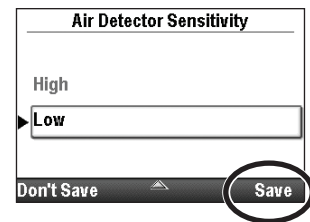
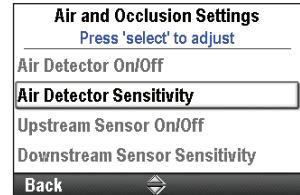
Air Detector Sensitivity

NOTE: If the air detector is set to “Off”, this feature does not appear in the menu.

If the air detector is set to “On”, the air detector sensitivity option allows you to customize how sensitive the pump should be to any air bubbles in the tubing. You may program the pump to have a high or low sensitivity.

To set the air detector sensitivity:

1. Make certain the pump is stopped.
2. In the Air and Occlusions Settings menu, press  or  to highlight **Air Detector Sensitivity** and press .
3. Unlock the keypad.
4. Press  or  to highlight **Low** or **High**, and select **Save**.








Upstream Sensor On/Off

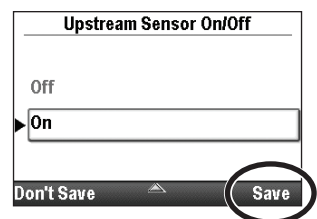
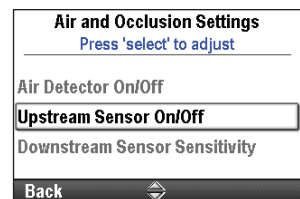
An upstream occlusion is a blockage in the tubing between the pump and the fluid container. The pump contains a sensor that determines if there is an upstream occlusion. This feature can be set to On or Off. If this sensor is set to On and an upstream occlusion is detected, an alarm sounds, delivery stops, and the display shows that there is an upstream occlusion. An upstream occlusion is detected only when administration sets and bags are being used.

NOTE: The pump does not test for an upstream occlusion when using a medication cassette reservoir, even if the upstream occlusion sensor is set to “On.”

WARNING: When the upstream occlusion sensor is turned off, the pump does not detect occlusions upstream between the pump and the reservoir. Periodic inspection of the fluid path for kinks, a closed clamp, or other upstream obstructions is recommended. Upstream occlusions may result in nondelivery of medication. If undetected, occlusions could lead to serious patient injury or death.

To turn the upstream sensor on or off:

1. Make certain the pump is stopped.
2. In the Air and Occlusions Settings menu, press  or  to highlight **Upstream Sensor On/Off** and press .
3. Unlock the keypad.
4. Press  or  to highlight **On** or **Off**, and select **Save**.








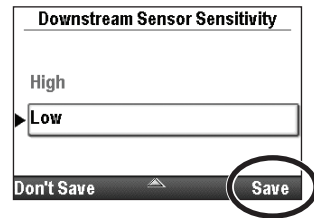
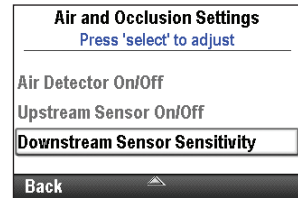
Downstream Sensor Sensitivity

A downstream occlusion is a blockage in the tubing between the pump and the patient. The pump contains a sensor that determines if there is a downstream occlusion. The downstream sensor sensitivity setting allows you to specify a high or low degree of response to pressure buildup in the downstream tubing. When a downstream occlusion is detected, an alarm sounds, delivery stops, and the display shows “High Pressure.”

NOTE: The downstream sensor is always On.

To set the downstream sensor sensitivity:

1. Make certain the pump is stopped.
2. In the Air and Occlusions Settings menu, press  or  to highlight **Downstream Sensor Sensitivity** and press .
3. Unlock the keypad.
4. Press  or  to highlight **Low** or **High**, and select **Save**.



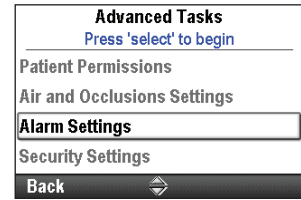
Alarm Settings

You may adjust various alarm settings through the Alarm Settings task in the Advanced Tasks menu. Some of the adjustable settings include whether or not an alarm is active, the type of alarm that sounds, and at what point the alarm sounds. See *Alarms and Messages* on page 110 for more information.

There are 4 types of alarms:

- High Priority
- Medium Priority
- Low Priority
- Informational

In the Advanced Tasks menu, press **▲** or **▼** to highlight **Alarm Settings** and press **select**.





Infusion Alerts On/Off

This alarm setting is available for Step delivery mode only. Turning this setting On results in infusion notification screens 5 minutes before each automatic step up occurs.

To set the infusion alerts:

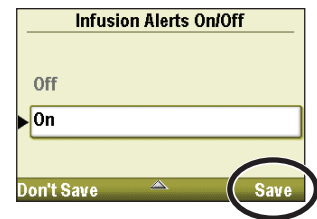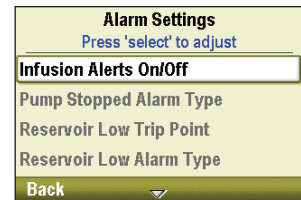1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press **▲** or **▼** to highlight **Infusion Alerts On/Off** and press **select**.
3. Unlock the keypad.
4. Press **▲** or **▼** to highlight **Off** or **On** and select **Save**.





Pump Stopped Alarm Type

The pump stopped alarm allows you to select either a high priority or informational alarm as a reminder that the pump is stopped.

To set the pump stopped alarm type:

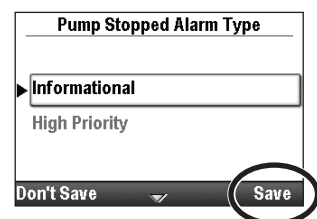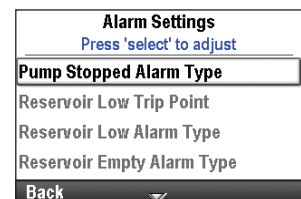1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press **▲** or **▼** to highlight **Pump Stopped Alarm Type** and press **select**.
3. Unlock the keypad.
4. Press **▲** or **▼** to highlight **Informational** or **High Priority**, and select **Save**.






Reservoir Low Trip Point

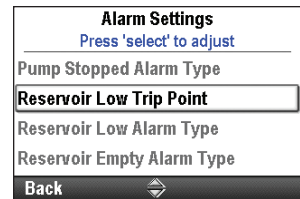
The reservoir low trip point allows you to program the pump to alarm when the reservoir volume reaches a specified level. It can be programmed from 0 to 999 mL in 1 mL increments. If you prefer not to sound an alarm until the reservoir is empty, program the trip point to zero.



NOTE: When programming the reservoir volume trip point, keep the programmed delivery rate of the pump in mind so that sufficient time is available to gather replacement supplies before the reservoir becomes empty.


Once activated, a medium or low priority alarm sounds, depending on the type of alarm selected in the reservoir low alarm type.

To set the reservoir low trip point:

1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press  or  to highlight **Reservoir Low Trip Point** and press .



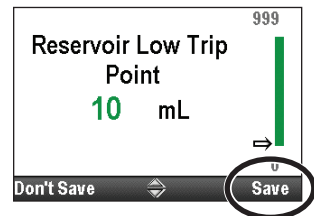3. Unlock the keypad.
4. Press  or  to set the trip point to the desired volume and select **Save**.






Reservoir Low Alarm Type

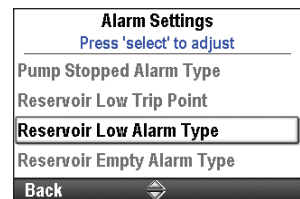
The reservoir low alarm type allows you to select either an “Insistent and One Time Only” alarm (medium priority) or a “Non-Insistent and Repeating” alarm (low priority) to inform the user that the reservoir volume is low.



Insistent and One Time Only: Alarm does not reoccur once it has been acknowledged with a push of the button on the pump.


Non-Insistent and Repeating: Alarm repeats at the 75%, 50%, and 25% marks of the amount determined in the reservoir low trip point.

To set the reservoir low alarm type:

1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press  or  to highlight **Reservoir Low Alarm Type** and press .



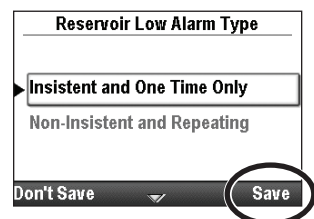3. Unlock the keypad.
4. Press  or  to highlight **Insistent and One Time Only** or **Non-Insistent and Repeating**, and select **Save**.



Reservoir Empty Alarm Type

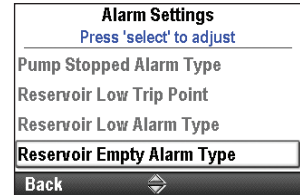
The reservoir empty alarm type option allows you to select between two types of alarms when the reservoir is empty:

One Time Only: When a reservoir empty alarm occurs, you must press a button to acknowledge or silence the alarm. This alarm does not reoccur once it has been acknowledged.

Repeat Until Removed or Reset: When a reservoir empty alarm occurs, you must press a button to acknowledge or silence the alarm. The alarm is repeated every 5 minutes until the reservoir volume has been set to a new value, the reservoir is removed, or the pump is powered down.

To set the reservoir empty alarm type:

1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press **▲** or **▼** to highlight **Reservoir Empty Alarm Type** and press **select**.
3. Unlock the keypad.
4. Press **▲** or **▼** to highlight **One Time Only** or **Repeat Until Removed or Reset**, and select **Save**.





PM Reminder On/Off

If your facility establishes a maintenance program for the pump, you can use the PM reminder on/off option to program the pump to sound an alarm when it is time for preventive maintenance. The message begins appearing on the date programmed and during every power up until the date is reset. Use the PM reminder on/off screen to turn the alarm on and then continue to the PM reminder screen to specify the interval at which the message should appear, or use it to reset the reminder.

To turn the PM reminder on or off:

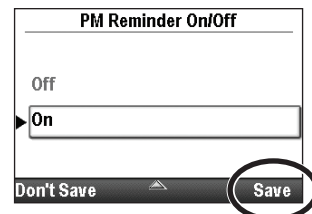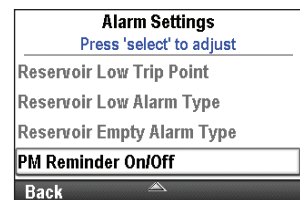1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press **▲** or **▼** to highlight **PM Reminder On/Off** and press **select**.
3. Unlock the keypad.
4. Press **▲** or **▼** to highlight **Off** or **On**, and select **Save**.






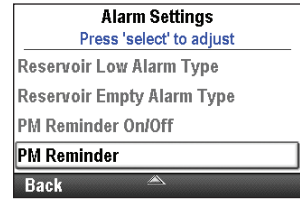
5. If you set the **PM Reminder On/Off** to On, press **▲** or **▼** to highlight **PM Reminder** and press **select**. Instructions for scheduling the PM Reminder are in the next section.

PM Reminder

The PM Reminder appears as an option only if the PM Reminder On/Off setting is set to “On”. Use the PM Reminder to schedule the amount of time that should pass before a preventive maintenance reminder displays (1 to 24 months).

To schedule the PM interval immediately after turning the PM Reminder On/Off setting to On:




1. In the Alarm Settings menu, press  or  to highlight **PM Reminder** and press .

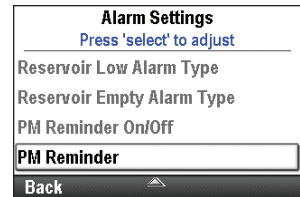





2. **PM Interval** is highlighted (see step 4 below to edit the settings).


PM Interval

To edit the PM interval:

1. Make certain the pump is stopped.
2. In the Alarm Settings menu, press  or  to highlight **PM Reminder** and press .



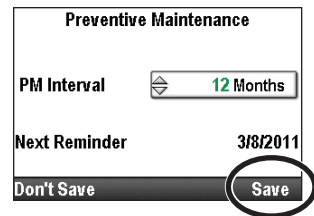3. Unlock the keypad, if necessary.
4. **PM Interval** is highlighted. Press  to edit it. Press  or  to choose the number of months that should pass before the next preventive maintenance reminder. The Next Reminder date displayed changes along with the PM Interval to reflect the actual date of the next reminder. Select **Save**.






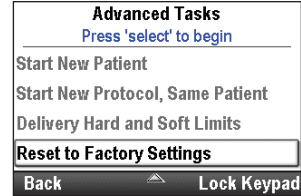
Reset to Factory Settings

Reset to Factory Settings allows you to erase the protocol library and pump settings, and restore the pump to the factory default settings.

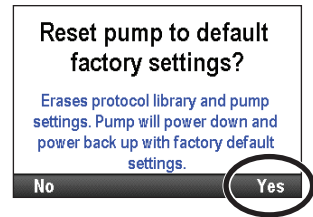
NOTE: When the pump is reset to factory settings, all protocol library information is lost.

To reset the pump to factory settings:

1. Make sure the pump is stopped.
2. In the Advanced Tasks menu, press  or  to highlight **Reset to Factory Settings** and press .



3. Unlock the keypad.
4. A screen appears that explains that you are about to erase the protocol library and pump settings. Select **Yes**.



5. The pump powers down and restarts in the continuous mode without a protocol library, and in the default factory settings.

References and Troubleshooting

Alarms and Messages

The pump can sound multiple alarms. For many of the alarms, you have the option to “acknowledge” or “silence.”

- **Acknowledge**—the alarm clears from the screen.
- **Silence**—the alarm stays on the screen, but is silenced for 2 minutes before it sounds again. The alarm continues until it is acknowledged or resolved.

The alarms may have different sounds depending on the sound theme selected. There are 3 different sound themes for the alarms and beeps that the pump makes: standard, intense, and distinctive. See *Sound Theme* on page 90 for more information on previewing and selecting the sound themes.

Types of Alarms

System Fault Alarm

If this screen appears, an unrecoverable error may have occurred, such as a hardware or software fault. The amber indicator light is continuously illuminated during these conditions and is accompanied by an audible two-tone alarm. If a system fault occurs, the fault should be reported to Smiths Medical Customer Service at 1 800.258.5361 (USA) or Smiths Medical International Ltd. at +44 (0) 1923 246434.

To clear this alarm, remove power from the pump by opening the battery door, and if necessary, removing the AC power. Close the battery door and turn the pump back on. If the error code does not repeat, Customer Service may suggest continued use of the pump. If the error is persistent, the pump must be returned for service.

NOTE: Document the error numbers displayed on the system fault screen to help Customer Service identify the problem.

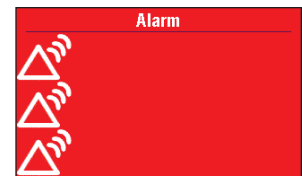


System fault alarm (red).
The numbers relate to the error that caused the system fault.

CAUTION: If the power up results in an error message indicating that the protocol library was lost, do not proceed with using the pump. Follow your facility’s procedures for downloading protocol libraries

High Priority Alarm


A high priority alarm always pauses or stops the pump if it is running. The alarm is accompanied by a **red** screen, and it persists until acknowledged by the press of a key on the pump or until the condition that triggered the alarm goes away (for example, high pressure). The alarm can be silenced with a key press and will sound again after 2 minutes if the alarm condition still exists.



High priority alarm (red)

Medium Priority Alarm

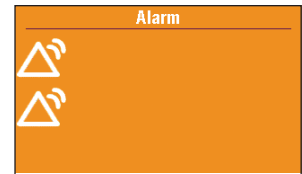
A medium priority alarm does not stop the pump if it is running. The alarm is accompanied by an **amber** screen, and it persists until acknowledged by the press of a key on the pump or until the condition that triggered the alarm goes away. The alarm can be silenced with a key press and will sound again after 2 minutes if the alarm condition still exists.



Medium priority alarm (amber)

Low Priority Alarm

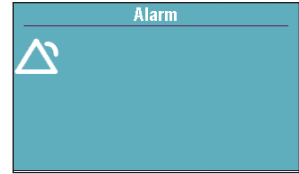A low priority alarm does not stop the pump if it is running. The alarm is accompanied by a **blue** screen, and it persists for 5 seconds unless it is acknowledged by the press of a key or the condition that triggered the alarm goes away before the 5 seconds have passed. (Some low priority alarms persist for longer than 5 seconds.)



Low priority alarm (blue)

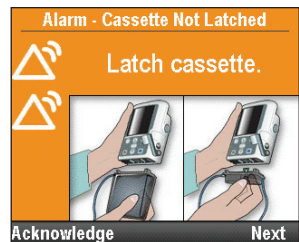
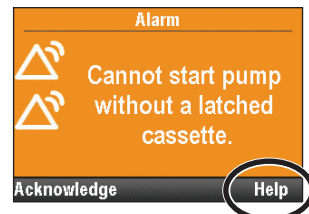
Informational Message

An informational priority message does not stop the pump if it is running. This message appears in the status bar, and does not display a new alarm screen. It persists for 5 seconds and may be silent, requiring no acknowledgement. Some informational examples are “Cassette Locked,” and “Cassette Unlocked”.

Alarm Help Screens

Additional information may be displayed when certain alarms occur. The help screens describe what you can do to try to solve the current problem that is causing the alarm.

1. When an alarm occurs, select **Silence** to quiet the alarm.
2. If help screens are available for the alarm, “Help” appears above the right soft key. To view the help screens, select **Help**.
3. Follow the applicable instructions provided on the help screen. To page through all available help screens, press **Next** repeatedly. Press **Acknowledge** at any time to exit Help.
4. “Retry Help” appears when no additional help steps are available. To page through the help screens again, press **Retry Help**. The Alarm screen reappears as in Step 2 above. Repeat Steps 2 and 3. When the alarm clears, the help screens will no longer be displayed.



Troubleshooting

Issue	Possible Solution
A continuous two-tone alarm is sounding, and the amber light is lit or flashing.	Delivery has stopped. Read the message on the display and refer to the list of messages in the table below. If the display is blank or contains random characters, the 4 AA batteries or the rechargeable battery pack may be depleted. Install 4 new AA batteries or a rechargeable battery pack.
The pump is sounding 2 beeps every 2 seconds, and the amber light is flashing.	Look at the message on the display and refer to the list of messages in the table below.
Three beeps sound every 5 minutes.	This is a reminder that the pump is stopped.
After installing 4 new AA batteries and powering up the pump, no screen appears and no beep sounds.	The batteries may be installed incorrectly. Review the procedure for installing batteries. Be sure to match the polarity (+ and -) markings inside the battery door with the markings on the batteries. If there is still no power, the batteries may be completely depleted. Install 4 new AA batteries.


Alarms and Messages, Alphabetical List


Alarm / Message	Alarm Priority	Description / Corrective Action
(Screen is blank and alarm is sounding)	High	The pump was delivering and the batteries were removed or the battery door was opened. The pump lost power and is no longer delivering. Clear the alarm by turning the pump back on, or the alarm will stop after the power has been off for a minimum of 2 minutes.
A setting was edited, but not saved, and the edit was lost.	Medium	A parameter was being manually edited, but it was not saved and the pump reverted to the home screen. Select Acknowledge to clear the alarm and if appropriate, edit the parameter and save.
AC adapter disconnected.	Low	The AC adapter was disconnected and the pump is being powered by the 4 AA batteries or the rechargeable battery pack. Select Acknowledge to clear the alarm, or the alarm will automatically clear after 5 seconds. If desired, reconnect the AC adapter.
Air in-line detected. Press 'Acknowledge' then prime tubing.	High	The air detector has detected air in the fluid path. The fluid path may contain air bubbles. The pump was delivering and is now stopped and will not run. Select Acknowledge to clear the alarm. If the fluid path contains air bubbles, close the clamps and disconnect the fluid path from the patient. Follow the instructions for removing air by priming (see page 83 for more information on priming). Restart the pump.
Battery depleted. Pump stopped.	High	The pump was delivering but is now stopped and the battery power is too low to operate the pump. If the AC adapter is attached, select Acknowledge to clear the alarm. Remove the batteries. Install 4 new AA batteries or a rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, restart the pump.
Battery depleted. Pump will not run.	Medium	The battery power is too low to operate the pump. If the AC adapter is attached, select Acknowledge to clear the alarm. Remove the batteries. Install 4 new AA batteries or a rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, restart the pump.

Alarm / Message	Alarm Priority	Description / Corrective Action
Battery low. Replace battery.	Low	The rechargeable battery pack or the 4 AA batteries are low but the pump is still operable. Select Acknowledge to clear the alarm, or the alarm will automatically clear after 5 seconds. Recharge or change the rechargeable battery pack or replace the 4 AA batteries soon.
Battery removed. Pump stopped.	High	The rechargeable battery pack or the 4 AA batteries were removed while the pump was running. The pump was delivering and is now stopped. Select Acknowledge to clear the alarm. Install 4 new AA batteries or a rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, restart the pump.
Battery removed. Pump will not run.	Medium	The pump is stopped and the rechargeable battery pack or the 4 AA batteries were removed, but the pump is still powered by the AC adapter. Select Acknowledge to clear the alarm. Install a fully charged rechargeable battery pack or 4 new AA batteries. To start delivery, good batteries must always be installed, even when an external source of power is connected.
Cannot start pump. Rechargeable battery has reached end of use.	Medium	The rechargeable battery pack is at the end of its life. It has been discharged and recharged so many times that it is no longer able to hold a good charge. Remove the battery from service. To start delivery, good batteries must always be installed even when an external source of power is connected. Select Acknowledge to clear the alarm. Install 4 new AA batteries or a new rechargeable battery pack. If appropriate, start the pump.
Cannot start pump with a depleted battery.	Medium	The battery power is too low to operate the pump. To start delivery, good batteries must always be installed even when an external source of power is connected. Select Acknowledge to clear the alarm. Install 4 new AA batteries or a fully charged rechargeable battery pack. If appropriate, start the pump.
Cannot start pump with a reservoir volume of zero.	Medium	The reservoir volume in the pump is set to zero. Select Acknowledge to clear the alarm. Edit or reset the reservoir volume to the correct value. If appropriate, start the pump. See page 86 for more information about resetting the reservoir volume.
Cannot start pump with air in-line. Prime tubing.	Medium	The air detector has detected air in the fluid path directly under the air detector. The fluid path may contain air bubbles. Select Acknowledge to clear the alarm, then: If the fluid path contains air bubbles, close the clamps and disconnect the fluid path from the patient. Follow the instructions for removing air by priming (see page 83 for more information). If appropriate, start the pump.
Cannot start pump with an unusable battery.	Medium	The batteries installed are either the wrong kind of AA batteries, or the rechargeable battery pack is not compatible with the pump. Select Acknowledge to clear the alarm. Remove the batteries. Install 4 new AA batteries or a rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, start the pump.
Cannot start pump without a battery.	Medium	The pump does not have any batteries installed, although it is still powered by the AC adapter. Select Acknowledge to clear the alarm. Install 4 new AA batteries. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, start the pump.

Alarm / Message	Alarm Priority	Description / Corrective Action
Cannot start pump without a latched and locked cassette.	Medium	The pump will not start without a cassette attached and locked. Select Acknowledge to clear the alarm. Make sure a cassette is properly attached and locked, then start the pump. NOTE: This alarm has associated help screens.
Cannot start pump without a latched cassette.	Medium	The pump will not start without a cassette attached. Select Acknowledge to clear the alarm. Make sure a cassette is properly attached, then start the pump. NOTE: This alarm has associated help screens.
Cassette detached. Pump stopped.	High	The cassette was detached while the pump was delivering and the pump is now stopped. Select Acknowledge to clear the alarm. Reattach the cassette and if appropriate, restart the pump.
Cassette locked, but not latched. Unlock and reattach cassette.	High	The cassette/keypad lock is locked, but there is no cassette attached. The pump is stopped and will not run. Select Acknowledge to clear the alarm. Unlock the cassette/keypad lock and reattach the cassette.
Cassette not attached properly. Reattach cassette.	High	The cassette is not properly attached. Close the tubing, remove the cassette, then reattach. If the alarm persists, replace the cassette. NOTE: You must remove the cassette to continue.
Cassette unlocked. Lock cassette.	Medium	This is a reminder that the cassette is not locked while the pump is delivering. Lock the cassette to clear the alarm.
Cassette was partially unlatched. Fully remove and reattach the cassette.	Medium	The cassette was not completely removed from the pump before it was reattached and the pump's sensors are not able to detect the cassette type. Remove the cassette and reattach it, then verify the cassette type in the pump display. If the alarm persists, replace the cassette. NOTE: You must remove the cassette to continue.
Check for empty tubing or reservoir. Pump stopped.	High	The tubing beneath the pump may not contain fluid, or the fluid container may be empty. The pump is stopped and will not run. Select Acknowledge to clear the alarm. Check to see if the fluid container is empty. If fluid is present in the reservoir, clamp the tubing, remove the cassette, and check for air in the tubing. If the alarm persists, remove the pump from service and contact Customer Service to return the pump for service.
Delayed start time adjusted to [date] [time].	Medium	A delayed start is programmed and a change to the time or date caused an adjustment to the start time in order to keep the interval unchanged from the current time and date.
Delivery limit reached. Partial PCA dose delivered. Running at KVO rate.	Low	The programmed delivery limit has been reached, a portion of the PCA dose was delivered, and the pump is delivering fluid at the KVO rate. This alarm occurs when the continuous rate is programmed to greater than 0 mL/hr, and either a PCA dose or the continuous rate has caused the delivery limit to be exceeded. Select Acknowledge to clear the alarm or the alarm will automatically clear after 5 seconds.
Delivery limit reached. Running at KVO rate.	Low	The programmed delivery limit has been reached, and the pump is delivering fluid at the KVO rate. This alarm occurs when the continuous rate is programmed to greater than 0 mL/hr, and either a PCA dose or the continuous rate has caused the delivery limit to be exceeded. Select Acknowledge to clear the alarm or the alarm will automatically clear after 5 seconds.

Alarm / Message	Alarm Priority	Description / Corrective Action
Delivery too slow. [Quantity] of fluid delivery delayed..	Medium	The pump is busy with too many activities and does not have sufficient resources to support the programmed delivery rate. Delivery has fallen behind by the quantity indicated. Select Acknowledge to clear the alarm. If this alarm occurs regularly, remove the pump from use and contact Customer Service to return the pump for service.
Depleted battery is charging.	Low	The rechargeable battery pack is depleted and is being recharged with the AC adapter. Select Acknowledge to clear the alarm or the alarm will automatically clear after 5 seconds.
Dose is now overdue. Pump is stopped.	Medium	The pump is stopped and a dose is overdue for its scheduled delivery. Select Acknowledge to clear the alarm, and then start the pump.
Dose scheduled to start in [time]. Pump is stopped.	Low	The pump is stopped and a dose is scheduled for delivery in the time indicated. Select Acknowledge to clear the alarm, and then start the pump.
Downstream occlusion. Clear occlusion between pump and patient.	High	The pump has detected high pressure, which may be resulting from a downstream blockage, a kink in the fluid path, or a closed tubing clamp. Delivery is paused and will resume if the occlusion is removed. Remove the obstruction to resume operation. Or select Stop Pump to stop the pump and silence the alarm for 2 minutes, then remove the obstruction and restart the pump. NOTE: To reduce the potential for a bolus delivery after an occlusion, perform the following: <ol style="list-style-type: none"> 1. Select Stop Pump to stop the pump. 2. Close the distal clamp. If the distal clamp is the cause of obstruction, keep it closed and continue with step 4. 3. Remove the obstruction. 4. Detach the CADD™ medication cassette reservoir or CADD® administration set from the pump. 5. Open the Flow Stop feature, if present. 6. Wait 10 seconds. 7. Close the Flow Stop feature, if present. 8. Reattach the CADD™ medication cassette reservoir or CADD® administration set to the pump. 9. Open the distal clamp. 10. Review the pump program. 11. Restart the pump. NOTE: This alarm has associated help screens.
External power source faulty. Change power source.	Medium	The AC adapter output voltage is too high. Select Acknowledge to clear the alarm. The AC adapter is faulty, remove from service.
High flow admin set required. Remove cassette.	High	The delivery-specific parameters are programmed to values that cause the maximum rate of delivery to exceed 250 mL/hr, which requires a high flow administration set. Replace the standard volume cassette with a high volume administration set to continue.
High volume admin set not allowed. Remove cassette.	High	The CADD® high volume administration set cannot be used. The pump is stopped and will not run. Replace the high volume administration set with a standard volume cassette to continue.
Infusion is now overdue. Pump is stopped.	Medium	The pump is stopped and the infusion is overdue for its scheduled delivery. Select Acknowledge to clear the alarm, and then start the pump.

Alarm / Message	Alarm Priority	Description / Corrective Action
Infusion scheduled to start in [time]. Pump is stopped.	Low	The pump is stopped and an infusion is scheduled for delivery in the time indicated. Select Acknowledge to clear the alarm, and then start the pump.
Key stuck. Release key or remove power. Pump stopped.	High	A key may be pressed down. The pump is stopped and will not run. Make sure there is nothing pressing on any of the keys. If the alarm persists, close the tubing clamp, remove the batteries to turn off the pump, and remove the pump from use. Contact Customer Service to return the pump for service.
Lock cassette to start pump.	Medium	PCA mode only. The cassette must be locked onto the pump before beginning delivery. Lock the cassette to clear the alarm and the pump will automatically start. NOTE: This alarm has associated help screens.
Loss of power occurred while pump was running.	Medium	The pump lost power while it was running. This alarm occurs when the pump restarts. Select Acknowledge to clear the alarm. To continue the infusion, restart the pump.
Motor service due.	Medium	The pump motor requires service. Select Acknowledge to clear the alarm. Remove the pump from use at the next cassette change and contact Customer Service to return the pump for service.
New pump settings downloaded. Press 'Acknowledge' then review.	Low	A new protocol was just sent to the pump from the CADD®-Solis Medication Safety Software. Select Acknowledge to clear the alarm. Review the protocol to ensure that it is correct.
Next dose start time adjusted to [date] [time].	Medium	A change to the time or date caused an adjustment to the next dose start time in order to keep the interval unchanged from the current date and time.
PCA dose cord button stuck. Release or remove cord.	High	The PCA remote dose cord button may be pressed down. The pump is stopped and will not run. Make sure nothing is pressing on the PCA remote dose cord button. If the alarm persists, remove the PCA remote dose cord to clear the alarm and contact Customer Service. Use the PCA dose  button on the pump, or use another PCA remote dose cord.
PCA dose cord disconnected.	Medium Low	Medium: The PCA remote dose cord was disconnected from the pump while the pump was delivering. Select Acknowledge to clear the alarm or reattach the PCA remote dose cord. Low: The PCA remote dose cord was disconnected from the pump while the pump was stopped. Select Acknowledge to clear the alarm or the alarm will automatically clear after 5 seconds.
Preventive maintenance due.	Medium	Your facility established a maintenance program for the pump and the pump is due for preventive maintenance. Select Acknowledge to clear the alarm and refer to your facility's policy for preventive maintenance.
Protocol library updating. Reselect protocol when update complete.	Medium	A new or updated protocol library is currently being sent to the pump. The pump will not allow the selection of any new protocols while this update is in process. Select Acknowledge to clear the alarm, or it will automatically clear when the update is complete.
Pump does not have a protocol library.	Medium	The pump had a protocol library the last time it was powered on, but now it does not. This may happen if the pump was manually reverted to the factory default, recently had a software update, or if an attempt to install a protocol library failed. Select Acknowledge to clear the alarm and contact the CADD®-Solis system administrator to download a new protocol library.

Alarm / Message	Alarm Priority	Description / Corrective Action
Pump settings and patient data lost.	Medium	The pump reverted to the factory default. The pump was either manually reverted to this default, recently had a software update, or has not been in use for some time. Select Acknowledge to clear the alarm and contact the CADD®-Solis system administrator to reprogram the pump.
Pump stopped by an alarm that has since cleared.	High	The pump was stopped by another high priority alarm. That alarm was not acknowledged, but the problem has since cleared. Select Acknowledge to clear the alarm and if appropriate, restart the pump. The Event Log recorded the alarm that stopped the pump. For information on accessing the Event Log, see page 98.
Pump stopped reminder. When ready, press the stop/start  key.	High	The pump was stopped and is not delivering. Select Acknowledge to clear the alarm. Start the pump, if appropriate. The alarm will repeat in 5 minutes if the pump is not restarted or powered down.
Rechargeable battery near end of use. Replace battery.	Medium	The rechargeable battery pack is near the end of its life. It has been discharged and recharged so many times that it will soon be at the end of its use. You may continue to use it in this state.
Rechargeable battery reached end of use. Pump will not run.	High	The rechargeable battery pack is at the end of its life. It has been discharged and recharged so many times that it is no longer able to hold a good charge. Remove the battery from service. Install 4 new AA batteries or a fully charged rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, restart the pump.
Remove and reattach cassette.	High	The pump has detected a damaged cassette. Close the tubing clamp and inspect the cassette for damage. The pump is stopped and will not run. Replace cassette if necessary. NOTE: You must remove the cassette to continue. NOTE: This alarm also occurs if a cassette is attached during one of the following situations: <ul style="list-style-type: none"> • The pump is set to the factory default, and then powered off, and then on again. • The pump was loaded with new software, and then powered off, and then on again.
Reservoir volume is zero. Pump stopped.	High	The reservoir volume is 0.0 mL. The pump was delivering, but is now stopped and will not run. Select Acknowledge to clear the alarm. Attach a new reservoir, and reset or edit the value of the reservoir volume, if appropriate.
Reservoir volume low.	Medium Low	Medium: The programmed reservoir volume trip point has been reached, indicating the level of fluid in the reservoir is low. Select Acknowledge to clear the alarm. Low: The programmed reservoir volume trip point has been reached, indicating the level of fluid in the reservoir is low. Select Acknowledge to clear the alarm, or the alarm will automatically clear after 5 seconds. Prepare to install a new reservoir, and reset or edit the value of the reservoir volume, if appropriate.
Unknown cassette type. Remove cassette.	High	A cassette that is incompatible with the pump is detected. The pump is stopped and will not run. Close the tubing clamp, remove and then reattach the cassette. If the alarm persists, replace the cassette. NOTE: You must remove the cassette to continue.

Alarm / Message	Alarm Priority	Description / Corrective Action
Unusable battery. Pump stopped.	High	<p>The batteries installed are either the wrong type of AA batteries, or the rechargeable battery pack is not compatible with the pump. Select Acknowledge to clear the alarm.</p> <p>Remove the battery or batteries from service. Install 4 new AA batteries or a fully charged rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, start the pump</p>
Unusable battery. Pump will not run.	Medium	<p>The batteries installed are either the wrong type of AA batteries, or the rechargeable battery pack is not compatible with the pump. Select Acknowledge to clear the alarm.</p> <p>Remove the battery or batteries from service. Install 4 new AA batteries or a fully charged rechargeable battery pack. To start delivery, good batteries must always be installed even when an external source of power is connected. If appropriate, start the pump.</p>
Upstream occlusion. Clear occlusion between pump and reservoir.	High	<p>Fluid is not flowing from the fluid container to the pump, which may be due to a kink, closed clamp, or air bubble in the tubing between the fluid container and pump. Delivery is paused and will resume if the occlusion is removed. Remove the obstruction to resume operation. The alarm will clear when the occlusion is removed. You must acknowledge this alarm after it clears if it has occurred and cleared more than 3 times within 15 minutes.</p> <p>NOTE: This alarm has associated help screens.</p>

Cleaning the Pump and Accessories

NOTE: Refer to the instructions for use for each accessory before proceeding with cleaning.

The following solutions may be used to clean the pump and accessories, unless otherwise specified:

- Soap solution
- Benzalkonium chloride concentrate (0.13%)
- Glutaral concentrate, USP (2%)
- 10% solution of household bleach (1 part household bleach to 9 parts water)
- Alcohol, USP (93%)
- Isopropyl alcohol, USP (99%)
- Chlorhexidine gluconate (4%)
- PDI® Super Sani-Cloth®
- Madacide (MADA Medical)
- Virex® II (Johnson Wax)
- Coverage® Spray and Coverage® HB Plus (Steris)
- CaviCide® (Metrex)
- Quik Fill Compac A-456-N (Airkem)

CAUTION:

- Do not immerse the pump in cleaning fluid or water. Do not allow solution to soak in to the pump, accumulate on the keypad, or enter the battery compartment, USB port, remote dose cord jack, or power jack areas. Moisture buildup inside the pump may damage the pump.
- Do not clean the pump with acetone, other plastic solvents, or abrasive cleaners, as damage to the pump may occur.

-
1. Ensure that debris is not allowed to build up on the pressure plate surface of the pumping mechanism. Consult the Technical Manual for a diagram of the pumping mechanism. Inspect the air detector sensor slot and remove any debris.
 2. Dampen a soft, lint-free cloth with cleaning solution. Apply the solution to the exterior surface of the pump or accessory according to manufacturer's instructions. ***Do not allow the solution to soak into the pump or accessory.***
 3. Allow the pump to dry completely before use.

Radiation and Magnetic Resonance Imaging (MRI)

CAUTION:

- The pump should *not* be directly irradiated by therapeutic levels of ionizing radiation due to the risk of permanent damage to the electronic circuitry. The best procedure to follow is to remove the pump from the patient during therapeutic radiation sessions or diagnostic levels of radiographic and fluoroscopic radiation. If the pump must remain in the vicinity during a diagnostic or therapy session, it should be shielded, and its ability to function properly should be confirmed following treatment.
 - Do not expose the pump directly to ultrasound, as permanent damage to the electronic circuitry may occur.
 - Magnetic fields produced by MRI equipment may adversely affect the operation of the pump. Remove the pump from the patient during MRI procedures and keep it at a safe distance from magnetic energy.
 - Use of this pump on patients monitored by electronic equipment may cause artifactual interference. As with all electronic equipment, electrical artifacts which affect the performance of other equipment, such as ECG monitors, can occur. The user should check the correct function of the equipment prior to use.
 - Do not use the pump in hyperbaric chambers as they affect how the pump works and may also cause damage to the pump.
-

Standards used in development of the pump

The following standards were used in whole or part in the development of the pump.

Medical Electrical Equipment

EN 60601-1 (1990), Medical Electrical Equipment, Part I: General Requirements for Safety. Amendment A1 (1993) Amendment A13 (1996) Amendment A2 (1995).

EN 60601-2-24 (1998), Medical Electrical Equipment, Part 2-24: particular Requirements for Safety of Infusion Pumps and Controllers.

EN 60601-1-4 (1996), Medical Electrical Equipment, Part 1-4: General Requirements for Safety—Collateral Standard: Programmable electrical medical systems. Amendment A1: 1999.

IEC 60601-1 (2nd Edition, 1988), Medical Electrical Equipment, Part 1: General Requirements for Safety. Amendment 1 (1991) Amendment 2 (1995).

IEC 60601-1-4 (2000), Medical Electrical Equipment, Part 104: General Requirements for Safety—Collateral Standard: Programmable electrical medical systems.

IEC 60601-2-24 (1998), Medical Electrical Equipment, Part 2-24: Particular Requirements for Safety of Infusion Pumps and Controllers.

CAN/CSA-C22.2 601.1-M90, Medical Electrical Equipment, Part 1: General Requirements for Safety - November 1990 (Canadian Deviations to IEC 60601-1) Update No. 2 (November 2003).

Electromagnetic Compatibility

RTCA/DO -160E (2004), Radiated Emissions Only, Category M Limit.

EN 60601-1-2 (2001), Medical Electrical Equipment, Part 1-2: General Requirements for Safety—Collateral Standard: Electromagnetic Compatibility—Requirements and Tests.

IEC 60601-1-2 (Edition 2.1, 2004-11), Medical Electrical Equipment, Part 1-2: General Requirements for Safety—Collateral Standard: Electromagnetic Compatibility—Requirements and Tests.

IEC 61000-4-2 (2001), Electromagnetic Compatibility (EMC), Part 4-2: Testing and measurement techniques. Electrostatic Discharge immunity test.

IEC 61000-4-3 (2006), Electromagnetic Compatibility (EMC), Part 4-3: Testing and measurement techniques. Radiated, radio frequency, electromagnetic field immunity test.

IEC 61000-4-4 (2004), Electromagnetic Compatibility (EMC), Part 4-4: Testing and measurement techniques. Electrical fast transient/burst immunity test.

IEC 61000-4-5 (2005), Electromagnetic Compatibility (EMC), Part 4-5: Testing and measurement techniques. Surge immunity test.

IEC 61000-4-6 (2004), Electromagnetic Compatibility (EMC), Part 4-6: Testing and measurement techniques. Immunity to conducted disturbances, induced by radio-frequency fields.

IEC 61000-4-8 (2001), Electromagnetic Compatibility (EMC), Part 4-8: Testing and measurement techniques. Power frequency magnetic field immunity test.

IEC 61000-4-11 (2004), Electromagnetic Compatibility (EMC), Part 4-11: Testing and measurement techniques. Voltage dips, short interruptions and voltage variations immunity test.

CISPR11 (2004), Limits and methods of measurement of electromagnetic disturbance characteristics of industrial, scientific and medical (ISM) radio frequency equipment. Amendment 1 (1999) Amendment 2 (2002).

EN 45502-1 (1998), Active implantable medical devices. Part 1. General requirements for safety, marking and information to be provided by the manufacturer.

EN 55011 (1998), Industrial, scientific and medical (ISM) radio frequency equipment—Radio disturbance characteristics—Limits and methods of measurement. Amendment 1 (1999), Amendment 2 (2002)

Miscellaneous Standards

USB 1.1, Universal Serial Bus (USB) Specification, Revision 1.1, September 23, 1998—USB.org

EN 1041 (1998), Information supplied by the manufacturer with medical devices.

IEC/TR 60878 (2003), Graphical symbols for electrical equipment in medical practice.

EN 980 (2003), Graphical symbols for use in the labeling of medical devices.

IEC 60529 (2001), Degrees of protection provided by enclosures (IP Code).

PCA Delivery Mode Scroll Ranges

PCA Continuous Rate Scroll Ranges

Units	Starting Value	Increment	Maximum
Milliliters	0	0.10	100.00
Milligrams only	10% of concentration	Values between 0.01 and 0.5: 0.01	Concentration x 100
Micrograms only	10% of concentration	Values between 0.1 and 0.5: 0.1	Concentration x 100
Milligrams and Micrograms	10% of concentration	Values between 0.5 and 100: 0.1 Values between 100 and 1000: 1.0 Values greater than 1000: 10.0	Concentration x 100

PCA Dose and Clinician Bolus Scroll Ranges: Milliliters

Milliliters	
Increment	Max.
0.05	50

PCA Dose and Clinician Bolus Scroll Ranges: Milligrams

Concentration mg/mL	Milligrams	
	Increment	Max.
0.1	0.01	5
0.2	0.02	10
0.3	0.03	15
0.4	0.04	20
0.5	0.05	25
1	0.05	50
2	0.10	100
3	0.15	150
4	0.20	200
5	0.25	250
6	0.30	300
7	0.35	350
8	0.40	400
9	0.45	450
10	0.50	500
11	0.55	550
12	0.60	600
13	0.65	650
14	0.70	700
15	0.75	750

Concentration mg/mL	Milligrams	
	Increment	Max.
20	1.00	1000
25	1.25	1250
30	1.50	1500
35	1.75	1750
40	2.00	2000
45	2.25	2250
50	2.50	2500
55	2.75	2750
60	3.00	3000
65	3.25	3250
70	3.50	3500
75	3.75	3750
80	4.00	4000
85	4.25	4250
90	4.50	4500
95	4.75	4750
100	5.00	5000

PCA Dose and Clinician Bolus Scroll Ranges: Micrograms

Concentration mcg/mL	Micrograms	
	Increment	Max.
1	0.05	50
2	0.10	100
3	0.15	150
4	0.20	200
5	0.25	250
6	0.30	300
7	0.35	350
8	0.40	400
9	0.45	450
10	0.50	500
11	0.55	550
12	0.60	600
13	0.65	650
14	0.70	700
15	0.75	750
20	1.00	1000
25	1.25	1250
30	1.50	1500

Concentration mcg/mL	Micrograms	
	Increment	Max.
35	1.75	1750
40	2.00	2000
45	2.25	2250
50	2.50	2500
55	2.75	2750
60	3.00	3000
65	3.25	3250
70	3.50	3500
75	3.75	3750
80	4.00	4000
85	4.25	4250
90	4.50	4500
95	4.75	4750
100	5.00	5000
200	10.00	10,000
300	15.00	15,000
400	20.00	20,000
500	25.00	25,000

Military Time

Military Time Conversion

12-Hour Time		24-Hour Time
12:00 AM	MIDNIGHT	00:00
1:00 AM		01:00
2:00 AM		02:00
3:00 AM		03:00
4:00 AM		04:00
5:00 AM		05:00
6:00 AM		06:00
7:00 AM		07:00
8:00 AM		08:00
9:00 AM		09:00
10:00 AM		10:00
11:00 AM		11:00

12-Hour Time		24-Hour Time
12:00 PM	NOON	12:00
1:00 PM		13:00
2:00 PM		14:00
3:00 PM		15:00
4:00 PM		16:00
5:00 PM		17:00
6:00 PM		18:00
7:00 PM		19:00
8:00 PM		20:00
9:00 PM		21:00
10:00 PM		22:00
11:00 PM		23:00

Specifications (Nominal)

General Pump Specifications

Used to test the pump:

- CADD™ medication cassette reservoirs, [REF](#) 21-7002
- CADD® extension sets, [REF](#) 21-7045
- CADD® administration sets, [REF](#) 21-7091
- CADD® high volume administration sets, [REF](#) 21-7055

Resolution

- CADD™ medication cassette reservoir: 0.050 mL per pump stroke nominal
- CADD® administration set: 0.050 mL per pump stroke nominal
- CADD® high volume administration set: 0.1 mL per pump stroke nominal

Size

4.1 cm × 10.2 cm × 12.7 cm (1.6 in × 4 in × 5 in), excluding cassette and accessories

Weight

595 g (21 oz) including 4 AA alkaline batteries, excluding other accessories

Pump alarms

High priority alarms:

- Air in line detected
- Battery depleted while delivering
- Battery removed while delivering
- Battery unusable while delivering
- Disposable attached improperly
- Disposable detached while delivering
- Disposable locked but not latched
- Disposable type high flow administration set not allowed
- Disposable type high flow administration set required
- Disposable type invalid
- Downstream occlusion
- Key stuck
- Pressure sensor faulty
- Pump automatically stopped
- Rechargeable battery end of life
- Remote dose cord key stuck
- Reservoir volume empty
- Stop mode reminder
- Upstream occlusion

Medium priority alarms: 23

Low priority alarms: 10

Informational messages/alerts: 23

Battery fallout alarm

Alarm sounds for 2 minutes if the pump has been powered up for a minimum of 2 minutes.

NOTE: Alarm enabled while pump is in run mode only.

Battery status

Battery State	CADD®-Solis Pump Status
25% to 100%	No alarm
Low battery	<ul style="list-style-type: none"> • Transition to low battery condition • Battery low message appears • Pump emits 3 beeps every 5 min • Low battery warning message appears on pump display • Pump is operable • LCD backlight flashes for 12 ms during each motor operation
Depleted battery	<ul style="list-style-type: none"> • Transition to depleted battery condition • Battery depleted message appears • Pump emits a continuous, variable-tone alarm • Depleted battery warning message appears on pump display • Battery power is too low to operate pump • Pump delivery operation stops
Shut down	Pump shuts off due to too low operating voltage.

Alkaline battery life with screen backlight intensity set to 3

These estimates are based on laboratory tests conducted at room temperature using new batteries. Actual battery life varies depending on the battery brand, shelf life, temperature conditions, delivery rate, and frequency of screen display and backlighting. It is recommended that new batteries be kept available for replacement.

PCA and Continuous delivery

Delivery Rate (mL/hr)	Operating Time (hrs)	Vol. Delivered (mL)
0.4	260	104
1	210	210
5	165	825
10	150	1500
30	70	2100
50	42	2100
100	30	3000
200	22	4400
350	15	5250
500	7	3500

Intermittent delivery

Volume (mL)	Duration (hr)	Cycle (hr)	KVO (mL/hr)	Operating Time (hrs)	Vol. Delivered (mL)
20.2	1	4	0.2	85	430
40.7	1	4	0.2	72	733
61	1	6	0.2	91	925
125	1	6	0.2	70	1460
200	1	12	0.2	52	1300

Step delivery

Volume to Deliver (mL)	Starting Rate (mL/hr)	Step Duration (min)	Step Rate Increase (mL/hr)	Max Rate (mL/hr)	No. of Steps	Operating Time (hrs)
900	45	15	45	315	7	21
1500	37.5	30	80	300	5	20
2500	30	30	90	300	4	18

Taper delivery

Volume (mL)	Period (hr)	Taper Up (hr)	Taper Down (hr)	KVO (mL/hr)	Operating Time (hrs)	Vol. Delivered (mL)
2000	12	1	1	5	29	4800
3000	12	1	1	5	17	4220

Rechargeable battery pack life with screen backlight intensity set to 3

These estimates are based on laboratory tests conducted at room temperature using a new CADD[®]-Solis rechargeable battery pack. Actual battery life varies depending on temperature conditions, delivery rate, and frequency of screen display and backlighting. It is recommended that new batteries be kept available for replacement.

PCA and Continuous delivery

Rate (mL/hr)	Life (hrs)	Volume (mL)
0.4	119	47.6
1	89	89
5	88	440
10	86	860
30	43	1290
50	30	1500
100	19	2375
200	16	3200
350	11	3850
500	8	8000

Intermittent delivery

Volume (mL)	Duration (hr)	Cycle (hr)	KVO (mL/hr)	Operating Time (hrs)	Vol. Delivered (mL)
20.2	1	4	0.2	54	273
40.7	1	4	0.2	51	520
61	1	6	0.2	54	550
125	1	6	0.2	43	896
200	1	12	0.2	48	800

Step delivery

Volume to Deliver (mL)	Starting Rate (mL/hr)	Step Duration (min)	Step Rate Increase (mL/hr)	Max Rate (mL/hr)	No. of Steps	Operating Time (hrs)
900	45	15	45	315	7	17
1500	37.5	30	80	300	5	14
2500	30	30	90	300	4	14

Taper delivery

Volume (mL)	Period (hr)	Taper Up (hr)	Taper Down (hr)	KVO (mL/hr)	Operating Time (hrs)	Vol. Delivered (mL)
2000	10	1	1	5	20	3360
3000	10	1	1	5	16	3950

ClassificationCF , Class II **Moisture protection**

Splashproof (IPX4) per IEC 60529

Maximum infusion pressure

1.86 bar (27.0 psi)

Maximum time to occlusion alarm and Bolus volume at occlusion alarm

Flow Rate (mL/hr)	Tubing Set	Maximum Time to Occlusion		Bolus at Occlusion	
		Raw Test Data (minutes)	Specification (minutes)	Raw Test Data (mL)	Specification (mL)
0.1	CADD™ medication cassette reservoir REF 21-7002 with CADD® extension set REF 21-7045	90	≤ 160	0.107	≤ 0.25
	CADD® administration set REF 21-7091	122	≤ 190	0.139	≤ 0.30
	CADD® high volume administration set REF 21-7055	1140	≤ 1200	1.250	≤ 1.40
Flow Rate (mL/hr)	Tubing Set	Maximum Time to Occlusion		Bolus at Occlusion	
		Raw Test Data (seconds)	Specification (seconds)	Raw Test Data (mL)	Specification (mL)
150	CADD™ medication cassette reservoir REF 21-7002 with CADD® extension set REF 21-7045	4	≤ 45	0.069	≤ 0.25
	CADD® administration set REF 21-7091	4	≤ 45	0.072	≤ 0.30
	CADD® high volume administration set REF 21-7055	32	≤ 90	1.059	≤ 1.40

Power sources

- AC adapter
- CADD®-Solis rechargeable battery pack
- Four AA (IEC LR6) alkaline batteries

Charging system for internal memory backup battery

The internal memory backup battery uses lithium manganese dioxide technology. It charges whenever the pump is powered on and has a 10-month memory capacity once it has been charged for 250 hours at 20°C (68°F).

System operating temperature

2°C to 40°C (36°F to 104°F)

System storage and transportation temperature

-20°C to 60°C (-4°F to 140°F)

Relative humidity

20% to 90% relative humidity, non-condensing

Atmospheric pressure

70 kPa to 106 kPa (10.2 psi to 15.4 psi)

System delivery accuracy**WARNING:**

- **Ensure that the $\pm 6\%$ system delivery accuracy specification is taken into account when programming the pump and/or filling the reservoir. Failure to do so may result in medication in the reservoir becoming depleted sooner than expected. If the pump is being used to deliver critical or life sustaining medication, the interruption in the delivery of medication could result in patient injury or death.**
- **System delivery inaccuracies may occur as a result of back pressure or fluid resistance, which depends upon drug viscosity, catheter size, and extension set tubing (for example, microbore tubing), and placing the infusion reservoir and/or pump above or below the level of the patient. System delivery inaccuracy may result in under or overdelivery of medication.**

$\pm 6\%$ (nominal). At low infusion rates, this accuracy may not be achieved for short periods. During the total infusion time, the accuracy averages out.

System definition

CADD®-Solis pump with:

- An attached medication cassette reservoir and CADD® extension set, *or*
- An attached medication cassette reservoir with Flow Stop feature and CADD® extension set, *or*
- An attached CADD® administration set *or*
- An attached CADD® administration set with Flow Stop feature

High pressure alarm threshold

18 ± 9 psi (1.24 bar \pm 0.62 bar)

Air detector alarm**Sensitivity:**

- Low: Single bubble > 400 μ L
- High: Single bubble > 150 μ L

Accumulated Air:

Greater than 1 mL air over 15 minutes (nominal)

Bolus accuracy specification: $\pm 6\%$

Actual test data for bolus accuracy at 0.05 mL:

Average	0.0508 mL
% Error	1.6%
Minimum Error %	-3.0%
Maximum Error %	4.2%

Actual test data for bolus accuracy at 20 mL:

Average	20.108 mL
% Error	0.5%
Minimum Error %	-2.9%
Maximum Error %	4.5%

Maximum volume infused under single-fault conditions

CADD[®] administration set: 0.15 mL

CADD[®] high volume administration set: 0.30 mL

Delivery rate during priming

Approx. 250 mL/hr with standard volume cassette tubing

500 mL/hr with high flow volume cassette tubing

Alarm disabled during priming

Air-In Line

Delivery Specifications

Common Delivery Specifications

Reservoir volume

0 to 9999

Programmable in 1 mL increments

Displayed in 0.1 mL increments

Given

0 to 99,999.99 in 0.01 unit increments

Res vol low trip point

1 to 999 mL in increments of 1 mL

Res vol empty alarm

- Insistent and one time only
- Non-insistent and repeating

Delayed start

0 to 96 hr in 5 min increments

Pump stopped alarm

Informational or high priority

Air detector

On or off

Air detector sensitivity

Low sensitivity: Single bubble > 400 μ L

High sensitivity: Single bubble > 150 μ L

Alarm volume

High, medium, or low

PM (preventive maintenance) reminder

Interval: 1 to 24 months in 1 month increments

Enable: On or off

Custom keypad code

001 to 899 in increments of 1

Custom clinician code

001 to 899 in increments of 1

Custom administrator code

001 to 899 in increments of 1

Date format

US standard (mm/dd/yy)

European standard (dd/mm/yy)

Japanese standard (yy/mm/dd)

Time format

00:00 to 23:59 military or 12-hour AM/PM

Downstream occlusion sensitivity

High Sensitivity: When the high pressure alarm threshold is reached, the downstream occlusion alarm is triggered immediately.

Low Sensitivity: When the high pressure alarm threshold is reached, the downstream occlusion alarm is delayed for 2 seconds. This allows for the pressure to stabilize before a possible alarm. If the pressure stabilizes below the high pressure alarm threshold before the 2 second delay is complete, the alarm will not occur.

Upstream occlusion sensor

On or off

NOTE: The upstream occlusion sensor is automatically disabled during use with medication cassette reservoirs.

PCA Delivery Specifications**Programming units**

If programming through manual mode. Otherwise the programming units are preset through the CADD[®]-Solis Medication Safety Software.

- Milliliters (mL)
- Milligrams (mg)
- Micrograms (mcg)

Concentration**mg/mL:**

- 0.1 to 0.5 mg/mL in increments of 0.1 mg/mL
- 0.5 to 1 mg/mL in increments of 0.5 mg/mL
- 1 to 15 mg/mL in increments of 1 mg/mL
- 15 to 100 mg/mL in increments of 5 mg/mL

mcg/mL:

- 1 to 15 mcg/mL in increments of 1 mcg/mL
- 15 to 100 mcg/mL in increments of 5 mcg/mL
- 100 to 500 mcg/mL in increments of 100 mcg/mL

Continuous rate

0 to 100 mL/hr (or the mg or mcg equivalent)

PCA dose

0 mL to 50 mL (or the mg or mcg equivalent)

Max delivery rate (continuous rate + PCA dose) programmable from 40 to 250 mL/hr

PCA dose lockout

1 minute to 24 hours in the following increments:

- 1 minute for values between 1 and 20 minutes
- 5 minutes between 20 minutes and 24 hours

Max doses per hour

1 to 60

Delivery limit amount

0.1 to 1900 mL (or the mg or mcg equivalent) in increments of:

- 0.01 mL from 0.01 to 0.5 mL
- 0.5 mL from 0.5 to 100 mL
- 1.0 mL from 100 to 1,000 mL

Clinician bolus

0 mL to 50 mL (or mg or mcg equivalent)

Delivery limit method

- Delivery limit
- Max doses per hour
- Not in use

Delivery limit period

1 to 12 hours in increments of 1 hour

Maximum delivery rate, combined bolus and continuous rate

40 to 250 mL/hr in increments of 1 mL

KVO rate

0 mL/hr or 0.1 mL/hr

Continuous Delivery Specifications

Continuous rate

0.1 to 500 mL/hr

KVO rate

0.1 to 10 mL/hr

Intermittent Delivery Specifications

Dose volume

0.1 to 1000 mL

Dose duration

1 min to 24 hr

Dose cycle

10 min to 96 hr

Next dose start time

0 to 96 hr in 5 min increments

KVO rate

0 to 10 mL/hr

Step Delivery Specifications**Initial rate**

0.4 to 499 mL/hr

Plateau rate

0.4 to 500 mL/hr

Rate increment

0.4 to 499 mL/hr

Infusion volume

1 to 9990 mL

Step duration

15 min to 24 hr

KVO rate

0 to 10 mL/hr

Step infusion alerts

On or Off

Taper Delivery Specifications**Infusion volume**

1 to 9990 mL

Taper up

0 min to 99:40 hr:min

Taper down

0 min to 99:40 hr:min

Infusion duration

10 min to 99:50 hr:min

KVO rate

0 to 10 mL/hr

Plateau rate upper limit

0.1 to 500 mL/hr

Electromagnetic Emissions and Immunity Declarations

Guidance and Manufacturer's Declaration—Electromagnetic Emissions		
The CADD®-Solis pump is intended for use in the electromagnetic environment specified below. The customer or the user of the CADD®-Solis pump should assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment—Guidance
RF emissions CISPR 11	Group 1	The Pump uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The Pump is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	


Compliance using:

- 100-240 VAC 50/60Hz to 7VDC switching desktop AC adapter, with an AC power cord length of 1.8 m (6 feet).
- Rechargeable battery pack.
- Remote dose cord with a length of 152 cm \pm 5 cm (60 in. \pm 2 in.).
- USB cable length of less than 2 m (6.5 feet).

WARNING:

- The use of power supplies and a remote dose cord other than those listed in the electromagnetic emissions declaration may result in increased emissions or decreased immunity of the pump.
- The pump should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, verify normal operation of the pump in the configuration in which it is to be used.

Guidance and Manufacturer's Declaration—Electromagnetic Immunity			
The CADD®-Solis pump is intended for use in the electromagnetic environment specified below. The customer or user of the CADD®-Solis pump should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment—Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	\pm 6 kV contact \pm 8 kV air	\pm 8 kV contact \pm 15 kV air (IEC 60601-2-24)	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	\pm 2 kV for power supply lines \pm 1 kV for input/output lines	\pm 2 kV for power supply lines Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	\pm 1 kV line(s) to line(s) \pm 2 kV line(s) to earth	\pm 1 kV line(s) to line(s) \pm 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations on power supply input lines IEC 61000-4-11	< 5% U_t (> 95% dip in U_t) for 0.5 cycle 40% U_t (60% dip in U_t) for 5 cycles 70% U_t (30% dip in U_t) for 25 cycles < 5% U_t (> 95% dip in U_t) for 5 sec	< 5% U_t (> 95% dip in U_t) for 0.5 cycle 40% U_t (60% dip in U_t) for 5 cycles 70% U_t (30% dip in U_t) for 25 cycles < 5% U_t (> 95% dip in U_t) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Pump requires continued operation during power mains interruptions, it is recommended that the Pump be powered from an uninterruptible power supply or a battery.
Power frequency 50/60 Hz magnetic field IEC 61000-4-8	3 A/m	400 A/m (IEC 60601-2-24)	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE: U_t is the a.c. mains voltage prior to application of the test level.			

Guidance and Manufacturer's Declaration—Electromagnetic Immunity			
The CADD®-Solis pump is intended for use in the electromagnetic environment specified below. The customer or user of the CADD®-Solis pump should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment—Guidance
			Portable and mobile RF communications equipment should be used no closer to any part of the Pump, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz outside ISM bands ^a	3 Vrms	Recommended separation distance $d = 1.2 \sqrt{P}$
Conducted RF IEC 61000-4-6	10 Vrms 150 kHz to 80 MHz in ISM bands ^b	10 Vrms	Recommended separation distance $d = 1.2 \sqrt{P}$
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.5 GHz	10 V/m	$d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz
			Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). ^c Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^c should be less than the compliance level in each frequency range. ^d Interference may occur in the vicinity of equipment marked with the following symbol: 
NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people			
^a The ISM (industrial, scientific, and medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.			
^b The compliance levels in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2.5 GHz are intended to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas. For this reason, an additional factor of 10/3 is used in calculating the recommended separation distance for transmitters in these frequency ranges.			
^c Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the CADD®-Solis pump is used exceeds the applicable RF compliance level above, the CADD®-Solis pump should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the CADD®-Solis pump.			
^d Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the CADD®-Solis pump				
The CADD®-Solis pump is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the CADD®-Solis pump can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the CADD®-Solis pump as recommended below, according to the maximum out put power of the communications equipment.				
Rated maximum output power or transmitter W	Separation distance according to frequency of transmitter m			
	150 kHz to 80 MHz outside ISM bands $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz in ISM bands $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.12	0.23
0.1	0.38	0.38	0.38	0.73
1	1.2	1.2	1.2	2.3
10	3.8	3.8	3.8	7.3
100	12	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be determined using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: The ISM (industrial, scientific, and medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.

NOTE 3: An additional factor of 10/3 is used in calculating the recommended separation distance for transmitters in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2.5 GHz to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas.

NOTE 4: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.



Collect Separately

This product contains electrical and electronic components (including batteries) that may contain materials, which if disposed of with general waste, could be damaging to the environment.

In accordance with Directive 2002/96/EC Waste Electrical and Electronic Equipment, residents of the European Union must follow specific disposal or recycling instructions for this product. Contact your local distributor, or visit the following web site for specific instructions:

<http://www.smiths-medical.com/recycle/index.html>.

Non-European Union residents must dispose of or recycle this product (including batteries) in accordance with the local laws or regulations that apply.

WARNING: There are potential health hazards associated with improper disposal of batteries, electronics, and contaminated (used) reservoirs and extension sets. Dispose of used batteries, reservoirs, extension sets, and other used accessories, or a pump that has reached the end of its useful life, in an environmentally safe manner, and according to any regulations that may apply.

Menu Maps

Tasks Menus

Tasks

Reset Reservoir Volume
Taper Down Now (Taper only)
Set Delayed Start (not Intermittent)
Set Next Dose Start Time (Intermittent only)
Prime Tubing
View Delivery Settings
Display and Sound Settings
Change Time and Date
View Reports
View Advanced Tasks

Delivery Settings

PCA	Continuous	Intermittent	Step	Taper
Continuous Rate	Continuous Rate	Dose Volume	Infusion Volume	Infusion Volume
PCA Dose		Dose Duration	Initial Rate	Taper Up
PCA Lockout		Dose Cycle	Rate Increment	Taper Down
Delivery Limit		Dose Rate	Plateau Rate	Infusion Duration
Max Doses/Hr		Next Dose	Step Duration	Plateau Rate
		Infusion Duration		
KVO Rate	KVO Rate	KVO Rate	KVO Rate	KVO Rate
Reservoir Vol.	Reservoir Vol.	Reservoir Vol.	Reservoir Vol.	Reservoir Vol.
Start Time	Start Time		Start Time	Start Time

Display and Sound Settings

Backlight Intensity
Alarm Volume
Sound Theme
Key Beep On/Off
Numeric Format

Time and Date

Current Time
Current Date
Time Format
Date Format

Reports

Total Given (<i>not</i> PCA)
Given and PCA Dose Counters (PCA only)
PCA Dose Graph (PCA only)
Delivery History and Pie Chart
Delivery Log
Event Log
Protocol Library Summary
Device Information

Advanced Tasks

Step Down (Step only)
Step Up (Step only)
Give Clinician Bolus (PCA only)
Patient Permissions
Air and Occlusions Settings
Alarm Settings
Security Settings
Start New Patient
Start New Protocol Same Patient
Delivery Hard and Soft Limits
Reset to Factory Settings

Advanced Tasks Menus

Patient Permissions

Priming Security On/Off
Delayed Start Security On/Off

Air and Occlusions Settings

Air Detector On/Off
Air Detector Sensitivity
Upstream Sensor On/Off
Downstream Sensor Sensitivity

Alarm Settings

Infusion Alerts (Step only)
Pump Stopped Alarm Type
Reservoir Low Trip Point
Reservoir Low Alarm Type
Reservoir Empty Alarm Type
PM Reminder On/Off
PM Reminder Interval

Security Settings

Keypad Security
Keypad Code
Clinician Code
Admin Code
Manual Programming Security

Delivery Hard and Soft Limits

Mode	Settings
PCA	Continuous Rate PCA Dose PCA Lockout Delivery Limit Method Delivery Limit Period Delivery Limit Amount Max Doses/Hr Max Delivery Rate Clinician Bolus Amount
Continuous	Continuous Rate
Intermittent	Dose Volume Dose Duration Dose Cycle
Step	Initial Rate Rate Increment Plateau Rate
Taper	Infusion Volume Max Plateau Rate

Default Factory Settings

The first time you use the pump, the protocol is set to the factory default, which is the **Continuous mode**. You may reset the pump to the factory default settings at any time (see *Reset to Factory Settings* on page 109).

Factory Default Pump Settings (Continuous Mode)

Submenu	Default Setting
Reset Reservoir Volume	100 mL
View Delivery Settings	(See <i>Manual Mode Initial Delivery Settings</i> , on page 141)
Display and Sound Settings	
Backlight Intensity	3
Alarm Volume	High
Sound Theme	Standard
Key Beep On/Off	On
Numeric Format	1,234.56
Change Time and Date	
Time Format	12:00:00 AM/PM
Date Format	Month/Day/Year
View Advanced Tasks	
Patient Permissions	
Priming Security On/Off	On
Delayed Start Security On/Off	On
Air and Occlusions Settings	
Air Detector On/Off	On
Air Detector Sensitivity	Low
Upstream Sensor On/Off	On
Downstream Sensor Sensitivity	Low
Alarm Settings	
Pump Stopped Alarm Type	Informational
Reservoir Low Trip Point	5 mL
Reservoir Low Alarm Type	Insistent and One Time Only
Reservoir Empty Alarm Type	One Time Only
PM Reminder On/Off	Off
Security Settings	
Keypad Security	Code only
Keypad Code	201
Clinician Code	997
Manual Programming Security	Admin Code
Continuous Mode Delivery Settings Hard and Soft Limits	(See <i>Manual Mode Initial Delivery Hard and Soft Limits</i> , on page 142)

Manual Mode Initial Settings

Manual Mode Initial Delivery Settings

Therapy	Default Setting
All Therapies	
Reservoir Vol.	0 mL
PCA	
Continuous Rate	0 mL/hr
PCA Dose	0 mL
PCA Lockout	1 Hrs 0 Min
KVO Rate	0.1 mL/hr
Give Clinician Bolus	0 mL/hr
Continuous	
Continuous Rate	0.1 mL/hr
KVO Rate	0 mL/hr
Intermittent	
Dose Volume	0.1 mL
Dose Duration	2 Hrs 0 Min
Dose Cycle	24 Hrs 0 Min
Dose Rate	0.1 mL/hr
Next Dose	[current time]
KVO Rate	0 mL/hr
Step	
Infusion Volume	1 mL
Initial Rate	0.4 mL/hr
Rate Increment	0.4 mL/hr
Plateau Rate	0.4 mL/hr
Step Duration	30 Min
Infusion Duration	2 Hrs 30 Min
KVO Rate	0 mL/hr
Taper	
Infusion Volume	24 mL
Taper Up	0 Min
Taper Down	0 Min
Infusion Duration	24 Hrs 0 Min
Plateau Rate	1 mL/hr
KVO Rate	1 mL/hr

Manual Mode Initial Delivery Hard and Soft Limits

Therapy	Default Settings			
PCA				
Continuous Rate	Hard Max:	100 mL/hr	Soft Min:	0 mL/hr
	Soft Max:	0 mL/hr	Hard Min:	0 mL/hr
PCA Dose	Hard Max:	50 mL/hr	Soft Min:	0 mL/hr
	Soft Max:	0 mL/hr	Hard Min:	0 mL/hr
PCA Lockout	Hard Max:	24 Hrs 0 Min	Soft Min:	24 Hrs 0 Min
	Soft Max:	24 Hrs 0 Min	Hard Min:	1 Min
Delivery Limit Method	Not in use			
Max Delivery Rate	125 mL/hr			
Clinician Bolus	Hard Max:	50 mL/hr	Soft Min:	0 mL/hr
	Soft Max:	0 mL/hr	Hard Min:	0 mL/hr
Continuous				
Continuous Rate	Hard Max:	500 mL/hr	Soft Min:	0.1 mL/hr
	Soft Max:	0.1 mL/hr	Hard Min:	0.1 mL/hr
Intermittent				
Dose Volume	Hard Max:	1,000 mL/hr	Soft Min:	0.1 mL/hr
	Soft Max:	0.1 mL/hr	Hard Min:	0.1 mL/hr
Dose Duration	Hard Max:	24 Hrs 0 Min	Soft Min:	24 Hrs 0 Min
	Soft Max:	24 Hrs 0 Min	Hard Min:	1 Min
Dose Cycle	Hard Max:	96 Hrs 0 Min	Soft Min:	96 Hrs 0 Min
	Soft Max:	96 Hrs 0 Min	Hard Min:	10 Min
Step				
Initial Rate	Hard Max:	499 mL/hr	Soft Min:	0.4 mL/hr
	Soft Max:	0.4 mL/hr	Hard Min:	0.4 mL/hr
Rate Increment	Hard Max:	499 mL/hr	Soft Min:	0.4 mL/hr
	Soft Max:	0.4 mL/hr	Hard Min:	0.4 mL/hr
Plateau Rate	Hard Max:	500 mL/hr	Soft Min:	0.4 mL/hr
	Soft Max:	0.4 mL/hr	Hard Min:	0.4 mL/hr
Taper				
Infusion Volume	Hard Max:	9,990 mL	Soft Min:	1 mL
	Soft Max:	1 mL	Hard Min:	1 mL
Max Plateau Rate	500 mL/hr			

Accuracy Test Results

In this device, as with all infusion pumps, the motion of the pumping mechanism and variations in individual disposables cause short-term fluctuations in rate accuracy. The following curves show typical performance of the pump system in two ways:

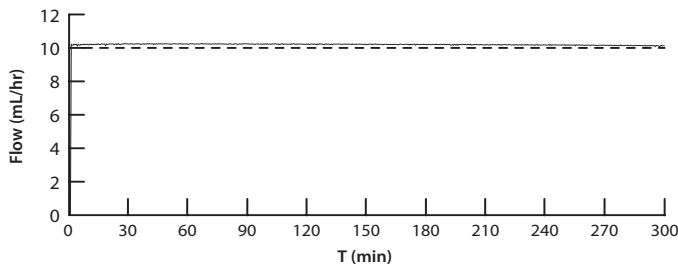
1. A flow versus time graph during the stabilization period (start-up curves).
2. The accuracy of fluid delivery of particular time periods, or observation windows, is measured (trumpet curves).

The start-up curve displays flow rate continuously from the start of the infusion. The curve visually represents flow rate uniformity. Trumpet curves are derived from the sixth hour of this data. The tests were performed per IEC 60601-2-24 standard.


Over long observation windows, short term fluctuations have minimal effect on accuracy as represented by the flat part of the curve. As the observation window is reduced, short term fluctuations have an increasing effect as represented by the “mouth” of the trumpet. Being aware of system accuracy over various observation windows may be of interest when certain drugs are being administered. Short term fluctuations in rate accuracy may have clinical impact depending on the half-life of the drug being infused, both the trumpet curve and drug half-life should be taken into consideration.

Start-up curve over the stabilization period

Low flow rate (10 mL/hr)




Intermediate flow rate (150 mL/hr)

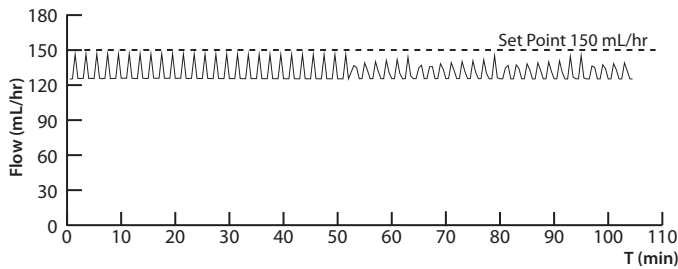


Trumpet curve over T(2) period

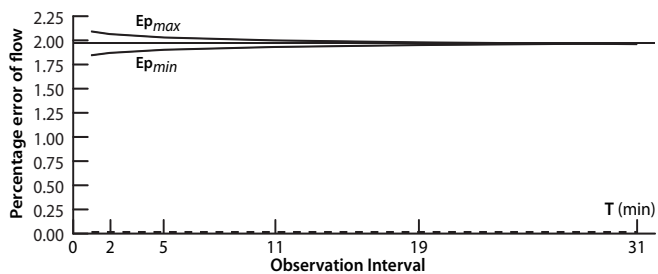Low rate (10 mL/hr)



Intermediate rate (150 mL/hr)

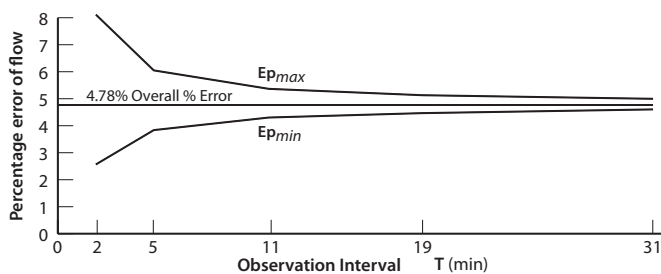

Limited Warranty

Smiths Medical ASD, Inc. (the “Manufacturer”) warrants to the Original Purchaser that the CADD®-Solis Ambulatory Infusion Pump (“Pump”), not including accessories, shall be free from defects in materials and workmanship under normal use, if used in accordance with this Operator’s Manual, for a period of two years from the actual date of sale to the Original Purchaser. THERE ARE NO OTHER WARRANTIES.

This warranty does not cover normal wear and tear and maintenance items, and specifically excludes batteries, administration sets, extension sets or any other accessory items or equipment used with the Pump.

Subject to the conditions of and upon compliance with this Limited Warranty, the Manufacturer will repair or replace at its option without charge (except for a minimal charge for postage and handling) any Pump (not including accessories) which is defective if a claim is made during such two-year period.

The following conditions, procedures, and limitations apply to the Manufacturer’s obligation under this warranty:

- A. Parties Covered by this Warranty:** This warranty extends only to the Original Purchaser of the Pump. This warranty does not extend to subsequent purchasers. The Original Purchaser may be a patient, medical personnel, a hospital, or institution which purchases the Pump for treatment of patients. The Original Purchaser should retain the invoice or sales receipt as proof as to the actual date of purchase.
- B. Warranty Performance Procedure:** Notice of the claimed defect must be made in writing or by telephone to the Manufacturer as follows: Customer Service Department, **Smiths Medical ASD, Inc.**, 1265 Grey Fox Road, St. Paul, MN 55112, (800) 258-5361 (USA, Canada) or **Smiths Medical International Ltd.** WD24 4LG, UK, +44 (0) 1923 246434. Notice to the Manufacturer must include date of purchase, model and serial number, and a description of the claimed defect in sufficient detail to allow the Manufacturer to determine and facilitate any repairs which may be necessary. Authorization must be obtained prior to returning the Pump. If authorized, the Pump must be properly and carefully packaged and returned to the Manufacturer, postage prepaid. Any loss or damage during shipment is at the risk of the sender.
- C. Conditions of Warranty:** The warranty is void if the Pump has been 1) repaired by someone other than the Manufacturer or its authorized agent; 2) altered so that its stability or reliability is affected; 3) misused; or, 4) damaged by negligence or accident. Misuse includes, but is not limited to, use not in compliance with the Operator’s Manual or use with nonapproved accessories. The Pump is a sealed unit, and the fact that the seal has been broken will be considered conclusive evidence that the Pump has been altered or misused. Removal or damage to the Pump’s serial number will invalidate this warranty.
- D. Limitations and Exclusions:** Repair or replacement of the Pump or any component part thereof is the EXCLUSIVE remedy offered by the Manufacturer. The following exclusions and limitations shall apply:
1. No agent, representative, or employee of the Manufacturer has authority to bind the Manufacturer to any representation or warranty, expressed or implied.
 2. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS OR USE OF THE PUMP FOR ANY PARTICULAR PURPOSE.
 3. The Pump can only be used under the supervision of medical personnel whose skill and judgment determine the suitability of the Pump for any particular medical treatment.
 4. All recommendations, information, and descriptive literature supplied by the Manufacturer or its agents are believed to be accurate and reliable, but do not constitute warranties.
- E. Computer Program License:**
1. The Pump is intended to be used in conjunction with a particular Licensed Computer Program supplied by Manufacturer and use of any other program or unauthorized modification of a Licensed Computer Program shall void Manufacturer’s warranty as set forth above.
 2. The Original Purchaser and any users authorized by the Original Purchaser are hereby granted a nonexclusive, nontransferable license to use the Licensed Computer Program only in conjunction with the single Pump supplied by Manufacturer. The Licensed Computer Program is supplied only in machine-readable object code form and is based upon Manufacturer’s proprietary confidential information. No rights are granted under this license or otherwise to decompile, produce humanly readable copies of, reverse engineer, modify or create any derivative works based upon the Licensed Computer Program.
 3. All other terms and conditions of this Limited Warranty shall apply to the Licensed Computer Program.
- The Manufacturer disclaims responsibility for the suitability of the Pump for any particular medical treatment or for any medical complications resulting from the use of the Pump. The Manufacturer shall not be responsible for any incidental damages or consequential damages to property, loss of profits, or loss of use caused by any defect or malfunction of the Pump.

This warranty gives the Original Purchaser specific legal rights, and the Original Purchaser may have other legal rights which may vary from state to state.

Index

A

AC adapter 15
 Accessories 15
 Accuracy 129, 130, 143
 AC power jack 12
 Admin code 18, 24, 131
 Advanced tasks 99, 138, 139
 Air and occlusion settings 102
 Alarm settings 105
 Delivery hard and soft limits 33
 Give clinician bolus 43
 Menu maps 138, 139
 Patient permissions 100
 Reset to factory settings 109
 Security settings 18, 22
 Start new patient 27
 Start new protocol same patient 27
 Step down 68
 Step up 67
 Air and occlusion settings 102, 139
 Air detector on/off 102
 Air detector sensitivity 103, 130
 Downstream sensor sensitivity 104, 131
 Upstream sensor on/off 103
 Air detector 14
 On/Off 102
 Sensitivity 103, 129
 Alarms and messages 110
 Alphabetical list 112
 Help screens 111
 High priority alarm 110, 125
 Informational message 111
 Low priority alarm 111
 Medium priority alarm 110
 System fault alarm 110
 Troubleshooting 112
 Alarm settings 105, 139
 Infusion alerts on/off 105
 PM interval 108
 PM reminder 108
 PM reminder on/off 107
 Pump stopped alarm type 105
 Reservoir empty alarm type 107
 Reservoir low alarm type 106
 Reservoir low trip point 106
 Alarm volume 89
 Amber indicator light 12, 13, 80, 86, 110
 Analgesics 9
 Anesthetics 9
 Autolock 21

B

Backlight intensity 89
 Batteries 78, 129
 AA batteries 14, 16, 78, 79
 AC adapter 12, 15, 78
 Battery compartment 12, 14, 79
 Battery life 126
 Battery status 16, 126
 Installing batteries 79
 Rechargeable battery pack 14, 15, 78, 79
 Symbols 10, 11, 16

C

CADD[®]-Solis Medication Safety Software 15, 17, 22, 26, 27, 33, 34, 80, 99
 Cassette 12, 81
 Attach a cassette 82, 87
 Remove a cassette 81
 Cassette/Keypad lock 10, 12, 14
 Cassette latch 12, 14
 Cautions 6
 Cleaning the pump 119
 Clinician bolus 43, 46
 Clinician code 18, 24, 131
 Color 17
 Amber 12, 13, 17, 34, 80, 86, 110
 Blue 12, 14, 17, 51, 80, 111
 Gray 17, 37, 47, 51, 59, 69, 80
 Green 12, 13, 17, 34, 37, 47, 51, 59, 69, 80, 85
 Olive 17, 59
 Purple 17, 37
 Red 17, 37, 47, 51, 59, 69, 86, 110
 Yellow 17
 Continuous rate 38, 40, 47, 49
 Continuous therapy 47
 Continuous rate 47
 Default values, manual program 29
 Delayed start 31
 Delivery specifications 132
 Home screen 47
 KVO rate 47
 Programming example 48
 Reservoir vol. 48
 Set delayed start 31
 Start time 48
 Contraindications 3, 9
 Counters, given and PCA dose 96
 Current date 93
 Current time 16, 92
 Current value, delivery limit 33

D

Date, current 93
 Date format 94, 131
 Daylight savings time 94
 Default factory settings 140
 Delayed start 31, 85, 130
 Delayed start security on/off 101
 Delivery hard and soft limits 33, 139, 142
 Delivery history and pie chart 97
 Delivery limits 28, 33, 34, 38, 139, 142
 Delivery log 97
 Delivery settings 34, 137, 141

- Continuous rate 38, 47
- Delivery limit 38
- Dose cycle 52
- Dose duration 52
- Dose rate 52
- Dose volume 52
- Edit delivery settings 35
- Infusion duration 61, 70
- Infusion volume 60, 70
- Initial rate 60
- KVO rate 39, 47, 53, 61, 71
- Max doses/hr 38
- Next dose 52
- PCA dose 38
- PCA lockout 38
- Plateau rate 60, 70
- Rate increment 60
- Reservoir vol. 39, 48, 53, 61, 71
- Review pump settings 36
- Specifications 130
- Start time 39, 48, 61, 71
- Step duration 60
- Taper down 70
- Taper up 70
- View delivery settings 35

 Device information 98
 Display and backlight 12, 13, 16
 Display and sound settings 89, 137

- Alarm volume 89
- Backlight intensity 89
- Key beep on/off 91
- Numeric format 91
- Sound theme 90

 Dose cycle 52, 55
 Dose duration 52, 55
 Dose rate 52, 55
 Dose volume 52, 54
 Downstream occlusion sensor 14

- Sensitivity 104, 131

 Drug 16, 26

E

Edit delivery settings 35
 Epidural administration 9
 Event log 98
 Examples *See* Programming examples

F

Factory settings, reset 109

G

Given and PCA dose counters 96
 Green indicator light 12, 13, 80, 85

H

Hard delivery limits 33, 142
 Help

- Screens 111
- Technical assistance 2

 High priority alarm 110, 125
 Home screen 16, 37, 47, 51, 59, 69

I

Indications 9
 Informational message 111
 Infusion alerts on/off 105
 Infusion duration 61, 65, 70, 74
 Infusion volume 60, 63, 70, 73
 Initial rate 60, 63
 Intermittent therapy 51

- Default values, manual program 29
- Delivery specifications 132
- Dose cycle 52
- Dose duration 52
- Dose rate 52
- Dose volume 52
- Home screen 51
- KVO rate 53
- Nest dose start time 32
- Next dose 52
- Programming example 53
- Reservoir vol. 53
- Reset a cycle 58
- Restart a dose 58
- Set delayed start 31
- Stop the pump during an infusion 58

K

Key *See* Pump key
 Key beep on/off 91
 Keypad 10, 12, 13
 Keypad code 18, 23, 131
 Keypad security 23
 KVO rate 39, 42, 47, 50, 53, 56, 61, 65, 71, 75

L

Last error code 98
 Latch 12, 14
 LCD *See* Display
 Libraries 26, 27
 Lights 12, 13
 Continuously on 13, 110
 Flashing 13, 80, 85, 86
 Limits *See* Delivery limits
 Lock 12, 14, 15, 16, 21
 Log
 Delivery 97
 Event 98
 Low priority alarm 111

M

Magnetic resonance imaging *See* MRI and radiation
 Manual mode, initial settings 141, 142
 Manual programming 28
 Protocol default values 29
 Security 25
 Max doses/hr 38, 42
 Medication Safety Software *See* CADD-Solis
 Medication Safety Software
 Medium priority alarm 110
 Menu maps 137
 Menu screen 16
 Messages 110
 Alphabetical list 112
 Help screens 111
 Informational message 111
 Troubleshooting 112
 MRI and radiation 120

N

New patient 27
 New protocol same patient 27
 New protocol selection 26
 Next dose 52, 56
 Next dose start time 32
 Numeric format 91

O

Occlusion and air settings 102, 139
 Air detector on/off 102
 Air detector sensitivity 103, 130
 Downstream sensor sensitivity 131
 Upstream sensor on/off 103
 Occlusion sensor
 Downstream 14
 Upstream 14

P

Patient permissions 100, 139
 Delayed start security on/off 101
 Priming security on/off 100
 PCA dose 10, 38, 41, 45, 46
 PCA dose counters 96
 PCA dose graph 96
 PCA lockout 38, 41
 PCA therapy 37
 Clinician bolus 43, 46
 Continuous rate 38, 40
 Default values, manual program 29
 Delayed start 31
 Delivery limit 38
 Delivery specifications 131
 Home screen 37
 KVO rate 39
 Max doses/hr 38
 PCA dose 38, 41, 45, 46
 PCA lockout 38, 41
 Programming example 39
 Programming screens 38
 Remote dose cord 15, 38, 45
 Reservoir vol. 39
 Scroll ranges 123
 Set delayed start 31
 Start time 39
 Stop clinician bolus 46
 Stop PCA dose 46
 Pie chart and delivery history 97
 Plateau rate 60, 64, 70, 75
 PM *See* Preventive maintenance (PM)
 Power down (off) 80
 Power jack 12, 14
 Power light 12
 Power switch 10, 12, 14
 Power up (on) 80
 Preventive maintenance (PM)
 PM interval 108
 PM reminder 108, 131
 PM reminder on/off 107
 Prime tubing 83, 130
 Priming security on/off 100

Programming examples
 Continuous therapy 48
 Intermittent therapy 53
 PCA therapy 39
 Step therapy 62
 Taper therapy 72

Programming screens
 Continuous therapy 47
 Intermittent therapy 52
 PCA therapy 38
 Step therapy 60
 Taper therapy 70

Programming the pump 26
 Manual programming 28
 Protocol selection from a library 27
 Single protocol download 26

Protocol library summary 98

Protocols 26
 Libraries 26, 98
 Manual programming 28
 New protocol selection 26
 Protocol selection from a library 27
 Single protocol download 26

Pump key 15, 22, 23

Pump screen 16

Pump stopped alarm type 105

Q

Qualifier 16, 26

R

Radiation and MRI 120

Rate increment 60, 64

Rechargeable battery pack 11, 14, 15, 16, 78, 79, 80, 127

Rechargeable Battery Pack 15

Remote dose cord 13, 15, 38, 45

Remote dose cord jack 12, 14

Reports 95, 138
 Delivery history and pie chart 97
 Delivery log 97
 Device information 98
 Event log 98
 Given and PCA dose counters 96
 PCA dose graph 96
 Protocol library summary 98
 Total given 95

Reservoir empty alarm type 107

Reservoir low alarm type 106

Reservoir low trip point 106, 130

Reservoir volume 39, 42, 48, 50, 53, 57, 61, 66, 71, 75, 87, 130

Reset reservoir volume 87

Reset to factory settings 109

Review pump settings 30, 36, 40, 49, 54, 62, 72

S

Same patient, new protocol 27

Screens and display 16

Security codes *See* Security settings

Security settings 18, 139
 Admin code 18, 24
 Autolock 21
 Clinician code 18, 24
 Customized settings 21, 22
 Entering security codes 21
 Keypad code 18, 23
 Keypad security 23
 Manual programming 25
 Tables 19, 20

Select drug 27, 40, 49, 54, 62, 72

Select qualifier 27, 30, 40, 48, 53, 62, 72

Select therapy 27, 30, 39, 48, 53, 62, 72

Select units 30

Serial number 98

Set delayed start 31

Set next dose start time 32

Single protocol download 26

Soft delivery limits 33, 142

Soft keys 10, 13

Software version 98

Sound theme 90

Specifications 125
 Delivery specifications 130
 General pump specifications 125

Start new patient 27, 30

Start new protocol same patient 27

Start the pump 85
 Daily infusion 76
 Delayed start 85
 Give a clinician bolus 43
 Next dose start time 85
 Power up 80
 Reset a cycle or infusion 58, 67
 Restart a dose or infusion 58, 67
 Start an infusion 67
 Start a PCA dose 45
 stop/start button 10

Start time 39, 48, 61, 71

Step down 68

Step duration 60, 65

Step therapy 59
 Default values, manual program 29
 Delayed start 31
 Delivery specifications 133
 Home screen 59
 Infusion duration 61
 Infusion volume 60
 Initial rate 60
 KVO rate 61
 Plateau rate 60
 Programming example 62

- Rate increment 60
- Reservoir vol. 61
- Restarting the pump during an infusion 67
- Set delayed start 31
- Starting each infusion 67
- Start time 61
- Step down 68
- Step duration 60
- Step up 67
- Step up 67
- Stopping during an infusion 67
- Stop the pump 86
 - Power down 80
 - Stop a clinician bolus 46, 86
 - Stop a PCA dose 46, 86
 - stop/start button 10
 - Stop the pump during an infusion 58
 - Taper down before stopping 77, 86
- Subarachnoid administration 9
- Symbols 10
 - Home screen 16, 37, 47, 51, 59, 69
- System fault alarm 110

T

- Taper down 70, 74, 77
- Taper down now 77
- Taper therapy 69
 - Default values, manual program 29
 - Delayed start 31
 - Delivery specifications 133
 - Home screen 69
 - Infusion duration 70
 - Infusion volume 70
 - KVO rate 71
 - Plateau rate 70
 - Programming example 72
 - Reservoir vol. 71
 - Set delayed start 31
 - Start time 71
 - Taper down 70
 - Taper up 70
- Taper up 70, 73
- Tasks menu 88
 - Change time and date 92
 - Display and sound settings 89
 - Menu maps 137
 - Next dose start time 32
 - Prime tubing 83
 - Reset reservoir volume 87
 - Set delayed start 31
 - Set next dose start time 32
 - Taper down now 77
 - View advanced tasks 99
 - View delivery settings 35
 - View reports 95
- Therapy 16, 26

- Time and date 92, 137
 - Current date 93
 - Current time 16, 92
 - Date format 94, 131
 - Daylight savings time 94
 - Military time 124
 - Time format 93, 131
- Total given 95
- Troubleshooting 112

U

- Ultrasound 120
- Unlock 21, 22
- Upstream occlusion sensor 14
 - On/Off 103
- USB port 12, 14

V

- View delivery settings 35
- Volume
 - Alarm 89, 130
 - Reservoir 39, 42, 48, 50, 53, 57, 61, 66, 71, 75, 87, 130

W

- Warnings 3

CADD™

Smiths Medical ASD, Inc.

St. Paul, MN 55112 USA
1 800.258.5361 (USA)
1 214.618.0218
www.smiths-medical.com

EC REP

Smiths Medical International Ltd.

WD24 4LG, UK
+44 (0) 1923 246434

CE
0473

AUS REP

Smiths Medical Australasia Pty. Ltd.

61 Brandl Street
Eight Mile Plains, QLD 4113, Australia
+61 (0) 7 3340 1300

CADD-Solis, CADD, and the CADD and Smiths Medical design marks are trademarks of the Smiths Medical family of companies. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are trade names, trademarks, or service marks of their respective owners.

© 2010 Smiths Medical family of companies. All rights reserved.

2010-05
40-5904-50A

smiths medical